



## Küschall® Compact 2.0

en User Manual

This manual MUST be given to the user of the product. BEFORE using this product, this manual MUST be read and saved for future reference

**Küschall®**  
UNLIMIT YOUR WORLD

# Contents

|                                                                              |           |
|------------------------------------------------------------------------------|-----------|
| <b>1 General</b>                                                             | <b>3</b>  |
| 1.1 Introduction                                                             | 3         |
| 1.2 Delivery Check                                                           | 3         |
| 1.3 Symbols in this Document                                                 | 3         |
| 1.3.1 Other Symbols                                                          | 3         |
| 1.4 Warranty Information                                                     | 3         |
| 1.5 Service Life                                                             | 3         |
| 1.6 Limitation of Liability                                                  | 3         |
| <b>2 Safety</b>                                                              | <b>4</b>  |
| 2.1 Safety Information                                                       | 4         |
| 2.1.1 Specific Safety Requirements for Carbon Line Models                    | 5         |
| 2.2 Safety Devices                                                           | 6         |
| 2.3 Labels and symbols on the product                                        | 6         |
| <b>3 Product overview</b>                                                    | <b>6</b>  |
| 3.1 Product description                                                      | 6         |
| 3.2 Intended use                                                             | 6         |
| 3.3 Main parts of the wheelchair                                             | 7         |
| 3.4 Parking brakes                                                           | 7         |
| 3.5 Backrest                                                                 | 8         |
| 3.5.1 Standard backrest cover                                                | 8         |
| 3.5.2 Adjusting the tension of the backrest cover                            | 8         |
| 3.5.3 Backrest height / backrest angle                                       | 8         |
| 3.5.4 Push handles                                                           | 8         |
| 3.5.5 Foldable push handles                                                  | 8         |
| 3.5.6 Height-adjustable push handles                                         | 9         |
| 3.6 Folding backrest (optional)                                              | 9         |
| 3.7 Stabilizing bar for the backrest                                         | 9         |
| 3.8 Armrests and side parts                                                  | 9         |
| 3.8.1 Tubular armrest, non-locking, height-adjustable, removable, swivelling | 9         |
| 3.8.2 Side rest, stepless height-adjustable, non-locking                     | 10        |
| 3.8.3 Kuschall armrest                                                       | 10        |
| 3.8.4 Mudguard                                                               | 10        |
| 3.8.5 Clothes-guard                                                          | 11        |
| 3.9 Leg rests                                                                | 11        |
| 3.10 Seat cushion                                                            | 11        |
| <b>4 Options</b>                                                             | <b>11</b> |
| 4.1 Anti-tip device                                                          | 11        |
| 4.2 Posture belt                                                             | 12        |
| 4.3 Attachment when using positioning belt                                   | 13        |
| 4.4 Headrest                                                                 | 13        |
| 4.5 Transit wheels                                                           | 13        |
| 4.6 Tipper aid                                                               | 14        |
| 4.7 Swing away table tray                                                    | 14        |
| 4.7.1 Swing away angle adjustable leg rests                                  | 15        |
| 4.8 Amputee support                                                          | 16        |
| 4.9 Dual handrim DHR                                                         | 16        |
| 4.10 Cane holder                                                             | 17        |
| 4.11 Calf strap / heel loop                                                  | 17        |
| 4.12 Passive illumination                                                    | 17        |
| 4.13 Tension adjustable seat upholstery                                      | 17        |

|                                                          |           |
|----------------------------------------------------------|-----------|
| <b>5 Setup</b>                                           | <b>17</b> |
| 5.1 Safety Information                                   | 17        |
| <b>6 Using the wheelchair</b>                            | <b>17</b> |
| 6.1 Safety Information                                   | 17        |
| 6.2 Braking during use                                   | 17        |
| 6.3 Getting in and out of the wheelchair                 | 18        |
| 6.4 Driving and steering the wheelchair                  | 18        |
| 6.5 Negotiating steps and kerbs                          | 19        |
| 6.6 Going up and down stairs                             | 20        |
| 6.7 Negotiating tamps and slopes                         | 20        |
| 6.8 Stability and balance when seated                    | 20        |
| 6.9 Transporting the occupied wheelchair in a vehicle    | 21        |
| 6.10 Transporting the manual wheelchair without occupant | 23        |
| <b>7 Transport</b>                                       | <b>23</b> |
| 7.1 Safety information                                   | 23        |
| 7.2 Folding and unfolding the wheelchair                 | 23        |
| 7.3 Removing and installing the rear wheels              | 24        |
| <b>8 Maintenance</b>                                     | <b>24</b> |
| 8.1 Safety Information                                   | 24        |
| 8.2 Maintenance schedule                                 | 24        |
| 8.3 Cleaning and disinfection                            | 25        |
| 8.3.1 General safety information                         | 25        |
| 8.3.2 Cleaning intervals                                 | 25        |
| 8.3.3 Cleaning                                           | 26        |
| 8.3.4 Disinfection                                       | 26        |
| <b>9 Troubleshooting</b>                                 | <b>27</b> |
| 9.1 Safety Information                                   | 27        |
| 9.2 Identifying and repairing faults                     | 27        |
| <b>10 After use</b>                                      | <b>28</b> |
| 10.1 Storage                                             | 28        |
| 10.2 Reconditioning                                      | 28        |
| 10.3 Disposal                                            | 28        |
| <b>11 Technical data</b>                                 | <b>28</b> |
| 11.1 Dimensions and weight                               | 28        |
| 11.2 Maximum weight of removable parts                   | 29        |
| 11.3 Tyres                                               | 30        |
| 11.4 Materials                                           | 30        |
| 11.5 Environmental conditions                            | 30        |

# 1 General

## 1.1 Introduction

This user manual contains important information about product handling. To ensure safety when using the product, read the user manual carefully and follow safety instructions.

Only use this product if you have read and understood this manual. Seek additional advice from a healthcare professional who is familiar with your medical condition and clarify any questions regarding the correct use and necessary adjustment with the healthcare professional.

There may be sections in this document, which are not relevant to your product. This document applies to all available models (on the date of printing). If not otherwise stated, each section in this document refers to all product models.

The models and configurations available in your country can be found in the country-specific sales documents. We reserve the right to alter product specifications without further notice.

Before reading this document, make sure you have the latest version. You find the latest version as a PDF on the Invacare website.

If you find that the font size in the printed document is difficult to read, you can download the PDF version from the website. The PDF can then be scaled on screen to a font size that is more comfortable for you.

For more information about the product, for example product safety notices and product recalls, contact your Invacare distributor. See addresses at the end of this document.

In case of a serious incident with the product, you must inform the manufacturer and competent authority in your country.

## 1.2 Delivery Check

Any transport damage must be reported immediately to the transport company. Remember to keep the packaging until the transport company has checked the goods and a settlement has been reached.

## 1.3 Symbols in this Document

Symbols and signal words are used in this document and apply to hazards or unsafe practices which could result in personal injury or property damage. This document is printed in greyscale. For your information, the safety messages have the following colour coding according to ANSI Z535.6: Danger (Red), Warning (Orange), Caution (Yellow) and Notice (Blue). See the information below for definitions of the signal words.

|                                                                                     |                                                                                                                      |
|-------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------|
|  | <b>WARNING</b><br>Indicates a hazardous situation that could result in serious injury or death if it is not avoided. |
|  | <b>CAUTION</b><br>Indicates a hazardous situation that could result in minor or slight injury if it is not avoided.  |
|  | <b>NOTICE</b><br>Indicates a hazardous situation that could result in damage to property if it is not avoided.       |



### Tips and Recommendations

Gives useful tips, recommendations and information for efficient, trouble-free use.

### 1.3.1 Other Symbols

(Not applicable for all manuals).



#### UK Responsible Person

Indicates if a product is not manufactured in the UK.



#### Triman

Indicates recycling and sorting rules (only relevant for France).

## 1.4 Warranty Information

We provide a manufacturer's warranty for the product in accordance with our General Terms and Conditions of Business in the respective countries.

Warranty claims can only be made through the provider from whom the product was obtained.

(missing or bad snippet)(missing or bad snippet)

## 1.5 Service Life

The expected service life of this product is five years when used daily and in accordance with the safety instructions, maintenance intervals and correct use, stated in this manual.

The effective service life can vary according to frequency and intensity of use.

## 1.6 Limitation of Liability

Invacare accepts no liability for damage arising from:

- Non-compliance with the user manual
- Incorrect use
- Natural wear and tear
- Incorrect assembly or set-up by the purchaser or a third party
- Unauthorised modifications, combinations and/or use of unsuitable spare parts

## 2 Safety

### 2.1 Safety Information

This section contains important safety information for the protection of the wheelchair user and assistant and for safe, trouble-free use of the wheelchair.



**WARNING!**

**Risk of serious injury or damage**

Improper use of this product may cause injury or damage.

- If you are unable to understand the warnings, cautions or instructions, contact a health care professional or provider before attempting to use this equipment.
- Do not use this product or any available optional equipment without first completely reading and understanding these instructions and any additional instructional material such as user manual, service manual or instruction sheet supplied with this product or optional equipment.



**WARNING!**

**Risk of death or serious injury**

In case of fire or smoke, wheelchair occupants are at particular risk of death or serious injury, when they are not able to move away from the source of fire or smoke. Lighted matches, lighter and cigarettes can cause an open flame in the wheelchair surroundings or on clothes.

- Avoid using or storing the wheelchair near open flames or combustible products.
- Do not smoke while using the wheelchair.



**WARNING!**

**Risk of accidents and serious injury**

Accidents with resulting serious injury can occur if the wheelchair is improperly adjusted

- Contact your wheelchair provider in order to a qualified technician does the adjustments needed.



Invacare expects that the qualified technician is familiar with the product, with good technical knowledge to understand and follow the steps of the described instructions in this manual and equipped with proper tools.



**WARNING!**

**Risk of overturning**

The longitudinal position of the rear wheels axis of the wheelchair compared to the backrest position can affect its stability.

- Moving the rear wheels axis frontwards makes the wheelchair less stable, increasing the risk of tipping backwards. But it improves the manoeuvrability with a better grip position of the handrim and a short turning radius.
- Conversely, by moving the rear wheels axis backwards, the wheelchair is more stable and tilts less easily, but its maneuverability is reduced.
- Depending on the user's abilities and its particular safety limits, the decrease in stability can be compensated for by installing an anti-tip device.



**WARNING!**

**Risk of tilting**

The rear wheel axle position and the angle of your wheelchair's backrest are two of the key adjustments that can affect your stability.

- The changes of rear/front wheels position and/or angle fork adjustments must only be performed by a qualified technician after the assessment of your healthcare professional.



**WARNING!**

**Risk due to driving style being unsuitable for environmental the conditions**

There is a risk of skidding on wet ground, gravel or uneven terrain.

- Always adjust your speed and driving style to the conditions (weather, surface, individual ability, etc.).



**WARNING!**

**Risk of injury**

In a collision you could sustain injury to parts of your body that extend beyond the wheelchair (e.g. feet or hands).

- Avoid an unbraked collision.
- Never drive into an object head-on.
- Drive carefully through narrow passages.



**WARNING!**

**Risk due to wheelchair being out of control**

At high speed you could lose control of your wheelchair and overturn.

- Always adjust your driving speed to the given conditions.
- Avoid collisions in general.

**WARNING!****Risk of accidents**

Safety devices that incorrectly set or no longer working (brakes, anti-tip devices) can cause accidents.

- Always check that the safety devices are working before you use the wheelchair and have them regularly checked by a qualified technician or your provider.

**CAUTION!****Risk of injuries**

- In case of pressure sore or injured skin, protect your injure to avoid a direct contact with the fabrics of the device. Refer to a healthcare professional for medical advice.

**CAUTION!****Risk of injuries**

The use of a manual wheelchair can lead to shoulder pain. These recommendations reduce this effect:

- Move the rear wheels further forwards so arms can move within the anatomical range.
- Make sure that the sits in an ergonomically seat position as the disability allows.

**CAUTION!****Risk of burning**

The wheelchair components can heat up when exposed to external sources of heat:

- Do not expose the wheelchair to strong sunlight before use.
- Before usage, check temperature of all components that come into contact with your skin.

**CAUTION!****Risk of catching body parts in moving parts**

You risk getting parts of the body e.g. fingers or arms, caught in the moving parts of the wheelchair.

- Take attention while activating and deactivating the mechanisms of moving parts, such as the removable axle of the rear wheel, folding backrest or anti-tip device.

**CAUTION!****Risk of tipping**

The hanging of additional load (back pack or similar items) onto your chair backposts can affect the rearward stability of your wheelchair.

- Invacare strongly recommend the use of anti-tip device (available as an option) when using your backposts with additional loa.

**CAUTION!****Risk of injury**

Non-original or wrong parts may affect the function and safety of the product.

- Only use original parts for the product in use.
- Due to regional differences, refer to your local Invacare catalogue or website for available options or contact your Invacare distributor. See addresses at the end of this document.

**NOTICE!**

When combining this product with another product, the restrictions of both products apply for the combination. E.g. the maximum user weight of one product in the combination may be lower.

- Only use combinations allowed by Invacare. Contact your Invacare distributor for information.
- Before use, read the user manual of each product and check the restrictions. The full compliance with the information, instructions and recommendations contained in both user manuals is necessary for a safe use.

**NOTICE!****Risk when using an add-on to support or substitute propelling the wheelchair**

For safety reasons, it must only be operated by persons who have the physical and mental ability to safely operate the wheelchair with the add-on attached to it in all possible situations (e.g. road traffic) and, in the event of the add-on failing to work or shutting down, are able to brake the wheelchair and stop safely. Any use contrary to the intended use involves a risk of injury.

At high speed you could lose control of your wheelchair and overturn.

- Drive carefully.
- Always adapt your speed and your driving style to the conditions (weather, surface, individual ability, etc.).
- Avoid collision.

An add-on to support or to substitute propelling the wheelchair increases the strain on the wheelchair structure.

- After each use of an add-on, check the wheelchair for frame cracks or loose screwing.
- If any sign of fatigue on the wheelchair structure, stop using the wheelchair immediately and contact a qualified technician.

### 2.1.1 Specific Safety Requirements for Carbon Line Models

**WARNING!****Risk of breakage****Risk of injury**

- You must not attach any add-on propulsion system to the carbon tubes comprising the wheelchair frame.
- Never drill holes in any of the carbon tubes.

**CAUTION!****Risk of damage**

There is a risk of deteriorating or damaging the carbon tubes of the wheelchair. We strongly recommend the following:

- Do not over tighten any attached objects (water bottle holder, telephone holder etc.) on a carbon tube.
- Avoid contact between the carbon and metallic parts. Consider using a soft material (e.g. rubber or silicone) between these 2 materials.

2.2 Safety Devices



**WARNING!**  
**Risk of accidents**

Safety devices that are incorrectly set or no longer working (brakes, anti-tip device) can cause accidents.


- Always check that the safety devices are working before you use the wheelchair and have them checked regularly by a qualified technician.

The functions of the safety devices are described in 3.1 Product description, page 6

2.3 Labels and symbols on the product

Identification label

The identification label is attached to the frame of your wheelchair and gives the following information:

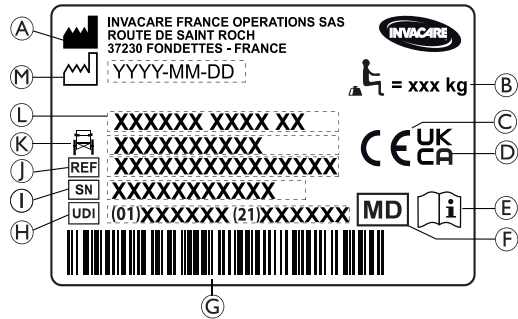

- A Manufacturer’s name and address
- B Maximum load of the wheelchair
- C European conformity mark
- D United Kingdom conformity mark
- E Consult user manual
- F Medical device
- G UDI barcode
- H Unique Device Identifier, with bar code
- I Serial Number
- J Reference number
- K Seat width
- L Product name
- M Date of manufacture

|  |                                                                                                    |
|--|----------------------------------------------------------------------------------------------------|
|  | Maximum load of the wheelchair or Maximum occupant mass without any additional load                |
|  | Consult user manual                                                                                |
|  | Indicates a hazardous situation that could result in serious injury or death if it is not avoided. |

Non-locking side parts warning label



Do not lift the wheelchair by non-locking side parts.

Posture belt information label



Do not lift the wheelchair by non-locking side parts.

Reclining mechanical backrest information label



Push on the back support tubes (1) before operating the levers (2), this is to release the auto-locked security system. This label is attached to the upper part of the backrest cane (right side).

Snap hook symbols

Depending on the configuration, some wheelchairs may be used as a seat in a motor vehicle, some may not.



Tie-down positions where the restraint system straps must be placed in case of transporting the occupied wheelchair in a motor vehicle. This symbol is only attached to the wheelchair when it is ordered with the transportation kit option.



**WARNING!**  
**Risk of serious injury or death**

Depending on the configuration, some wheelchairs may be used as a seat in a motor vehicle, some may not.

- Do not use the wheelchair for seat in a motor vehicle in.

3 Product overview

3.1 Product description

The **Küschall Compact SA** is a folding active wheelchair with cross struts and swing-away leg rests.

The **Küschall Compact FF** is a folding active wheelchair with cross struts.



**NOTICE!**

The wheelchair is manufactured and configured individually to the specifications in the order. The specification must be performed by a healthcare professional/ clinician according to the occupant’s requirements and health condition.

- Consult a healthcare professional if you intend to adapt the wheelchair configuration.
- Any adaptation should be performed by a qualified technician.

3.2 Intended use

The wheelchair may be used indoors and outdoors on level ground and accessible terrain.

The wheelchair is intended to provide mobility to persons limited to a sitting position, who are propelling the wheelchair themselves most of the time.

The wheelchair is not intended for abusive or excessive activities, such as sports, that are not part of the intended (daily) use.



### Intended user

The wheelchair is intended for persons, ages 12 and up (adolescents and adults). The weight of the wheelchair occupant should not exceed the maximum user weight as stated in the Technical Data section and on the identification label.

The intended user is the wheelchair occupant and/or an assistant. The user should physically and mentally be able to use the wheelchair safely (e.g. to propel, steer, brake).


### Indications

- Considerable to complete impairment of mobility due to structural and/or functional damage to the lower extremities.
- Sufficient strength and gripping function of arms and hands to propel the wheelchair.

### Contraindications

There are no contra-indications known, when using the wheelchair as intended.

## 3.3 Main parts of the wheelchair

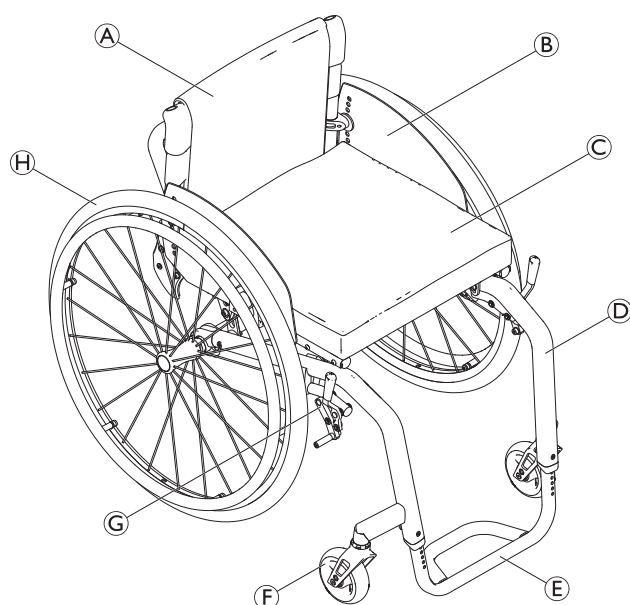

|     |                                        |
|-----|----------------------------------------|
| (A) | Backrest                               |
| (B) | Clothes-guard                          |
| (C) | Seat with cushion                      |
| (D) | Frame                                  |
| (E) | Footrest                               |
| (F) | Castor fork with castor wheel          |
| (G) | Parking brake                          |
| (H) | Rear wheel with handrim and quick-axle |



The equipment of your wheelchair may differ from the diagram as each wheelchair is manufactured individually to the specifications in the order.

## 3.4 Parking brakes

The parking brakes are used to immobilize the wheelchair when it is stationary to prevent it from rolling away.



### WARNING!

#### Risk of overturning if you brake sharply

If you engage the parking brakes while you are moving, the direction of movement can become uncontrollable and the wheelchair may stop suddenly, which can lead to a collision or to you falling out

- Never engage the parking brakes while you are moving.



### WARNING!

#### Risk due to wheelchair being out of control

- The parking brakes must be operated simultaneously.
- Do not engage the parking brakes to slow down the wheelchair.
- Do not lean on the parking brakes for support or transfer.



### WARNING!

#### Risk of overturning

- The parking brakes will not operate correctly unless there is sufficient air in the tyres. Ensure the correct tyre pressure, *11.3 Tyres, page 30*.



### CAUTION!

#### Risk of pinching or crushing


There may be a very small gap between the rear wheel or the armrest front socket and the parking brake with the risk that you could trap your fingers.

- Keep your fingers away from movable parts when using the brake, always keep your hand on the brake lever.



The distance between the brake shoe and the tyre can be adjusted. The adjustment must be carried out by a qualified technician.

### Standard brake

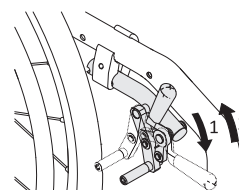


1. To engage the brake, push the brake lever forwards with the palm as far as possible.
2. To disengage the brake, pull the brake lever backwards with the low fingers.



The lever of the push handle can be folded back to facilitate transfers. To do so, pull up the lever and fold it backwards.


### Active brake

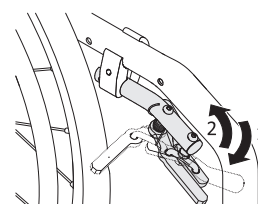

1. To engage the brake, pull the brake lever forwards beside or between your knees as far as possible.

- To disengage the brake, push the brake lever backwards beside or between your knees.

#### Pull brake

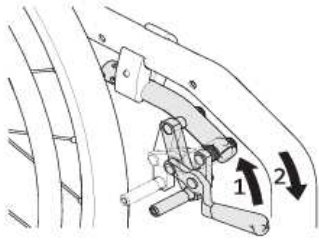

- To engage the brake, push the brake lever forwards with the palm as far as possible.
- To disengage the brake, pull the brake lever backwards with the low fingers.

#### Light active brake




##### CAUTION!

##### Risk of pinching or crushing

Do not grasp the brake mechanism with the full hand when engage or disengage the brake.



The hardness of the locking/unlocking could be adapted to your needs. The adjustment must be carried out by a qualified technician.

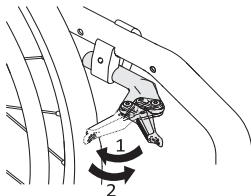

- To engage the brake, pull the brake lever backwards with the palm against the tyre as far as possible.
- To disengage the brake, push the brake lever forwards with the hand edge.

### 3.5 Backrest



##### WARNING!

##### Risk of injury to the assistant and damage to the wheelchair

Tilting the chair with a heavy occupant can hurt the assistant's back and damage the chair.

- Make sure to be able to safely control the wheelchair with a heavy occupant before performing a tilting maneuver.

#### 3.5.1 Standard backrest cover



A standard backrest cover is available for all backrest types. These covers are not adjustable.



##### WARNING!


##### Risk of tipping

If the standard backrest cover become slack, the tipping point of your wheelchair is changed for the worse.

- Slack standard backrest cover must be replaced by a qualified technician.

In addition to the standard backrest cover, there are tension-adjustable backrest covers available.

#### 3.5.2 Adjusting the tension of the backrest cover

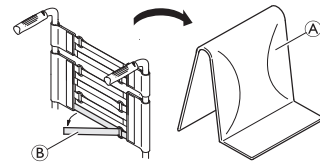

- Remove the backrest cushion (A).
- Loosen the hook-and-loop straps (B) on the rear of the backrest cover by simply pulling on them.
- Tighten or loosen the bands as required then reattach them.



##### WARNING!

##### Risk of tipping

If the bands are set very loose, the tipping point of your wheelchair is changed for the worse.

- Ensure with your healthcare professional that the bands are set correctly.



##### NOTICE!

##### Risk when using an add-on to support or substitute propelling the wheelchair

There is a risk that the geometry of the wheelchair is changed when the bands are overtightened.

- Do not overtighten the bands.
- Only ever tighten the bands when the wheelchair is unfolded.

#### 3.5.3 Backrest height / backrest angle

The height and the angle of the backrest can be adjusted. The adjustment must be carried out by a qualified technician.



##### WARNING!

##### Risk of tipping

The hanging of additional load (back pack or similar items) onto your wheelchair back posts can affect the rearward stability of your wheelchair. This can cause the chair to tip backwards causing injury.

- Try to avoid putting additional load in the back of the wheelchair.
- We strongly recommend the use of anti-tip devices when using your back posts with additional load.


#### 3.5.4 Push handles



##### NOTICE!

- Always check the push handles before using the wheelchair, as to whether the hand grips are secure, cannot be turned and cannot be pulled off.

#### 3.5.5 Foldable push handles

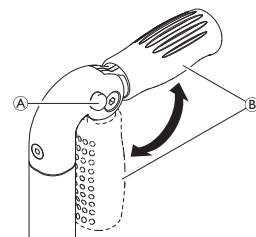



1. Fully depress button (A) and fold up respectively fold down the push handle (B) until it audibly engages.


**CAUTION!**

If not correctly engaged, the push handle could accidentally fold down while pushing the wheelchair.

- Make sure that the push handle is correctly engaged.

### 3.5.6 Height-adjustable push handles

The stepless height-adjustable push handles allow any assistant to adjust the handles to a comfortable level.

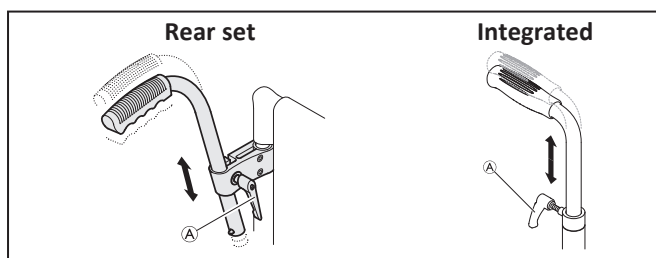

1. To set the height of the push handles, loosen the bolt (A), push the handle to the desired position and tighten the bolt again.

**WARNING!****Risk of injury**

- Do not use the height-adjustable push handles to pull the wheelchair.

### 3.6 Folding backrest (optional)

The backrest can be folded down and straightened up when the wheelchair is folded or unfolded.

#### Folding down the backrest

1. To fold down the backrest, grip the strap on the backrest of the wheelchair and pull it upwards.
2. At the same time, push the upper edge of the backrest and fold it upwards.

#### Folding up the backrest

1. To fold down the backrest, grip the strap on the backrest of the wheelchair and pull it upwards.

**NOTICE!**

- Do not fold and unfold the wheelchair when the backrest is folded down.

### 3.7 Stabilizing bar for the backrest

An optional stabilizing bar is available to increase the rigidity of the backrest on wheelchairs with large back heights or widths.


**NOTICE!**

- Do not use the stabilizing bar to lift or push the wheelchair. It could break.
- Release the stabilizing bar before folding the wheelchair.



Use a seat cushion with anti-slip underlay or hook fastening tapes to avoid slipping of the seat cushion. A loop fastening tape is pre-attached to the seat cover.

### Release the stabilizing bar

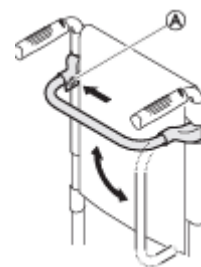

1. Push the spring (A) on the fastening on the left-hand side and swivel the stabilizing bar downwards.

The stabilizing bar can also be completely removed by pushing the spring pins on both sides. The stabilizing bar can then be taken off.

### 3.8 Armrests and side parts

**WARNING!****Risk of injury**


The side rests are not locked and can be easily pulled out upwards.

- Do not lift or transfer the wheelchair using the side rests.
- Do not use the side rests for transportation when carrying the wheelchair up- or downstairs.

**CAUTION!****Risk of pinching**

- Keep your fingers away from movable parts during removing, fitting or adjusting the armrest.

#### 3.8.1 Tubular armrest, non-locking, height-adjustable, removable, swivelling

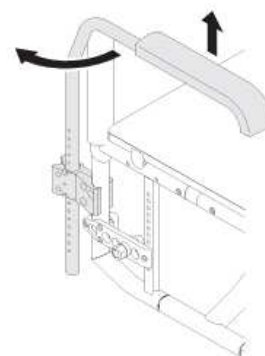

#### Removing

1. Raise the armrest and pull it upwards out of the holder.

#### Installing

1. Push the armrest into the holder.


#### Adjusting the height

1. Pull the armrest out of the holder.
2. Loosen the screw on the tube of the armrest and screw it back in at the desired height.
3. Push the armrest back into the holder.
4. Perform the setting on both sides.

## Swivelling

1. Slightly raise the armrest and swivel it outwards.

### 3.8.2 Side rest, stepless height-adjustable, non-locking

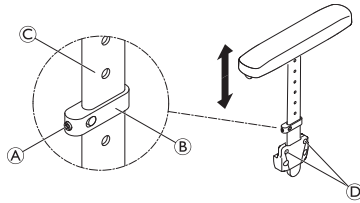

#### Removing

1. Pull the side rest by the arm cushion out of the holder.
2. To adjust how easy or difficult it is to pull the side rest out of the holder, change the degree to which the screws ⑤ are tightened.

#### Installing

1. Push the side rest into the holder.


#### Adjusting the height

1. Slightly loosen the grub screw ① on the clamping part ②.
2. Move the clamping part upwards or downwards on the profile ③ to the desired height.
3. Re-tighten the grub screw.
4. Perform the setting on both sides.

### 3.8.3 Kûschall armrest

This armrest option is locking in place and can be flipped backwards, swivelled outwards and removed. The T-armrest pad ① is height adjustable and also locks in place.

#### Removing the armrest

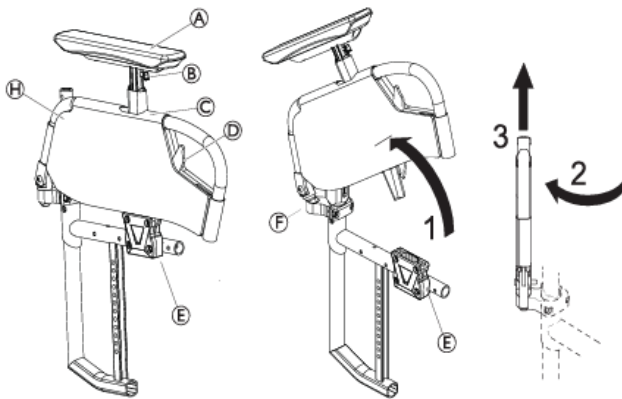

1. Push release lever ④ and swivel the armrest assembly upwards out of the armrest holder ⑤.
2. Pull the armrest assembly upwards out of the joint ⑥.
3. Swivel the armrest assembly 45° outwards.



Reverse the current procedure to reinstall the armrest.


#### Removing the cover

The cover ⑧ can be removed to be machine washed at 30°C.

1. Set the T-armrest in the lowest position pressing the ②.
2. Remove partially the cover ⑧.
3. Unscrew the screw ①.
4. Remove the T-armrest by pressing ③.



Reverse the current procedure to reinstall the cover.

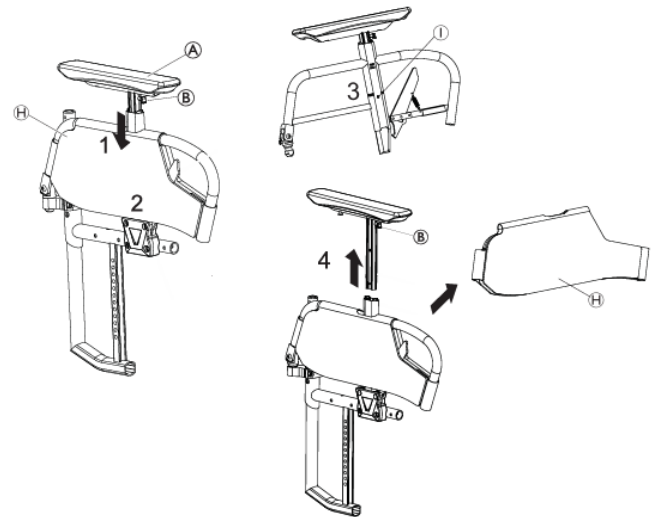

#### Adjusting the height

1. Press the ② and move the armrest pad ① to the required height.
2. Release the ② and allow it to engage with an audible "click" in the nearest hole.




#### WARNING! Risk of injury

Even if the armrests are locking, they are not designed to carry the weight of the occupied wheelchair.

- Never use the armrests to lift the wheelchair.
- Do not use the armrests for transportation when carrying the wheelchair up- or downstairs.

### 3.8.4 Mudguard

In order to prevent dirt from the wheels spraying upwards, a removable mudguard can be fitted. The position of the mudguard can be adjusted. The adjustment must be carried out by a qualified technician.

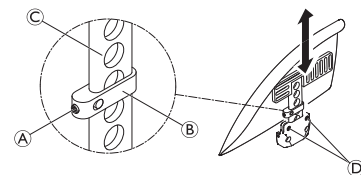

#### WARNING! Risk of injury

The mudguards are not locked and can be easily pulled out upwards.

- Do not lift or transfer the wheelchair using the mudguards.
- Do not use the mudguards for transportation when carrying the wheelchair up- or downstairs.

#### Removing

1. Pull the mudguard upwards out of the holder.
2. To adjust how easy or difficult it is to pull the mudguard out of the holder, change the degree to which the screws ⑤ are tightened.

## Installing


1. Push the mudguard into the holder.

## Adjusting the height

1. Slightly loosen the grub screw ① on the clamping part ②.
2. Move the clamping part upwards or downwards on the profile ③ to the desired height.
3. Re-tighten the grub screw.
4. Perform the setting on both sides.

## 3.8.5 Clothes-guard

A clothes-guard can be installed to protect your clothes. Before folding the back of the wheelchair, the clothes-guard must first be folded inwards against the backrest.

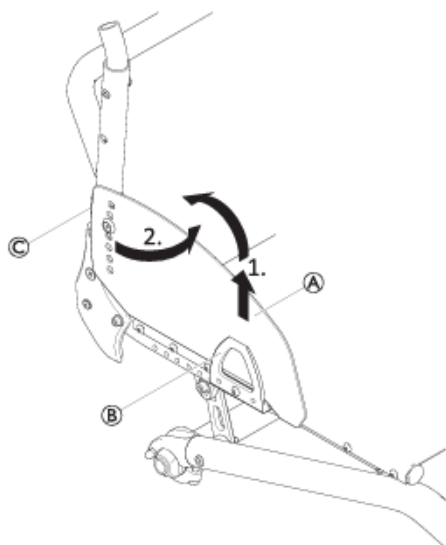

## Folding in


1. Raise the clothes-guard ① at the front out of the holder ② and swivel it inwards around the joint ③ in front of the backrest.

## 3.9 Leg rests



### WARNING! Risk of injury

- Never use the legrests or foot plates to lift the wheelchair.

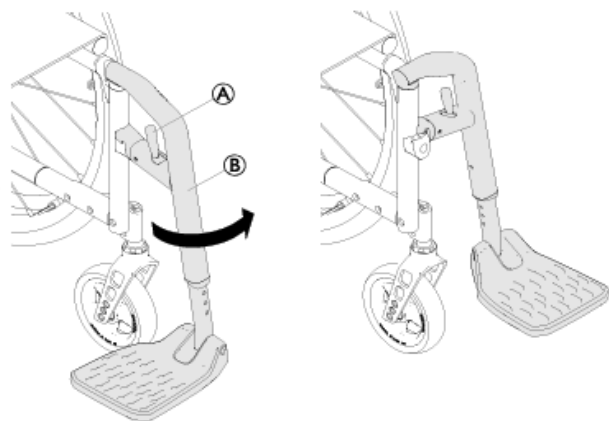

## Swivelling

1. Activate the release lever ① and swivel the leg rest ② to the outside.

## Swivelling forwards

1. Swivel the leg rest forwards until it engages.

## Unhinging

1. Activate the release lever A.
2. Pull the leg rest upwards.

## Hinging

1. Hinge the leg rest on the front of the frame and swivel it forwards until it engages.

## Adjusting the height

1. Remove screw, washer and nut from the frame B with 4 mm Allen key and 8 mm spanner.
2. Slide the footrests to the desired height.
3. Reinsert screw, washer and nut on the frame B and tighten them (4 Nm).
4. Perform the setting on both sides.
5. The distance between the lowest part of the footrest and the ground must be at least 50 mm. Invacare recommends that foot plate adjustment is carried out by a qualified technician.

## 3.10 Seat cushion

A suitable cushion is recommended to provide an even pressure distribution on the seat.



### CAUTION!

#### Risk of instability

Adding a cushion to the seat will raise your height above the ground and can affect your stability in all directions. If a cushion is changed it may also change the user's stability.

- If the thickness of the cushion is changed, a complete set up of the wheelchair need to be done by a qualified technician.
- We recommend using an Invacare or Matrx cushion with anti-slip base cover in order to prevent sliding.



Use a seat cushion with anti-slip underlay or hook fastening tapes to avoid slipping of the seat cushion. A loop fastening tape is pre-attached to the seat cover.

## 4 Options

### 4.1 Anti-tip device

An anti-tip device prevents the wheelchair from tipping backwards.



### WARNING!

#### Risk of overturning

Incorrectly set or malfunctioning anti-tip device can lead to overturning.

- Always check that the anti-tip device is working before using the wheelchair and have it set or readjusted by a qualified technician when required.
- In some configurations, the static stability of the wheelchair may be lower than 10°; Invacare strongly recommends the use of anti-tip device (available as an option).


**WARNING!****Risk of overturning**

On uneven or soft ground, the anti-tip device can sink into potholes or directly into the ground, thereby curtailing or eliminating its safety function.

- Only use the anti-tip device when travelling on even and firm ground.

**NOTICE!****Risk of damage**

- Never step on the anti-tip device or use it as a tipping aid.

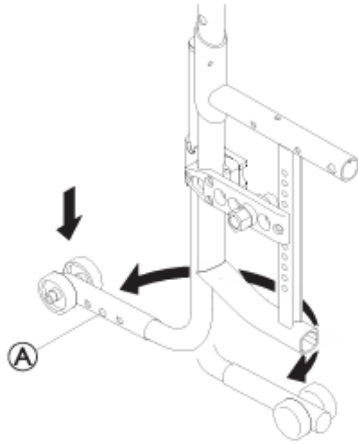**Activating the anti-tip device**

1. Push the anti-tip device completely downwards and swivel it backwards until it engages.

**WARNING!****Risk of tipping**

If the wheelchair is equipped with a retractable anti-tip device, it must always remain in the active position for user safety.

- It can be temporarily deactivated to overcome an obstacle that it cannot pass if it is in an active position.
- Never leave a user alone in their wheelchair with the anti-tip device deactivated.

**Deactivating the anti-tip device**

1. Push the anti-tip device completely downwards and swivel it forwards until it engages.

**Setting the position**

1. Press the spring pin (A) on the anti-tip device and push the outer section to the position required. Allow the spring pin to engage in the nearest hole.



If further adjustments of height and/or position are necessary, refer to a qualified technician.



**4.2 Posture belt**

The wheelchair can be equipped with a posture belt. It prevents the user from sliding downward in the wheelchair or from falling out of the wheelchair. The posture belt is not a positioning device.

**WARNING!****Risk of serious injury / strangulation**


A loose belt can cause the user to slip down and create a risk of strangulation.

- The posture belt should be mounted by a qualified technician and fitted by a healthcare professional.
- Always make sure that the posture belt is tightly fitted across the lower pelvis.
- Each time the posture belt is used, check if it fits properly. Changing the seat and / or backrest angle, the cushion and even your clothes influence the fit of the belt.

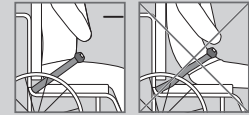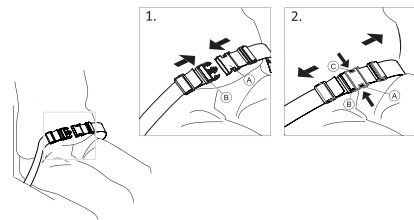
**Closing and opening the posture belt**

Ensure that you are sitting fully back in the seat and that the pelvis is as upright and symmetrical as possible.


1. To close, push the catch (A) into the buckle clasp (B).
2. To open, push the PRESS button (C) and pull the catch (A) out of the buckle clasp (B).

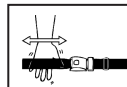**Adjusting the length**

The posture belt has the good length, when there's just sufficient space for a flat hand between body and belt.

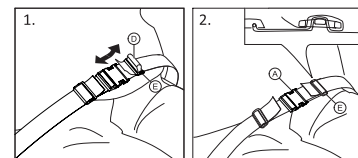

1. Shorten or extend the loop (D) as required.
2. Thread loop (D) through catch (A) and plastic buckle (E) until the loop is flat.
3. Ensure loop (D) is perfectly fitted in plastic buckle (E).
4. Secure adjustment with the strap extremity fitted in the buckle (E).

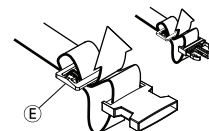

**WARNING!****Risk of sliding and strangulation / Risk of falling**

- Do the adjustments on both sides equally, so that the buckle clasp remains in a central position.
- Ensure that both strap extremities are threaded two times into the buckle ⑤ to avoid the belt from loosening.
- Make sure that the belts do not get caught in the spokes of a rear wheel.

### 4.3 Attachment when using positioning belt

**CAUTION!****Risk of sliding**

- The belt must be placed tight to the body.
- If the length of the belt has been adjusted to fit one type of seat cushion height, it is important that this cushion is used. If the seat cushion is replaced, the length of the belt needs to be adjusted.



Positioning belts which are CE-marked for the purpose of using on wheelchairs, can be mounted on the chair with preserved CE-marking. The positioning belt should be fitted by a healthcare professional and be mounted by a qualified technician.




The purpose of using a belt for positional, is to position the user and to give him / her a better posture.

### 4.4 Headrest



The height and the removal are operated by the handwheel. The bar is equipped with an adjustable stop block.

#### Height adjustment

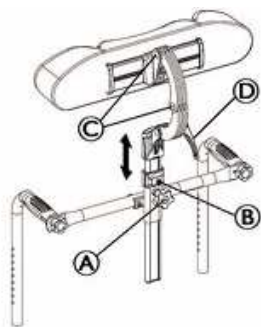


1. Loosen the screw ⑥ with 5 mm Allen key in the stop block.
2. Loosen the handwheel ①.
3. Adjust the headrest to the desired position.
4. Re-tighten the handwheel ① down to the top of the headrest attachment.
5. Slide the stop block.
6. Re-tighten the screw ⑥ (2–3 Nm).



You can now remove the headrest and reinsert it in the desired position without further adjustments.

#### Depth / angle adjustment

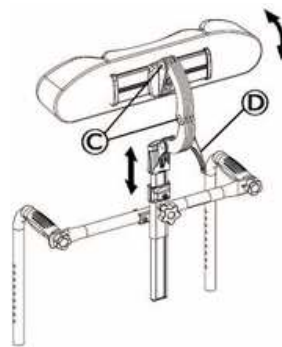


1. Loosen the handle ②.
2. Loosen the screws ③ with 5 mm Allen key.
3. Adjust the depth and angle of the headrest.
4. Re-tighten firmly the handle ② and the screws ③ (7 Nm).

### 4.5 Transit wheels

If your wheelchair is too wide for certain applications, e.g. trains / planes, narrow passageways, narrow doors, the transit wheels can be used.

**WARNING!****Risk of severe accidents**

When using transit wheels, the parking brakes no longer have any effect and you cannot control your wheelchair via the handrims.

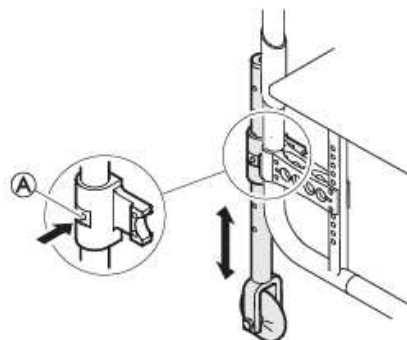

#### Switching to transit wheels

1. Lower the transit wheels into the operating position. To do this, press the pin ① inwards until you can pull the tube down.
2. Press the second pin in.
3. Pull the tube further downwards until the pin engages into the slot.
4. Move the wheelchair to a stable object (door frame, heavy table), which you can use for support and hold onto it firmly with one hand.
5. Tilt the wheelchair until the rear wheel on the opposite side rises off the floor.
  -  Prior to removing the rear wheels, the distance between the lowered transit wheels and the floor should not be more than 2 cm.
6. Remove the rear wheel → 7.3 Removing and installing the rear wheels, page 24
7. Repeat the process on the other side.




### Switching from transit wheels to rear wheels

1. Move the wheelchair to a stable object, which you can use for support and hold onto it firmly with one hand.
2. Tilt the wheelchair until the transit wheel on the opposite side rises slightly off the floor.
3. Install the rear wheel back onto the removable axle → 7.3 *Removing and installing the rear wheels, page 24*
4. Bring the transit wheels back into the inactive position.
5. Repeat the process on the other side.

## 4.6 Tipper aid

With the tipper aid, an assistant can tip the wheelchair more easily, in order to navigate steps for example.

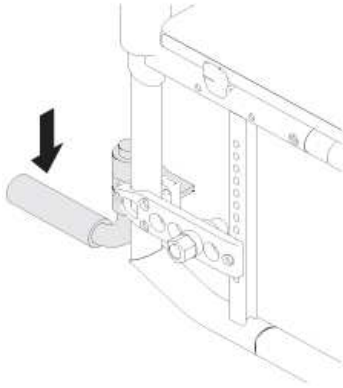

#### **WARNING!** Risk of overturning

- Ensure that the tipper aid does not project beyond the external diameter of the rear wheel.

1. Hold the wheelchair by the push handles.
2. Press the tipper aid with your foot and hold the wheelchair in the tipped position until you have overcome the obstacle.

## 4.7 Swing away table tray



#### **WARNING!** Risk of falling/injury

- The table tray must never be used as a replacement for the posture belt.



#### **WARNING!** Risk of tipping / injury

- Max. load on the table tray: 5 kg



#### **CAUTION!** Risk of discomfort/minor bruises

- Make sure that the elbows of the user are placed on the table tray when pushing the wheelchair. If the elbows protrude from the table tray while pushing the wheelchair, there is a risk of discomfort or minor bruises.




#### **WARNING!** Risk of injury/ bruises

- Take care not to hit the user whilst swinging the table tray.



#### **WARNING!** Risk of table tipping forwards


- Do not adjust the depth of the table tray too much forwards due to the risk of tipping the table when adding weight on it.

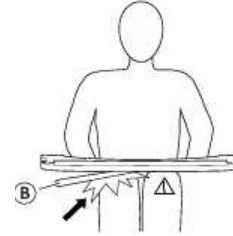

#### **CAUTION!** Risk of discomfort/minor bruises


- Make sure that the user has enough space for the legs. The table tube must not rest on the legs of the user.

### Depth adjustment

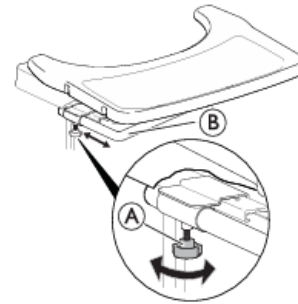


1. Loosen the 2 handwheels ①.
2. Adjust the table in depth with the table tube ②.
3. Re-tighten the 2 handwheels ①.

### Swing away function

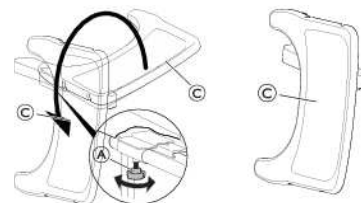

1. Loosen the 2 handwheels ①.
2. Swing the table tray ③ sideways.
3. Re-tighten the 2 handwheels ①.

### 4.7.1 Swing away angle adjustable leg rests

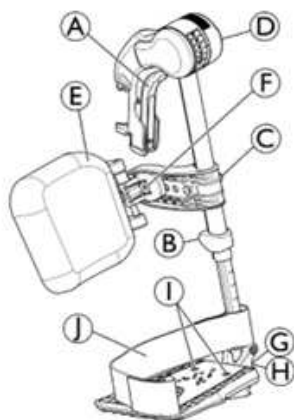

#### Swinging outwards

1. Activate the release lever A and swing the leg rests outwards.

#### Swinging forwards

1. Swing the leg rest forwards until it engages.


#### Unhinging

1. Activate the release lever A, swing the leg rest outwards.
2. Pull the leg rest upwards.

#### Hinging

1. Push the leg rest down into its receiver and swing it forwards until it engages.


#### Adjusting angle

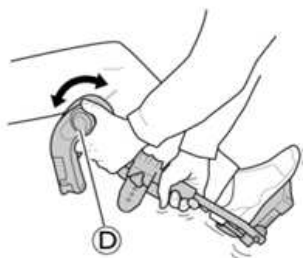

There are six preset positions available for angle adjustment.

1. Rotate the knob D with one hand while supporting the leg rest with your other hand.
2. When a suitable angle is obtained, let go off the knob and the leg rest will lock into the desired position.

#### Swinging calf pad and adjustments

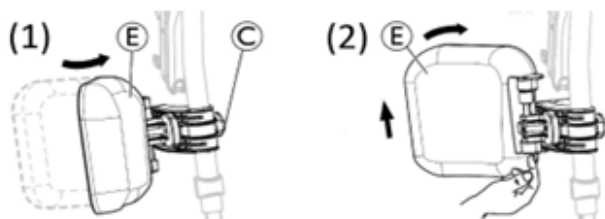


- Swinging calf pad forwards (1) or backwards (2):

1. Swing calf pad E forwards (1).
2. Pull up the calf pad E and swing it backwards (2).

#### Adjusting height or calf pad:

1. Loosen knob C.
2. Adjust to the desired height and firmly tighten the knob.

#### Adjusting depth of calf pad

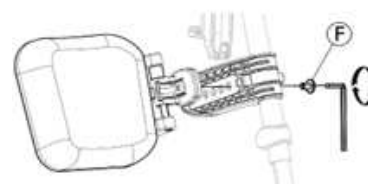

The calf pad has four depth adjustment options:

1. Remove the fixing screw F with 5 mm Allen key.
2. Adjust to one of the four positions and firmly tighten the fixing screw F (3–5 Nm).

#### Adjustable footplates

There are three different fold-up foot plates available:

##### Standard height-adjustable footplates:




1. Loosen the fixing screw B with 5 mm Allen key.
2. Adjust the height and let the screw catch one of the recesses on the foot plate tube.
3. Tighten the fixing screw B (3–5 Nm) in the desired position.




The distance between the lowest part of the footrest and the ground must be at least 50 mm.

##### Depth-, angle- (1) and width- (2) adjustable foot plates:

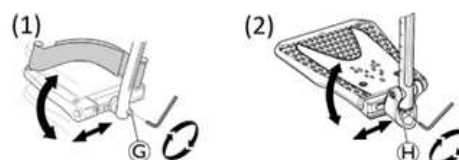

1. Loosen the fixing screw G or H with 5 mm Allen key.
2. Adjust the depth and angle positions and firmly tighten the fixing screw G (12 Nm) or H (8–9 Nm) in the desired position.

##### Width-adjustable foot plates (2):

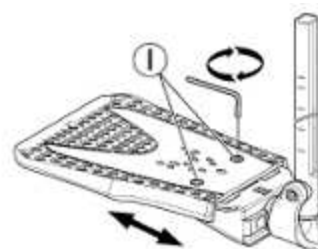

1. Loosen the two fixing screws I with 5 mm Allen key.
2. Adjust the width position and firmly tighten the two fixing screws I (3–5 Nm) in the desired position.

We recommend that foot plate adjustment is carried out by a qualified technician.



To ensure a good position of the feet, two types of straps can be provided; the heel strap J (serial) and the calf strap (optional) attached to the leg rest support. Both adjustable by hook and loop fasteners.

## 4.8 Amputee support



### WARNING!

#### Risk of tipping

There is an increased risk of tipping for users with amputated legs.

- Use anti-tip device and/or rebalance the wheelchair when amputee supports are used.




### WARNING!

#### Risk of falling

The wheelchair can tip forwards when getting into the wheelchair and applying load on the amputee support.


- Never sit or lean on the amputee support.

### Adjusting the angle / height on amputee supports

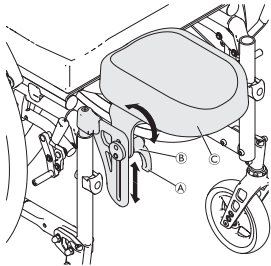

1. Loosen the lever (A) by turning it counter clockwise.
2. Loosen the nut (B) by turning it counter clockwise.
3. Adjust the desired angle and height of the amputee support (C).
4. Re-tighten the lever and nut.

### Removing / installing amputee support

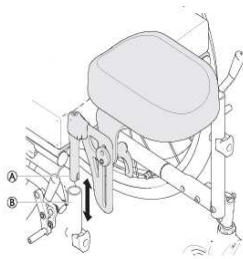

1. Remove the amputee support (A) by simply pulling it out of the frame tube (B).
2. Insert the amputee support by pushing it completely into the holder.

## 4.9 Dual handrim DHR

This quick release dual handrim system (left or right hand) allows the user to drive the wheelchair with a single arm using two handrims on the same wheel with 2 possible side positions of the small handrim.

The user catches both handrims to propel the wheelchair forward or backwards; turns on the LH or RH side by gripping the large or the small handrim only. The DHR system is always delivered in passive rear wheel position for a better stability.



The handrims can either be placed in line with each other (internal position) or with the inner, smaller, handrim outside the larger handrim (external position). Make sure that the 3 screws (B) are locked into position.



### CAUTION!

#### Risk of injury

- Its recommended that the change of control side and position of the small handrim must be carried out by a qualified technician.



### WARNING!

#### Risk of injury

- Always remove the telescopic rear wheel shaft (C) before removing the rear wheels.



### WARNING!

#### Risk of injury

Wheels can come off

- Check that the rear wheel is securely locked in position! It should not be possible to remove the wheels when the removable axle button (A) is inactivated. Test this by trying to pull the wheel off.



### CAUTION!

#### Risk of injury

Minor bruises, pain

- Be careful not to get your fingers trapped between the rear wheel spokes and the three brackets of the outer handrim.

### Disengaging and engaging the rear wheel shaft

To fold the wheelchair it is necessary to disengage the telescopic rear wheel shaft (C).

1. Disengage the telescopic shaft (C) by sliding one part inside the other.
2. To facilitate this operation move the small handrim forwards to backwards.
3. Reverse procedure to engage the telescopic rear wheel shaft.



### CAUTION!

#### Risk of accident

- The telescopic rear wheel shaft (C) is an integral part of the wheelchair and the user will be unable to propel the wheelchair without it.


### Removing and installing the rear wheel

The removable and interchangeable rear wheels facilitate transport as well as change of the control side.

1. Remove the telescopic rear wheel shaft (C).
2. Press the removable axle button (A). Keep it pressed and pull the wheel outwards.

3. Remove the rear wheel.
4. Repeat procedure on the opposite side.
5. Reverse procedure to reinstall the rear wheels.

## 4.10 Cane holder

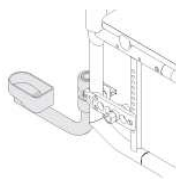

### WARNING!

#### Risk of overturning

- Ensure that the cane holder does not project beyond the external diameter of the rear wheel.



### NOTICE!


#### Risk of damage

- Never on the cane holder or use it as a tipping aid.

1. Place the cane in the holder.
2. Fix the top of the cane to the backrest

## 4.11 Calf strap / heel loop

A calf strap and/or heel loops are installed to the leg rests of the wheelchair to provide a better position of the users' legs.

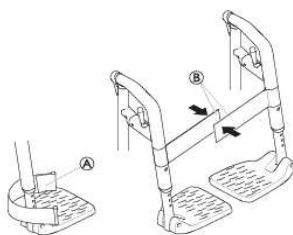

### Adjusting the heel loops

1. Adjust the heel loops (A) to the desired depth using the hook-and-loop fasteners.

### Installing / adjusting the calf strap

1. Completely unfold the wheelchair.
2. Close the calf strap by placing both hook-and-loop straps (B) on top of each other.

## 4.12 Passive illumination

You can attach reflectors to the rear wheels.

## 4.13 Tension adjustable seat upholstery

The tension of the straps must not disrupt the good operation of the folding mechanism, as excessive tension may prevent the locking of the mechanism. Your therapist must verify the adjustment.

# 5 Setup

## 5.1 Safety Information



### CAUTION!

#### Risk of injury

- Before using the wheelchair, check its general condition and its main functions in 8.2 *Maintenance schedule*, page 24

Your provider will supply your wheelchair ready for use. Your provider will explain the main functions and ensure that the wheelchair meets your needs and requirements.

Adjustments of the axle position and the castor supporters must be carried out by a qualified technician.

If you receive your wheelchair folded, read section 1 *Folding and unfolding the backrest*, page 1.

# 6 Using the wheelchair

## 6.1 Safety Information



### WARNING!

#### Risk of accidents

- Uneven tyre pressure can have a huge effect on handling.
- Check the tyre pressure before each journey.



### WARNING!

#### Risk of falling out of the wheelchair

- When using front wheels that are too small, the wheelchair could get stuck at curbs or in floor grooves.
- Make sure that the front wheels are suitable for the surface you're driving on.



### CAUTION!

#### Risk of crushing

- There may be a very small gap between the rear wheel and the mudguard or armpad with the risk that you could trap your fingers.
- Ensure that you always propel your wheelchair using the handrims only.



### CAUTION!

#### Risk of crushing

- There may be a very small gap between the rear wheel and the mudguard or armpad with the risk that you could trap your fingers.
- Ensure that you always propel your wheelchair using the handrims only.

### Theft and metal detection systems

In seldom cases the materials used in the wheelchair may activate theft and metal detection systems.

## 6.2 Braking during use

Whilst you are moving, you brake by transferring force to the handrim with your hands.

**WARNING!****Risk of accidents**

If you engage the parking brakes while you are moving, the direction of movement can become uncontrollable and the wheelchair may stop suddenly, which can lead to a collision or to you falling out.

- Never engage the parking brakes while you are moving.

**WARNING!****Risk of falling of the wheelchair**

If the wheelchair is rapidly decelerated by an assistant pulling at the push handles, the user may fall out of the wheelchair.

- Always apply the posture belt if present
- Make sure your assistant is qualified in transferring occupied wheelchairs.
- Operate the two brakes simultaneously and do not take slopes exceeding 7° to ensure perfect control of the wheelchair steering.
- Always use the parking brakes to secure the wheelchair during short or long stops.

**CAUTION!****Risk of burning your hands**


If you brake for a long time, a lot of frictional heat is produced at the handrims (especially anti slip handrims)

- Wear suitable gloves.

1. Hold the handrims and press evenly with both hands until the wheelchair stops.

**Braking with an attendant**

With an attendant's brake (drum brake) an attendant can brake the wheelchair while it is moving. The attendant's brake can also be used as a parking brake.

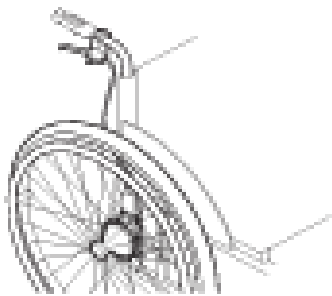

1. To brake, pull the brake lever on the push handle.
2. To park, pull the brake lever firmly and move the safety lever upwards.
3. To release, pull the brake lever until the safety lever disengages.

## 6.3 Getting in and out of the wheelchair

**WARNING!****Risk of overturning**

There is a high risk of overturning during the transfer.

- Only get in and out without assistance if you are physically able to do so.
- When transferring, position yourself as far back as possible in the seat. This will prevent damaged upholstery and the possibility of the wheelchair tipping forward.
- Make sure that both castors are facing straight to the front.

**WARNING!****Risk of overturning**

The wheelchair could tip forwards if you stand on the footrest.

- Never stand on the footrest when getting in and out.

**CAUTION!**

A heavy load hanging on the backrest can affect the wheelchair's centre of gravity.

- Change your driving style accordingly.

**CAUTION!**

If you disengage or damage the brakes the wheelchair could roll away out of control.

- Do not support yourself on the brakes when getting in and out.

**NOTICE!**

The mudguards and armrests could become damaged.

- Never sit on the mudguards or armrests when getting in and out.



1. Propel the wheelchair as near as possible to the seat that you want to move to.
2. Engage the parking brakes.
3. Remove the armrests or move them upwards out of the way.
4. Detach the leg rests or swing them outwards.
5. Place your feet on the ground.
6. Hold the wheelchair and, if necessary, also hold a fixed object in the vicinity.
7. Move slowly to chair.

## 6.4 Driving and steering the wheelchair

You drive and steer the wheelchair using the handrims. Before driving without an assistant you must find your wheelchair's tipping point.



**WARNING!****Risk of tipping**

The wheelchair can tip backwards.

- When finding the tipping point, an assistant must stand immediately behind the wheelchair to catch it if it tips over.
- To prevent tipping, install anti-tip device.

**WARNING!****Risk of tipping**

The wheelchair can tip forward.

- When setting up your wheelchair, test its behaviour in terms of tipping forward and adjust your driving style accordingly.

**CAUTION!**

A heavy load hanging on the backrest can affect the wheelchair's centre of gravity.

- Change your driving style accordingly.

**Finding the tipping point**

1. Release the brake.
2. Roll backwards a short distance, grasp both handrims firmly and push forwards with a slight kick.
3. The shift in weight and steering in opposite directions with the handrims will enable you to identify the tipping point.

**6.5 Negotiating steps and kerbs****WARNING!****Risk of overturning**

When negotiating steps you could lose your balance and tip the wheelchair over.

- Always approach steps and kerbs slowly and carefully.
- Do not go up or down steps that are higher than 25 cm.

**CAUTION!**

An activated anti-tip device prevents the wheelchair from tipping backwards.

- Deactivate the anti-tip device before going up or down steps or kerbs.

**WARNING!****Risk of injury to the assistant and damage to the wheelchair**

Tilting the chair with a heavy user can hurt the assistant's back and damage the chair.

- Make sure to be able to safely control the wheelchair with a heavy user before performing a tilting maneuver.

**Going down a step with an assistant**


1. Move the wheelchair right up to the kerb and hold the handrims.
2. The assistant should hold both push handles, place one foot on the tipper aid (if installed) and tilt the wheelchair backwards so that the front wheels lift off the ground.
3. The assistant should then hold the wheelchair in this position, push it carefully down the step and then tilt it forwards until the front wheels are back on the ground.

**Going up a step with an assistant****WARNING!****Risk of serious injuries**

Going up steps and kerbs often than the normal use can cause a fatigue break of the wheelchair backrest and the user might fall out of the wheelchair.

- Make sure that the wheelchair can be controlled by the user or the assistant before performing a tilting maneuver.

1. Move the wheelchair backwards until the rear wheels touch the kerb.
2. The assistant should tilt the wheelchair using both push handles so that the front wheels lift off the ground, then pull the rear wheels over the kerb until the front wheels can be placed back on the ground.

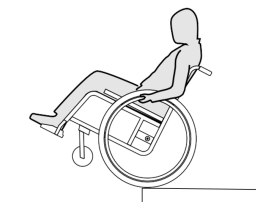**Going down a step without an assistant****WARNING!****Risk of tipping**

When going down a step without an assistant you could tip over if you cannot control your wheelchair.

- First learn how to go down a step with an assistant.
- Learn how finding the tipping point to balance on the rear wheels, *6.4 Driving and steering the wheelchair*, page 18

1. Move the wheelchair right to the kerb, lift the front wheels and keep the wheelchair balanced.
2. Now slowly roll both rear wheels over the kerb. While doing this, hold the handrims firmly with both hands until the front wheels are back on the ground.

## 6.6 Going up and down stairs



### WARNING!

#### Risk of overturning

When negotiating stairs you could lose your balance and overturn your wheelchair.

- Always use two assistants when negotiating stairs with more than one step.



1. You can go up and down stairs by taking them one step at a time, as described above. The first assistant should stand behind the wheelchair holding the push handles. The second assistant should hold a solid part of the front frame to steady the wheelchair from the front.

## 6.7 Negotiating tamps and slopes



### WARNING!

#### Risk due to wheelchair being out of control

When negotiating slopes or gradients your wheelchair could tip backwards, forwards or sideways.

- Always have an assistant behind the wheelchair when approaching long slopes.
- Avoid lateral slopes.
- Avoid slopes of more than 7°.
- Avoid jerking when changing direction on a slope.




### CAUTION!

Your wheelchair could run away even on slightly sloping ground if you do not control it using the handrims.

- Use the parking brakes if your wheelchair is stationary on sloping ground.

### Going up slopes


To go up a slope, you must create some momentum, keep up the momentum and control the direction at the same time.

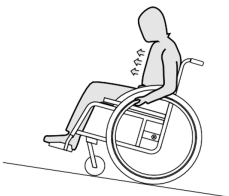

1. Bend your upper body forwards and propel the wheelchair with quick, powerful strokes on both handrims.

### Going down slopes

When going down slopes, it is important to control your direction and particularly your speed.

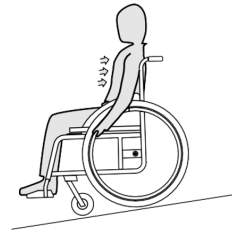

1. Lean back and carefully allow the handrims to run through your hands. You should be able to stop the wheelchair at any time by gripping the handrims.

## 6.8 Stability and balance when seated

Some everyday activities and actions require you to lean forwards, sideways or backwards out of the wheelchair. This has a major effect on the wheelchair's stability. To keep your balance at all times, proceed as follows:

### Leaning forwards




### WARNING!

#### Risk of falling

If you lean forwards out of the wheelchair you could fall out of it.

- Never bend too far forwards and do not shift forwards in your seat to reach an object.
- Do not bend forwards between your knees to pick up something off the floor

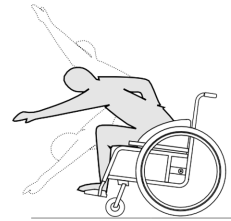

1. Point the front wheels forwards. (To do this, move your wheelchair forwards slightly then back again.)
2. Apply both parking brakes.
3. When leaning forwards your upper body must remain over the front wheels.

### Reaching backwards




### WARNING!

#### Risk of falling

If you lean too far backwards you could tip your wheelchair over.

- Do not lean out over the backrest.
- Use an anti-tip device.

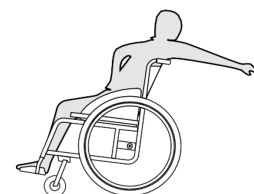

1. Point the front wheels forwards. (To do this, move your wheelchair forwards slightly then back again.)
2. Do not apply the parking brakes.
3. When reaching backwards do not reach so far that you have to change your sitting position.

## 6.9 Transporting the occupied wheelchair in a vehicle

Even when the wheelchair is properly secured and the following rules are met, injuries to passengers may occur in a collision or sudden stop. Therefore Invacare strongly recommends transferring the wheelchair user to the vehicle seat with the 3-point passenger restraint system on.

Do not make alterations or substitutions to points of the wheelchair (structure, frame or parts) without the written consent of Invacare. The wheelchair has been successfully tested according to the requirements of ISO 7176-19 (Frontal impact).



### **WARNING!**

#### **Risk of serious injury or death**

To use the wheelchair as a seat in a vehicle the backrest height minimum must be at least 400 mm.

To transport the wheelchair with user in a vehicle, a restraint system must be installed in the vehicle. Wheelchair tie-downs and occupant restraint systems must be approved according to ISO 10542-1.

Contact your local Invacare authorized provider for more information on getting and installing an approved and compatible restraint system.



### **WARNING!**

If, for some reason, it is impossible to transfer the wheelchair user to a vehicle seat, the wheelchair can be used as a seat in a vehicle if the following procedures and regulations are followed. A transportation kit (optional) have to be fitted on the chair for such a purpose.

- The wheelchair must be secured in the vehicle with a 4-point wheelchair restraint system.
- The user must wear a 3-point passenger restraint system secured to the vehicle.
- The user can be additionally secured with a posture belt in the wheelchair.



### **WARNING!**

Safety restraint devices must only be used when the wheelchair user's weight is 22 kg or more (ISO-7176-19).

- Do not use the wheelchair as a seat in a vehicle when the user weight is lower than 22 kg.



### **WARNING!**

- Before journey contact transporter and request information about the capability for the below required equipment.
- Make sure there is sufficient free space around the wheelchair and user to avoid the user making contact with other vehicle occupants, unpadded parts of the vehicle, wheelchair options or anchor points of the restraint system.



### **WARNING!**

- Make sure the tie-down points at the wheelchair are not damaged and that the parking brakes are fully functional.
- It's recommended to use puncture-proof tyres during transport to avoid brake problems due to reduced tyre pressure.



### **WARNING!**

Injury or damage may occur from wheelchair components or options loosened during a collision or sudden stop.

- Ensure all removable or detachable components and options are removed from the wheelchair and securely stored in the vehicle.
- It is essential to have your wheelchair checked by a qualified technician after an accident, collision etc.



We strongly recommend transferring the wheelchair user to the vehicle with the posture belt on.




### **NOTICE!**

- Refer to the user manuals supplied with the restraint systems.
- The following illustrations may differ depending on the restraint system supplier.



The choice of wheelchair configuration (seat width and depth, wheelbase) influences maneuverability and access to motor vehicles.

### **Securing the wheelchair with a 4-point restraint system**

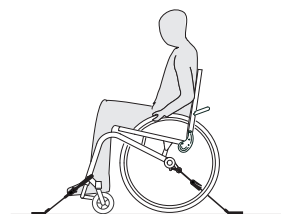


### **WARNING!**

- Place the wheelchair with the user forward-facing in direction of vehicle travel.
- Engage the wheelchair parking brakes.
- Activate the anti-tip device (if installed).

The wheelchair tie-down positions where the 4-point restraint system straps must be placed are marked with snap hook symbols, see following figures and section .


1. Using the front and rear straps of the 4-point restraint system, secure the wheelchair to the vehicle mounted rails. Refer to the user manual supplied with 4-point restraint system.

### Front side tie-down positions for belt straps:

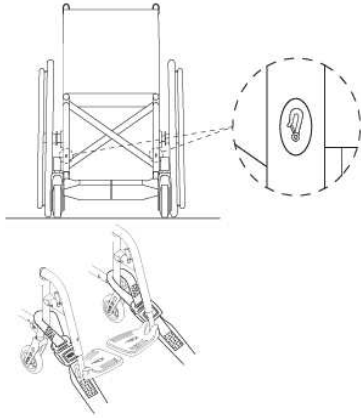

1. Attach the front straps above the castor supporters as shown in the two figures above (see location of the attachment labels).
2. Attach front straps to the rail system referring to best practice recommended instructions from the safety belt manufacturer.
3. Release the parking brakes and apply tension front straps by pulling the wheelchair backwards from the rear.
4. Readjust the parking brakes.

### Rear side tie-down positions for snap hooks:

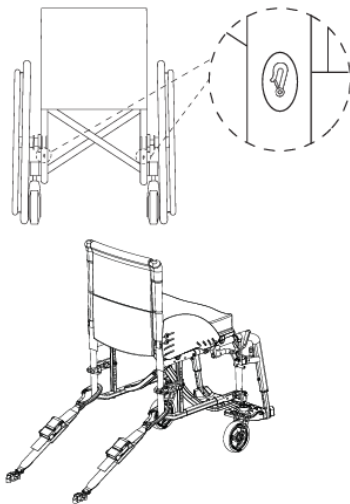

1. Attach rear straps to the rail system referring to best practice recommended instructions from the safety belt manufacturer.
2. Tighten the straps.



#### NOTICE!

- Make sure that the snap hooks are covered with slip resistant material to avoid laterally slipping on the axle.



#### NOTICE!

- Check the plungers are fully engaged on both sides and located in the same position of the cut-out section of the rail.
- Make sure that the angle between the rails and the straps lies between 40° and 45°.

### Fit posture belt




#### WARNING!

- The posture belt can be used in addition to but never as a substitute for an approved 3-point passenger restraint system.

1. Adjust posture belt to fit the wheelchair occupant, see section 4.2 *Posture belt*, page 12.

### Fastening the 3-point passenger restraint system

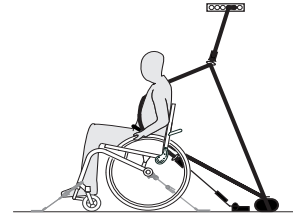

#### NOTICE!

- The illustration above may differ depending on the restraint system supplier.



#### WARNING!

- Ensure the 3-point passenger restraint system fits as tightly across the user's body as possible without discomfort and no part is twisted.
- Ensure the 3-point passenger restraint system is not held away from the user's body by parts of the wheelchair such as armrests or wheels etc.
- Ensure the 3-point passenger restraint system has a clear path from the user to the anchor point without interference by any part of the vehicle, wheelchair, seating or option.
- Ensure the pelvic belt fits snugly over the occupant's pelvis and is not allowed to ride up into the abdominal area.
- Ensure the user is able to reach the release mechanism unaided.




1. The 3-point belt restraint should fit as tightly across the user's body as possible without discomfort.



#### NOTICE!

- Ensure that the belt webbing is not twisted when in use.

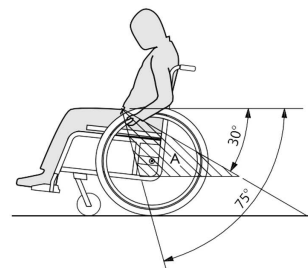

**! NOTICE!**

- Apply the pelvic belt of the 3-point passenger restraint system low across the pelvis so that the angle of the pelvic belt is within the preferred zone **A** of 30° to 75° to the horizontal. A steeper angle is preferred, but never exceeding 75°.



The image shows an incorrect placement of three-point belt restraint.

- The belt restraint must not be held away from the user's body by parts of the wheelchair such as armrests or wheels etc.
- If the vehicle is equipped with a head restraint and the wheelchair with the head/neck support, then assess with a clinician which one should be used.



The Invacare head and neck supports are not restraint devices. Nevertheless, they have been crash tested to verify their attachment points.

## 6.10 Transporting the manual wheelchair without occupant

**WARNING!**

**Risk of injury if the wheelchair is not properly secured**

In the event of an accident, braking manoeuvre, etc. you may suffer serious injuries from flying wheelchair parts.

- Firmly secure all wheelchair components in the means of transport to prevent them from coming loose during the journey.



You can remove the rear wheels when transporting the wheelchair.

**NOTICE!**

**Excessive wear and abrasion could affect the strength of load-bearing parts.**

- Do not pull your wheelchair across abrasive surfaces without the wheels fitted (e.g. pulling the frame over tarmac).

**CAUTION!**

**Risk of injury**

- If you are unable to fasten your mobility device securely in a transport vehicle, Invacare recommends that you do not transport it.

Your mobility device may be transported without restrictions, whether by road, rail or by air. Individual transport companies have, however, guidelines which can possibly restrict or forbid certain transport procedures. Please ask the transport company regarding each individual case.

- We strongly recommend securing the manual wheelchair to the floor of the transporting vehicle.

## 7 Transport

### 7.1 Safety information

**WARNING!**

**Risk of injury if the wheelchair is not properly secured**

In the event of an accident, braking maneuver, etc. you may suffer serious injuries from flying wheelchair parts.

- Always remove the rear wheels when transporting the wheelchair.
- Firmly secure all wheelchair components in the means of transport to prevent them from coming loose during the journey.
- When transporting the unoccupied wheelchair in a car or an airplane make sure that it is folded and secured.

**NOTICE!**

Excessive wear and abrasion could affect the strength of load-bearing parts.

- Do not pull your wheelchair across abrasive surfaces without the wheels fitted (e.g. pulling the frame over tarmac).

### 7.2 Folding and unfolding the wheelchair

**NOTICE!**

- Do not fold and unfold the wheelchair when the back is folded down.

#### Folding the wheelchair

**NOTICE!**

- Release the stabilizing bar, if present, before folding the wheelchair.

1. Remove the seat cushion, if present.
2. Fold the footrests upwards.
3. Pull the seat cover upwards.
4. Now you can fold the backrest down, if your wheelchair is fitted with a folding back.

#### Unfolding the wheelchair

**CAUTION!**

**Risk of getting fingers caught**

There is a risk of getting your fingers caught between the seat edge and the frame.

- Never clasp the seat edge with your fingers.

1. Position the wheelchair next to you.
2. If the back is down, fold it up.
3. Tilt the wheelchair, if possible, gently towards you.
4. With your hand flat, press the seat edge closest to you, so that the seat surface flattens out.
5. Put the wheelchair back completely on the floor and check that the seat edges on both sides are in the guide components.
6. Fold the footrests down.



## 7.3 Removing and installing the rear wheels




### WARNING!

#### Risk of overturning

If the removable axle of a rear wheel is not fully engaged, the wheel can become loose during use. This can lead to overturning.

- Always ensure that the removable axles are fully engaged whenever you fit a wheel.

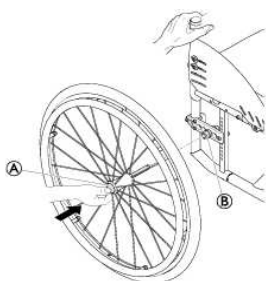

### Removing the rear wheels

1. Disengage the wheel locks.
2. With one hand, hold the wheelchair upright.
3. With the other, hold the wheel through the external spoke rim around the wheel hub.
4. Using your thumb, press the removable axle button (A). Keep it pressed and pull the wheel out of the adapter sleeve (B).

### Installing the rear wheels

1. Disengage the wheel locks.
2. With one hand, hold the wheelchair upright.
3. With the other, hold the wheel through the external spoke rim around the wheel hub.
4. Using your thumb, press the removable axle button (A). Keep it pressed and pull the wheel out of the adapter sleeve (B).
5. Push the axle into the adapter sleeve (B) up to the stop.
6. Release the removable axle button and make sure that the wheel is secure.

## 8 Maintenance

### 8.1 Safety Information



### WARNING!

#### Some materials deteriorate naturally over time.

This could result in damage to wheelchair components.

- Your wheelchair should be checked by a qualified technician at least once a year or if it has not been used for a long period.



### Packaging for return

The wheelchair shall be sent to the authorized provider in an appropriate packaging to avoid damage during the shipment.

### 8.2 Maintenance schedule

To ensure safe and reliable operation, carry out the following visual checks and maintenance regularly or have it carried out by another person

|                                                   | Weekly | Monthly | Annually |
|---------------------------------------------------|--------|---------|----------|
| Check the tyre pressure                           | x      |         |          |
| Check that rear wheels are seated correctly       | x      |         |          |
| Check the adjustment of the backrest joint        | x      |         |          |
| Check posture belt                                | x      |         |          |
| Visual Check                                      |        | x       |          |
| Check the folding mechanism                       |        | x       |          |
| Check the castors and its fixation                |        | x       |          |
| Check bolts                                       |        | x       |          |
| Check spokes                                      |        | x       |          |
| Check parking brakes (mechanism and cable)        |        | x       |          |
| Have wheelchair checked by a qualified technician |        |         | x        |

### Check the tyre pressure

1. Check the tyre pressure, see section 11.3 Tyres, page 30.
2. Inflate the tyres to the required pressure.
3. Check the tyre tread at the same time.
4. If necessary, change the tyres.

### Check that rear wheels are seated correctly

1. Pull on the rear wheel to check that the removable axle is seated correctly. The wheel should not come off.
2. If the rear wheels are not engaged properly, remove any dirt or deposits. If the problem persists, have the removable axles re-fitted by a qualified technician.

### Check the adjustment of the backrest joint

1. Sit into the wheelchair and lean back on the backrest. The backrest must be correctly engaged.
2. Pull the release cord. The pins must be free-moving on both sides. Otherwise, the backrest joint will have to be readjusted by your provider.

### Check posture belt

1. Check that the posture belt is adjusted correctly.



### NOTICE!


- Loose posture belts must be adjusted by an authorized provider.
- Damaged posture belts must be replaced by a qualified technician.

### General visual check

1. Examine your wheelchair for loose parts, cracks or other defects.
2. If you find anything, stop using your wheelchair immediately and contact an authorized provider.

### Check the castors and its fixation

1. Check that the castors turn freely.
2. Push and pull in the 3 directions (front and rear; left and right; up and down) to observe that there is no free play and movement, also check that there is no visual damage.

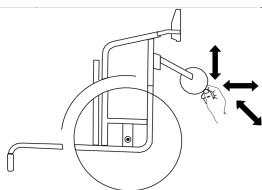

3. Remove any dirt or hair from the castor bearings.
4. Faulty or worn out castors fixation must be replaced by a qualified technician.

### Check the bolts

Bolts can work loose through constant use.

1. Check that the bolts are tight (on the footrest, seat cover, sides, backrest, frame, seat module).
2. Tighten any loose bolts with the suitable torque. Therefore refer to the service manual, available on our web site.



#### NOTICE!

Self-locking screws/nuts or thread-locking adhesive are used for several connections. If these are loosened, they must be replaced by new self-locking screws/nuts or they must be secured using new thread-locking adhesive.

- Self-locking screws/nuts must be replaced by a qualified technician

### Check the spoke tension

1. The spokes should not be loose or distorted.
2. Broken spokes must be replaced by a qualified technician.
3. Loose spokes must be tightened by a qualified technician.

### Check the parking brakes

1. Check that the parking brakes are positioned correctly. The brake is set correctly if the brake shoe depresses the tyre by a few millimetres when the brake is engaged.
2. If you find that the setting is not correct, have the brakes correctly adjusted by a qualified technician



#### NOTICE!

The parking brakes must be reset after replacing the rear wheels or changing their position.

### Checking after a heavy collision or blow



#### NOTICE!

The wheelchair can sustain visibly undetectable damage as a result of a heavy collision or hard blow.

- It is essential to have your wheelchair checked by a qualified technician after a heavy collision or hard blow.

### Repairing or changing an inner tube

1. Remove the rear wheel and release any air from the inner tube.
2. Lift one tyre wall away from the rim using a bicycle tyre lever. Do not use sharp objects such as a screwdriver which could damage the inner tube.
3. Pull the inner tube out of the tyre.
4. Repair the inner tube using a bicycle repair kit or, if necessary, replace the tube.
5. Inflate the tube slightly until it becomes round.
6. Insert the valve into the valve hole on the rim and place the tube inside the tyre (the tube should lie right round the tyre with no creases).
7. Lift the tyre wall over the edge of the rim. Start close to the valve and use a bicycle tyre lever. When doing this, check all the way round to ensure that the inner tube is not trapped between the tyre and the rim.
8. Inflate the tyre to the maximum operating pressure. Check that no air is escaping from the tyre.



#### Spare parts

All spare parts may be obtained from an Invacare authorized provider.

## 8.3 Cleaning and disinfection

### 8.3.1 General safety information



#### CAUTION!

##### Risk of Contamination

Take precautions for yourself and use appropriate protective equipment.



#### NOTICE!

Wrong fluids or methods can harm or damage the product.

All cleaning agents and disinfectants used must be effective, compatible with one another and must protect the materials they are used to clean.

- Never use corrosive fluids (alkalines, acid etc.) or abrasive cleaning agents. We recommend an ordinary household cleaning agent such as dishwashing liquid, if not specified otherwise in the cleaning instructions.
- Never use a solvent (cellulose thinner, acetone etc.) that changes the structure of the plastic or dissolves the attached labels.
- Always make sure that the product is completely dried before taking it into use again.



For cleaning and disinfection in clinical or long-term care environments, follow your in-house procedures.

### 8.3.2 Cleaning intervals



#### NOTICE!

Regular cleaning and disinfection enhance smooth operation, increases the service life and prevents contamination.

- Clean and disinfect the product:
- Regularly while in use,
- Before and after any service procedure,
- When it has been in contact with any body fluids,
- Before using it for a new user.

### 8.3.3 Cleaning


**NOTICE!**

- The product does not tolerate cleaning in automatic washing plants, with high-pressure cleaning equipment or steam.


**NOTICE!**

- Dirt, sand and seawater can damage the bearings and steel parts can rust if the surface is damaged.
- Only expose the wheelchair to sand and sea water for short periods and clean it after every trip to the beach.
  - If the wheelchair is dirty, wipe off the dirt as soon as possible with a damp cloth and dry it carefully.

1. Remove any options fitted (only options which do not require tools).
2. Wipe down the individual parts using a cloth or soft brush, ordinary household cleaning agents (pH = 6 - 8) and warm water.
3. Rinse the parts with warm water
4. Thoroughly dry the parts with a dry cloth.



Car polish and soft wax can be used on painted metal surfaces to remove abrasions and restore gloss.

#### Cleaning upholstery

For cleaning upholstery refer to the instructions on the labels of the seat, cushion and backrest cover.



If possible, always overlap hook and loop strips (the self-gripping parts) when washing, to minimize lint and thread build-up on hook strips and prevent damage to upholstery fabric by these.

### 8.3.4 Disinfection

The wheelchair may be disinfected by spraying or wiping with tested, approved disinfectants.



Spray a soft cleaning and disinfecting product (bactericidal and fungicide meeting the EN1040/ EN1276 / EN1650 standards) and follow the instructions given by the manufacturer.

1. Wipe down all generally accessible surfaces with a soft cloth and ordinary household disinfectant.
2. Allow the product to air-dry.

## 9 Troubleshooting

### 9.1 Safety Information

Faults may arise as a result of daily use, adjustments or changing demands on the wheelchair. The table below shows how to identify and repair faults. Some of the actions listed must be carried out by a qualified technician. These are indicated. We recommend that all adjustments are carried out by a qualified technician.


**CAUTION!**

– If you notice a fault with your wheelchair, e.g. a significant change in handling, stop using your wheelchair immediately and contact your provider.

### 9.2 Identifying and repairing faults

| Fault                                             | Possible cause                            | Action                                                                                      |
|---------------------------------------------------|-------------------------------------------|---------------------------------------------------------------------------------------------|
| The wheelchair does not travel in a straight line | Incorrect tyre pressure on one rear wheel | Correct tyre pressure → 11.3 Tyres, page 30                                                 |
|                                                   | One or more spokes broken                 | Replace faulty spoke(s) → qualified technician                                              |
|                                                   | Spokes tightened unevenly                 | Tighten loose spokes → qualified technician                                                 |
|                                                   | Castor supporter is not aligned correctly | Correct position and alignment of the castor supporter → qualified technician               |
|                                                   | Castor bearings are dirty or damaged      | Clean the bearings or replace the castor → qualified technician                             |
|                                                   | Rear wheels are off track                 | Change track → qualified technician                                                         |
|                                                   | Castors not adjusted to same height       | Position the castors such that they touch the floor at the same time → qualified technician |
| The wheelchair tips backwards too easily          | Rear wheels are mounted too far forwards  | Mount the rear wheels further back → qualified technician                                   |
|                                                   | Back angle too large                      | Reduce the backrest angle → qualified technician                                            |
| Removable axles cannot be inserted correctly      | Removable axles dirty                     | Clean removable axles                                                                       |
|                                                   | Removable axles misaligned                | Adjust removable axles → qualified technician                                               |

| Fault                                                    | Possible cause                                    | Action                                                                                |
|----------------------------------------------------------|---------------------------------------------------|---------------------------------------------------------------------------------------|
| The parking brakes are gripping poorly or asymmetrically | Incorrect tyre pressure in one or both rear tyres | Correct tyre pressure → 11.3 Tyres, page 30                                           |
|                                                          | Brake setting incorrect                           | Correct the brake setting → qualified technician                                      |
| Rolling resistance is very high                          | Tyre pressure in rear tyres is too low            | Correct tyre pressure → 11.3 Tyres, page 30                                           |
|                                                          | Tyre is punctured                                 |                                                                                       |
|                                                          | Rear wheels not parallel                          | Ensure the rear wheels are parallel → qualified technician                            |
| The castors wobble when moving fast                      | Too little tension in castor bearing block        | Tighten the nut on the bearing block axle slightly → qualified technician             |
|                                                          | Castor has worn smooth                            | Change castor → qualified technician                                                  |
| The castor is stiff or stuck                             | Bearings are dirty or faulty                      | Clean or replace the bearings → qualified technician                                  |
| The backrest is very difficult to unfold                 | Backrest cover is too tight                       | Loosen the hook-and-loop straps of the backrest cover a little → 3.5 Backrest, page 8 |

## 10 After use

### 10.1 Storage


**NOTICE!**
**Risk of damage to the product**

Do not store the product near heat sources.

Never store other items on top of the wheelchair.

Store the wheelchair indoors in a dry environment.

Refer to temperature limitation in section 11.5 *Environmental conditions*, page 30.

After long-term storage (more than four months) the wheelchair must be inspected in accordance to section 8.2 *Maintenance schedule*, page 24.

### 10.2 Reconditioning

This product is suitable for reuse. To recondition the product for a new user, carry out the following actions:

- Inspection
- Cleaning and disinfection
- Adaptation to the new user

For detailed information, please see section 8.2 *Maintenance schedule*, page 24 and the service manual for this product. Make sure that the user manual is handed over with the product.

If any damage or malfunction is detected, do not reuse the product.

### 10.3 Disposal

Be environmentally responsible and recycle this product through your recycling facility at its end of life.

Disassemble the product and its components, so the different materials can be separated and recycled individually.

The disposal and recycling of used products and packaging must comply with the laws and regulations for waste handling in each country. Contact your local waste management company for information.


## 11 Technical data

### 11.1 Dimensions and weight


All dimension and weight specifications refer to a wide range of the wheelchair in a standard configuration. Dimension and weight (based on ISO 7176–1/5/7) may alter according to different configurations. In some configurations, the wheelchair exceeds the recommended dimensions for a wheelchair.


**NOTICE!**

- In some configurations, the overall dimensions of the wheelchair when it is ready for use exceeds the authorised limits and the access to emergency escape routes is not possible.
- In some configurations, the wheelchair exceeds the size recommended for travelling by train in the EU.

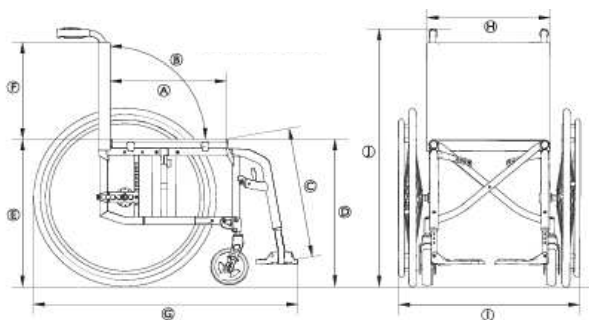

**Compact SA**

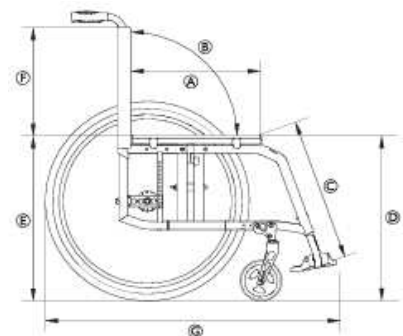

**Compact FF**



|   |                                     |                                                                                  |
|---|-------------------------------------|----------------------------------------------------------------------------------|
|   | Maximum load of the wheelchair (*)  | 130 kg                                                                           |
| Ⓐ | Seat depth                          | 350 – 500 mm, in increments of 20 mm                                             |
| Ⓑ | Backrest angle                      | 82°/86°/90°/94°/98°/102°                                                         |
| Ⓒ | Lower leg length                    | 200 – 510 mm, in increments of 10 mm                                             |
| Ⓓ | Front seat-to-floor height          | 370 – 530 mm, in increments of 10 mm                                             |
| Ⓔ | Rear seat-to-floor height           | 370 – 500 mm, in increments of 10 mm                                             |
| Ⓕ | Backrest height                     | 300 – 510 mm, in increments of 15 mm                                             |
| Ⓖ | Overall length with legrest         | 70° frame angle: approx. 870 – 1220 mm<br>80° frame angle: approx. 850 – 1200 mm |
| Ⓗ | Seat width                          | 280 – 500 mm, in increments of 20 mm                                             |
| Ⓘ | Overall width                       | Seat width + 180 mm                                                              |
|   | Overall width, folded               | approx. 290 mm                                                                   |
| ⓫ | Overall height                      | approx. 660 – 1060 mm                                                            |
|   | Folded height                       | approx. 660 – 1060 mm                                                            |
|   | Total weight                        |                                                                                  |
|   | Compact SA                          | approx. 7.8 kg                                                                   |
|   | Compact FF                          | approx. 6.7 kg                                                                   |
|   | Pivot width                         | 1420 mm                                                                          |
|   | Minimum turning radius              | 950 mm                                                                           |
|   | Max. safe slope                     | 7°                                                                               |
|   | Static stability downhill           | 20°                                                                              |
|   | Static stability uphill             | 7°                                                                               |
|   | Static stability sideways           | 20°                                                                              |
|   | Seat plane angle                    | 0° – 12°                                                                         |
|   | Leg to seat surface angle           | 100° / 110°                                                                      |
|   | Armrest to seat distance            | 120 – 280 mm                                                                     |
|   | Front location of armrest structure | 220 – 330 mm                                                                     |
|   | Handrim diameter                    | 450 – 590 mm                                                                     |
|   | Horizontal location of axle         | aprox. –60 to +90 mm                                                             |

(\*) The maximum load of the wheelchair is the total of the maximum occupant mass and the maximum mass of any other items intended to be carried by the wheelchair.

The maximum occupant mass is the maximum load of the wheelchair without additional load.

**Example:** for a wheelchair with a back pack of 10kg:

Maximum Occupant Mass = Maximum Load of the wheelchair = 10kg

## 11.2 Maximum weight of removable parts

The components used to manufacture Küschall wheelchairs consist of following materials:

| Maximum weight or removable parts                          |                |
|------------------------------------------------------------|----------------|
| Part                                                       | Maximum weight |
| Leg rest angle adjustable with calf pad and footrest       | 3.2 kg         |
| Armrest                                                    | 2 kg           |
| Rear wheel 600 mm (24") solid with handrim and spoke guard | 2.5 kg         |
| Headrest / Neckrest / Cheek support                        | 1.4 kg         |
| Trunk support                                              | 0.8 kg         |
| Backrest cushion                                           | 2 kg           |
| Seat cushion                                               | 2 kg           |
| Table tray                                                 | 3.9 kg         |

## 11.3 Tyres

The ideal tyre pressure depends on the tyre type.

The table below is an indication. In case the tyre differs from the list below, check the side of the tyre, the maximum pressure is often listed there.

| Tyre         | Diameter                                       | Max. pressure |          |         |
|--------------|------------------------------------------------|---------------|----------|---------|
| Profile tyre | 610 mm (24")                                   | 7 bar         | 700 kPa  | 101 psi |
| Other tyres  | 610 mm (24");<br>635 mm (25");<br>660 mm (26") | 10 bar        | 1000 kPa | 145 psi |

 The compatibility of the tyres listed above depends on the configuration and/or model of your wheelchair.

 In case of a tyre puncture consult a suitable workshop (e.g. bike repair shop, bicycle dealer ...) to have the tube replaced by a skilled person.

 The size of the tyre is mentioned on the sidewall of the tyre. The change of appropriate tyres must be carried out by a qualified technician.



### CAUTION!

— The tyre pressure has to be equal in both wheels to avoid a less driving comfort, to keep the brakes efficiency and an easy propelling of the wheelchair.

## 11.4 Materials

The components used to manufacture Küschall wheelchairs consist of following materials:

| Component                      | Material                |
|--------------------------------|-------------------------|
| Frame tubes                    | Aluminum                |
| Backrest tubes                 | Aluminum                |
| Axle tube                      | Aluminum                |
| Seat cover / Backrest cover    | PA / PE / PVC           |
| Push handles                   | Aluminum / TPE          |
| Clothes guard / Mudguard       | Carbon Fiber or Plastic |
| Castor forks                   | Aluminum                |
| Footrest tube                  | Aluminum                |
| Footplate                      | Carbon Fiber or Plastic |
| Supporting parts / Attachments | Steel / Aluminum        |
| Screws, washers and nuts       | Steel                   |

 All materials used are protected against corrosion. We use only REACH compliant materials and components.

 Theft and metal detection systems: in seldom cases the materials used in the wheelchair may activate theft and metal detection systems.

## 11.5 Environmental conditions

|                      | Storage and transportation            | Operation |
|----------------------|---------------------------------------|-----------|
| Temperature          | -20 °C to 40 °C                       |           |
| Relative humidity    | 20 % to 90 % at 30 °C, not condensing |           |
| Atmospheric pressure | 800 hPa to 1060 hPa                   |           |

 Be aware that when a wheelchair has been stored under low temperatures, it must be adjusted before use. Please see 8 *Maintenance*, page 24.





**United Kingdom & Ireland:**

Invacare Ltd.

Pencoed Technology Park, Pencoed

Bridgend CF35 5AQ

Tel: (44) (0) 1656 776 200

uk@invacare.com

www.invacare.co.uk



Invacare UK Operations Limited  
Unit 4, Pencoed Technology Park,  
Pencoed  
Bridgend CF35 5AQ  
United Kingdom



Invacare France Operations SAS  
Route de St Roch  
F-37230 Fondettes  
France

