
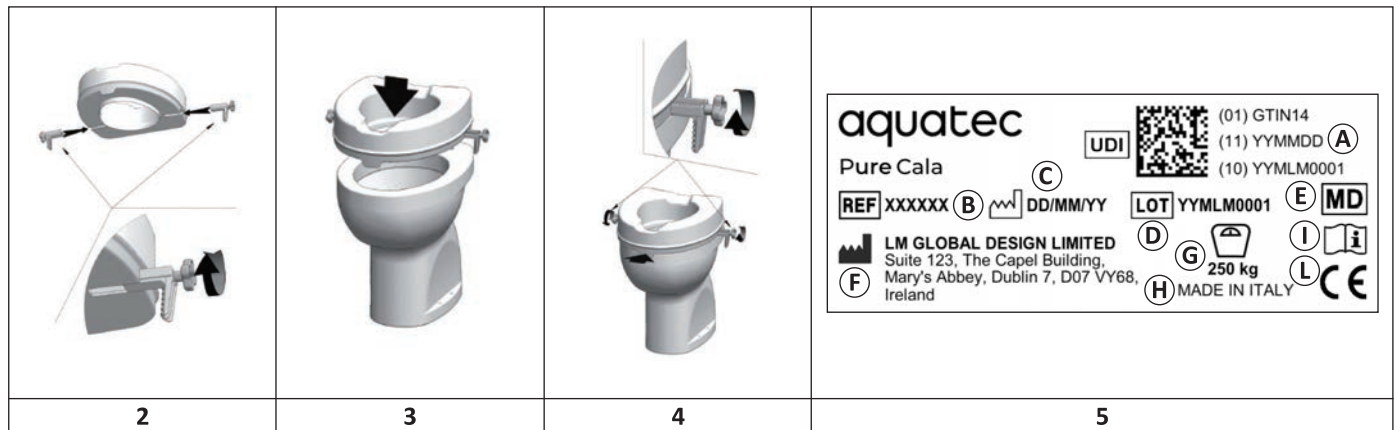
Aquatec Pure Cala

en	Raised Toilet Seat User Manual	1
de	Toilettensitzerhöhung Gebrauchsanweisung	3
fr	Réhausse WC Manual d'utilisation	5
es	Asiento elevador de indoor Manual del usuario	7
pt	Elevador de sanita Manual de utilizacao	9
nl	Verhoogde toiletzitting Gebruiksaanwijzing	11
it	Rialzo WC Manuale 'uso	13

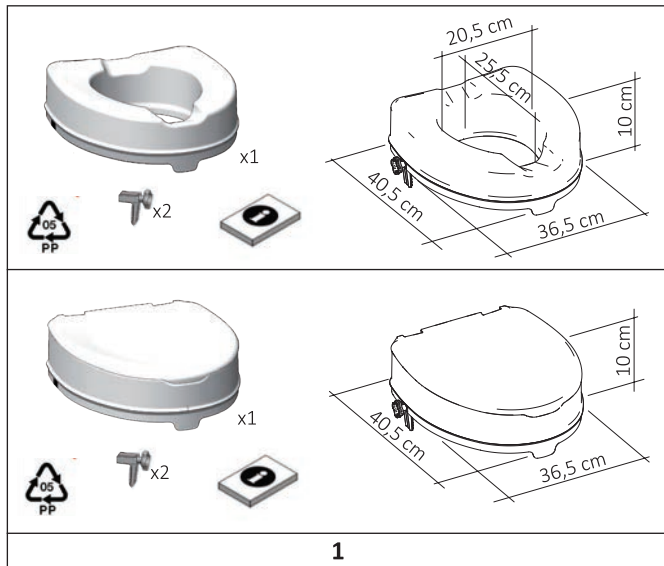1 General

1.1 Introduction

This user manual contains important information about the handling of the product. To ensure safety when using the product, read the user manual carefully and follow the safety instructions. If you find that the font size in the printed document is difficult to read, you can download the PDF version from the website.

Invacare reserves the right to alter product specifications without further notice.

In case of a serious incident with the product, you should inform the manufacturer and the competent authority in your country.

1.2 Warranty Information

We provide a manufacturer's warranty for the product in accordance with our General Terms and Conditions of Business in the respective countries.

Warranty claims can only be made through the provider from whom the product was obtained.

1.3 Compliance

This CE-marked medical device conforms to: EU Medical Device Regulation, ISO 13485, ISO 14971, REACH, ISO 17966 and related standards.

1.4 Service Life

The expected service life of this product is five years when used in accordance with safety instructions and intended use, stated in this manual. The effective service life can vary depending on frequency and intensity of use.

This manual **MUST** be given to the user of the product. **BEFORE** using this product, this manual **MUST** be read and saved for future reference.

aquatec

1.5 Limitation of Liability

Invacare accepts no liability for damage arising from: Non-compliance with the user manual, Incorrect use, Natural wear and tear, Incorrect assembly or set-up, Unauthorised modifications and/or use of unsuitable spare parts.

2 Safety

2.1 Safety Information

WARNING!

Risk of Serious Injury or Damage

Improper use of this product may cause injury or damage.

- If you are unable to understand the warnings, cautions or instructions, contact a health care professional or provider before attempting to use this equipment.
- Do not use this product or any available optional equipment without first completely reading and understanding these instructions and any additional instructional material.

WARNING!

Risk of Injury

- The maximum user weight indicated in the fitting instructions must not be exceeded.
- Before every use, be sure the product is in a stable position on the WC and correctly fixed on it by the locking system supplied.
- During the assembly, strictly adhere to what is shown in the user manual.
- Do not overtighten the locking system.
- Do not modify the product and its components in any way.
- Do not use the product if damaged or altered in any way.
- Be sure the product height is such as making it easy and safe for the user
- During use, avoid positions that keep the body weight too much forward. Wrong positions may lead to instability and overturning the product
- The product must not be stored outdoors.
- The risk of finger injuries may increase during the assembly process.
- Do not use the product for different purposes from the intended one.
- Avoid direct contact with open wounds or broken skin to prevent infection.

WARNING!

- Before first use, clean the product thoroughly.
- Before every use, check the product for damage and ensure that it is correctly positioned.

CAUTION!

Risk of Injury

- Non-original or wrong parts may affect the function and safety of the product.
- Only use original parts for the product in use.

2.2 Labels and Symbols on the Product (Fig. 5)

(A)	Unique Device Identifier	(F)	Manufacturer
(B)	Reference number	(G)	Maximum user weight
(C)	Production date	(H)	Made in Italy
(D)	Lot number	(I)	Read the user manual
(E)	Medical Device	(L)	European conformity

The identification label is attached to the back of the toilet seat.

3 Product Overview

Intended Use

The Aquatec Pure Cala is a raised toilet seat that elevates the seating height on toilets for 10 cm for a better seating position. It is fixed with a clamp system at the rear of the toilet bowl.

The product is intended to be used as an aid for self-hygiene care and toilet purposes. The product provides support for sitting down or standing up from a toilet.

Intended Users

Users are any adults and adolescents with limited mobility or limited muscle strength. The user must be able to keep own balance or ensures the availability of additional support i.e. a wall mounted grab bar.

The maximum user weight is 250kg.

Indications

- Limited mobility (e.g. legs and/or hips) but adequate body stability or
- having muscle strength issues to stand up/sit down due to transitional or permanent disabilities.

Contraindications

The use of the product is not recommended in cases of perception disorders or sitting balance impairments.

Material

Product is made from polypropylene. All product components are corrosion resistant and do not contain natural rubber latex.

4 Setup

1. Slide the fixing elements into the allocated cavity on the corpus. Twist the star screw just enough until the fixing element takes hold. (Fig. 2)
2. Flip up your toilet seat. (Alternative: Remove toilet seat)
3. Place the toilet seat raiser on the toilet bowl. (Fig. 3)
4. Push backward until the stop lug aligns with the toilet bowl. (Fig. 4)
5. Tighten the star screws. (Fig. 4)

5 Maintenance

5.1 General Maintenance Information

This product is maintenance-free, provided that the cleaning and disinfection instructions are adhered to.

5.2 Cleaning and Disinfection

General Safety Information

CAUTION!

Risk of Contamination

- Take precautions for yourself and use appropriate protective equipment.

NOTICE!

- Wrong fluids or methods can harm or damage the product.
- All cleaning agents and disinfectants used must be effective, compatible with one another and must protect the materials they are used to clean.
 - Never use corrosive fluids (alkalines, acid etc.) or abrasive cleaning agents. We recommend an ordinary household cleaning agent such as dishwashing liquid.
 - Never use a solvent (cellulose thinner, acetone etc.) that changes the structure of the plastic or dissolves the attached labels.

For cleaning and disinfection in clinical or long-term care environments, follow your in-house procedures.

Cleaning Intervals

NOTICE!

- Regular cleaning and disinfection enhances smooth operation, increases the service life and prevents contamination.
- Clean and disinfect the product:
- regularly while it is in use,
 - before and after any service procedure
 - when it has been in contact with any body fluids.

Cleaning Instructions

NOTICE!

- Maximum cleaning temperature:
- 60 °C by hand
 - 85 °C for 3 minutes by automatic washing systems

1. Clean the product by wiping down all generally accessible surfaces.
2. Dry the parts with a cloth.

Disinfection Instructions

Method: Follow the application notes for the used disinfectant and wipe-disinfect all accessible surfaces.
Disinfectant: Ordinary household disinfectant.
Drying: Allow the product to air-dry.

6 After Use

6.1 Reconditioning

This product is suitable for reuse. To recondition the product for a new user, carry out the following actions:

- Inspection
- Cleaning and disinfection

For detailed information, see 5 Maintenance, page 2.

Make sure that the user manual is handed over with the product.

6.2 Disposal

Disassemble the product and its components, so the different materials can be separated and recycled individually.

The disposal and recycling of used products and packaging must comply with the laws and regulations for waste handling in each country. Contact your local waste management company for information.

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Invacare Austria GmbH
Invacare Poirier SAS
LM Global Design Limited