

Invacare® | **Comet® , Comet® Alpine and Comet® HD**



Invacare **Comet**



Invacare **Comet Alpine**



Invacare **Comet HD**

Rugged, Robust, Reliable: The *Invacare Comet* Family

The **Comet** Family is packed with features that enable a safe, fast and enjoyable drive. Reliability remains the hallmark of this secure and speedy scooter. Whether stocking up at the supermarket or challenging its outdoor performance, you can rely on the **Comet** to get you where you want to go, quickly, in comfort and with peace of mind. By combining rugged outdoor capabilities with ergonomics and good looks, the *Invacare Comet* Family is the solution for your transportation needs.

Trusted Solutions, Passionate People



Yes, you can.®



Features

For a safe, controlled drive

With its large wheels, a powerful motor and adjustable suspension system in both the front and rear of the chassis, the **Comet** simply invites adventure. The **Comet** is equipped with an electronic speed reduction feature, which automatically decreases speed when going round bends: this ensures safe, comfortable manoeuvres at all times. The lighting system includes an easily visible DUAL-bulb headlight, indicators and integrated heavy-duty brake light for safety.

Compact for convenience

Ergonomic handles help you execute precise manoeuvres in tight spaces and offer relief to hands and arms during longer trips. The adjustable tiller lever is positioned opposite the handbrake and is within easy reach. A two-step disengaging lever protects the scooter from free-wheeling if accidentally knocked out of drive mode.



Finite tiller movement

Users can easily adjust the tiller with a lever to suit their needs.



Hand brake

Ensures immediate braking if required.



Brake light

Heavy duty brake lighting warns people that the scooter is slowing down.

Technical data

3 wheels	510 mm	480 mm	440 - 515 mm	500 mm	660 mm	1450 mm	133 kg
4 wheels	510 mm	480 mm	440 - 515 mm	500 mm	660 mm	1475 mm	136 kg
3 wheels	81 kg	160 kg	2 x 75 Ah	3130 mm	100 mm	135 mm	10°
4 wheels	84 kg	160 kg	2 x 75 Ah	3260 mm	100 mm	135 mm	10°
				* 10 km/h ** 12,8 km/h / 15 km/h *** Peak performance			
3 wheels	60* / 55** km	10 / 12,8 / 15 km/h	10 km/h : 250 W / 1200 W***				
4 wheels	60* / 55** km	10 / 12,8 / 15 km/h	12,8 km/h : 250 W / 1200 W***				
			15 km/h : 250 W / 1200 W***				

Colour of shrouds



Moonstone Silver



Onyx Blue



Ruby red

Upholstery



Black vinyl

Speed reduction functionality



Features and Accessories

Extra power

Beyond its stylish appearance, the **Comet[®] Alpine** proves itself a real winner. The scooter will win you over with its powerful driving characteristics and impressive manoeuvrability going uphill and even when driving downhill or on loose gravel. The stability of the **Comet[®] Alpine** in challenging surroundings leaves no doubt that this scooter is out of the ordinary.

Challenge your Comet[®] Alpine

With the **Comet[®] Alpine** you can leave limitations behind, enjoying outings across uneven or wet ground. The **Comet[®] Alpine** has built-in anti-splash guards, which protect vital components from water and dirt. If you want to explore unfamiliar terrain, just keep the optional "automatic speed-reduction" active: this will prevent you from going too fast when turning.



Seat suspension module

Offer yourself increased comfort by adding a seat suspension.



Seat raiser


















Electronically raises the seat to the preferred height.



Off-board charging

An off-board charger kit is available, which allows you to conveniently charge the batteries away from the scooter.

Technical data

Comet Alpine 4 wheels		510 mm		480 mm		440 - 515 mm		500 mm		660 mm		1475 mm		143 kg
Comet Alpine 4 wheels		91 kg		160 kg		2 x 75 Ah		3260 mm		100 mm		135 mm		10°
Comet Alpine 4 wheels		50* km		10 km/h		650 W / 1800 W**	<p>* 10 km/h ** Peak performance</p>							

Colour of shrouds



Slate Grey

Upholstery



Black vinyl



Features

Durable design and power in one

Even for users with a weight of up to 220 kg, safety is never a concern with the **Comet^{HD}**. All necessary precautions for guaranteeing the stability and safety of the scooter have been taken into consideration during the design process.


That comfortable feeling

The XL contoured seat has been developed in conjunction with therapists to ensure optimum seating comfort. The seat base and height- and width-adjustable armrests have been reinforced to allow for additional stresses and can easily be used as supports when getting on and off the scooter. Seating comfort can be increased by using the pneumatically adjustable lumbar support.



Unique suspension system

It has proven itself in the professional mountain biking industry and allows for simple independent user adjustment.

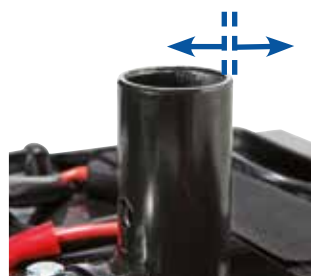

Reinforced seat column

It ensures stability and safety while seated.



Reinforced armrest

It allows for the additional stress of getting on and off the **Comet^{HD}**.

Technical data

Comet^{HD} 4 wheels	 660 mm	 530 mm	 440 - 515 mm	 610 mm	 660 mm	 1475 mm	 148 kg
Comet^{HD} 4 wheels	 96 kg	 220 kg	 2 x 75 Ah	 3260 mm	 100 mm		
Comet^{HD} 4 wheels	 135 mm	 10°	 50* km	 10 km/h	 650 W / 1800 W**		

* 10 km/h

** Peak performance

Invacare International Sàrl
Route de Cité-Ouest 2 - 1196 Gland
Switzerland
Tel: +41 (0)22 354 60 10
Fax: +41 (0)22 354 60 11
hqeurope@invacare.com
www.invacare.eu.com

© 2015 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

Colour of shrouds

Upholstery



Slate Grey



Black vinyl



Weight
capacity 220 kg



Yes, you can.®

Comet Family - EU - 08/2015