



Azalea[®] Max



A tailor-made solution for plus
size users



The **Rea Azalea Max** is designed to meet the specific needs of plus size, disabled users. The **Rea Azalea Max** is a niche product and compliments the existing **Azalea** family from Invacare.

Compatibility

The **Rea Azalea Max** includes features such as tilt-in-space with weight-shifting and DSS (Dual Stability System) mechanism from the **Azalea** range. Stability in this wheelchair originates from body length compensation in the back and leg rests, beneficial for users with limited movement and that have difficulty changing position.

High ergonomic push handle

The ergonomically shaped push handles are positioned to optimise the pushing for the carer. An alternate version with higher handles offers a different position for the hands and is optional.



Two-way backrest

The lower and upper parts of the backrest are reversible, giving users the option of more space in the seat or chest area.



Footplate converter

Divided into three parts, the footplate comes with a removable mid section, to facilitate transfers.



Power tilt and recline

With the hand control the user or the carer can easily change seating position.

Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

Azalea Max	550; 610 + 50 mm	530 - 600 mm	400 / 450 mm	500 + 80 mm	320 - 420 mm	330 - 627 mm	1° - 20°	0°-30°
Azalea Max Assist	550; 610 + 100 mm	530 - 600 mm	400 / 450 mm	500 + 80 mm	320 - 420 mm	330 - 627 mm	1° - 20°	0°-30°

Azalea Max	SW + 250 mm	1020 - 1240 mm	1050 - 1135 mm	54 kg	max. 180 kg	33 kg ¹	7°
Azalea Max Assist	SW + 220 mm	1020 - 1240 mm	1050 - 1135 mm	54 kg	max. 180 kg	37 kg ¹	7°


1. Transport weight increases with 4,2 Kg if Power tilt & recline is chosen



Frame colour ▼



Upholstery ▼

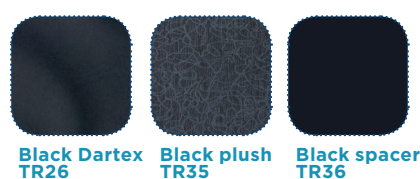

Tel: +41 61 487 70 70
hqeurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Benkenstrasse 260
4108 Witterswil
Switzerland

Rea, Rea design and DSS (Dual Stability System) design are registered trademarks of Invacare International. The **Rea Azalea** is designed with DSS (Dual Stability system). It provides the user with enhanced stability and safety when being tilted or the backrest is reclined.



Yes, you can.®