

Sonata[®]



Engineered for comfort and safety



Whatever a patient's individual needs, the *Invacare Sonata* has been engineered to offer peace of mind, with a flexible range of options designed for maximum comfort.

The **Sonata** is a highly modular and stylish care bed that can be personalised to suit any environment, from a private home to nursing centre. With the **Sonata**, individual care needs and comfort come first. Users can customise the bed by choosing and combining all elements - from the mattress support and legs, to the castors and bed end.

As expected with Invacare, the **Sonata** is sturdy, reliable and built to the highest quality standards.



Easy assembly

The **Sonata** could not be easier to assemble, as all parts are lightweight and portable. The easy-fit legs and simple-assembly bed ends ensure that the **Sonata** can be mounted and transported by one person.



Modular

Different bed designs can be chosen to suit individual tastes and environments. Choose from a wide range of options (bed ends, legs, base, operating method) to ensure their unique and individual care requirements are perfectly catered for.

The flexibility of the *Invacare Sonata* also allows the bed to be upgraded at a later stage, e.g. from a manual to electric operating system.

Personalised comfort

To cater to each individual's need for comfort and support, **Sonata** comes with a choice between two-, three-, or four-sectioned mattress support. Opt for steel slats for easier cleaning or wooden slats for a more flexible level surface. Finally, choose between manual operation or an electric operating system for a more independent and active use of the bed.

Features and options



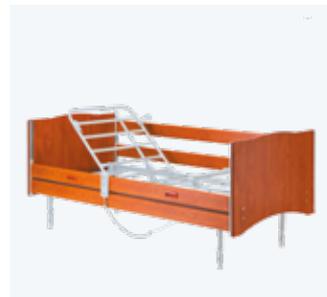
Sonata with Valeria bed ends and Verso II side rails.



Sonata with Sabrina bed ends, Diana side rail and height adjustable legs.



Sonata with electrically operated mattress support, Vibeke bed ends and Verso II side rails.



Sonata with Camila bed ends and Flora side rail.

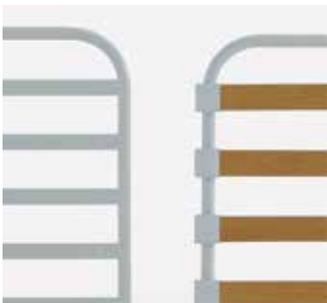
Features and options



Two-sectioned profiling mattress support



Four-sectioned profiling mattress support



Wooden or steel slats.



Movement of lower leg section - Synchronized

Operates together with the upper leg section in a fixed movement.



Handles - Manual

Operated by one handle (two-sectioned support) or two handles (four-sectioned support).



Hand Control - Electric

One electrical function (two-sectioned support) or two electrical functions (four-sectioned base).



Lifting pole on base



Castors

Aluminium 100 mm or 125 mm.



Verso II side rail

A collapsible metal side rail. Compliant with IEC 60601-2-52.



Lifting pole and I.V. Drip Rod

The catching handle for the lifting pole can be adjusted in both height and depth.



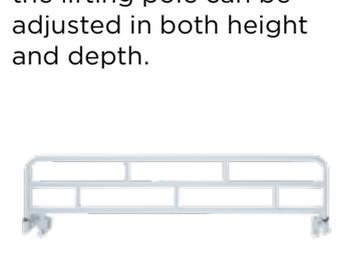
Legs at fixed height of 50 cm. With or without castors.



Legs - telescopic height adjustable (37,5 - 57,5 cm) in 3 steps. With or without castors.



Legs and bed end in one. With or without castors.



Diana side rail






A collapsible metal side rail. Compliant with IEC 60601-2-52.



	<i>Valeria</i>	<i>Nadine</i>	<i>Sabrina</i>	<i>VibekeNG</i>	<i>Camila</i>
Standard wood / colour	Beech laminate	Beech laminate	Beech / Cherry	Cherry	Cherry
Fit in with Verso II	✓	✓	✓	✓	
Fit in with <i>Diana</i>	✓	✓	✓	✓	
Fit in with <i>Flora</i>					✓

Technical data

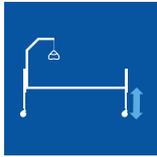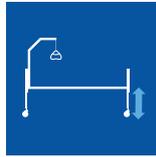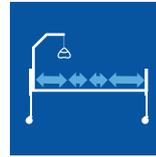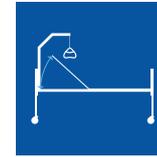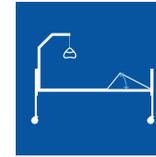For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Width outside! Width inside	Length outside! Length inside	Height adjustment	Height adjustment	Mattress support dimensions	Backrest angle	Thigh angle
Sonata	 900 mm 890 mm Mattress: 900 mm	 2130 mm 2000 mm Mattress: 1900 - 2000 mm	 Height with fixed legs (without castors): 500 mm	 Height with adjustable legs: Without castors: 360, 430, 500 mm With castors ² : 490, 560, 630 mm	 2-sectioned: 800, 1190 mm 4-sectioned: 800, 230, 300, 610 mm	 0 - 72°	 0 - 32°

	Total weight	Max. user weight
Sonata	 2-sectioned: 41 kg 3 and 4-sectioned: 48 kg	 SWL 170 kg

Complies with IEC 60601-2-52

1. With Valeria bed end
2. Castors 100 mm

NB
The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour and wood work



Please notice colors may vary slightly from those displayed above.

Tel: +41 61 487 70 70
hqeurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Benkenstrasse 260
4108 Witterswil
Switzerland

QR
Code

