

SB[®] 910



Designed with safety, ease
of use and comfort in mind




Yes, you can.™

The *Invacare SB 910* bed combines a high specification profiling platform with enhanced patient support, comfort and safety. Its ergonomic design provides ease of use and good working conditions for the carer.

SB 910 has a range of user-friendly features including a synchroniser function operated by a "green button" on the hand control, which allows the backrest and knee break to be operated simultaneously reducing shear and friction whilst profiling.

Caters for individual needs

The *Invacare SB 910* has been developed in close collaboration with health authorities to ensure it meets the requirements and demands of the modern hospital environment. The variable height-adjustment range of the bed from 38 to 83 cm ensures safe patient transfer. The side rails accommodate a range of mattress heights and can be raised and lowered quickly and easily for added safety. The integrated mattress support allows the bed to be adjusted for taller patients.

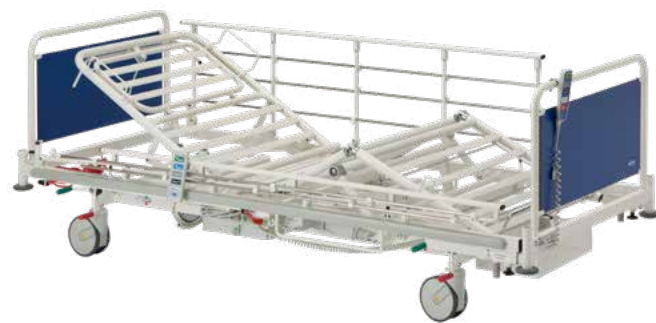

Independence and a high level of patient comfort

The 5-section mattress support with reclining backrest and "green button" function allows the patient to adjust their own sitting and lying position, achieving a high level of comfort along with an optimal reclining position without assistance. This minimises shear and friction forces, thereby preventing skin damage. By using the electrically operated Anti-Trendelenburg feature, which is integrated in the "green button" function, it is possible to achieve an optimal sitting position, providing both comfort and support for the patient.

Flexibility and safety

Security and safety were high on the list of priorities when designing the **SB 910**. The bed comes as standard with an Attendant Control Unit, which enables care and nursing staff to lock one or more of the bed's motor functions. The CPR function immediately brings the bed to the lowest position after placing the mattress platform sections in a horizontal position. The **SB 910** also has a fully welded frame allowing it to be hand or machine washed, meeting every hospital's cleaning and decontamination policy.

Features and options



Hand control
"Green button" function raises or lowers the mattress support to provide an optimum sitting position.



Attendant Control Unit
Allows the carer to lock each of the bed's functions as required. Activates the Trendelenburg/anti-Trendelenburg, the CPR function and other functions of the hand control.



Fender beads
Fender beads on the side of the bed and fender rolls on the corners of the bed for protection.

Features and options



Side rail
Accommodates mattresses from **12-18 cm**, can be raised and lowered quickly and easy for safety.



Side rail height extender
Foldable side rail height extender for use with mattresses from 19 - 30 cm in height.



Integrated mattress support extension
Invacare SB 910 can be extended by up to 20 cm.



IV rack/ telescopic IV rack
IV rack and telescopic IV rack come with three movable hooks.



Basket
Basket for personal belongings, to be placed on the base frame of the bed.



CPR quick release
Quick release handle lowers the backrest to a neutral position in case of an emergency.



Bed hoist - Lifting Pole
The height and depth of the grab handle on the hoist (lifting pole) can be adjusted.



Bed end
The beds ends can be easily and quickly removed to aid care, especially during cardiac arrest.



Neck support
Increased comfort and reduced manual handling for care and nursing staff.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Width outside Width inside	Length outside Length inside	Height adjustment	Mattress support dimensions	Under bed clearance	Backrest angle	Thigh angle
SB 910	800/940 mm 850/990 mm	2000/2160 mm 2050/2210 mm	380 - 830 mm Total height: 1590 - 2040 mm	835 - 280 - 335 - 610 mm	170 mm	0 - 74°	0 - 33°

	Legrest lift	Minimum angle below horizontal	Angle below/ above frame	Total weight	Max. user weight
SB 910	0 - 22°	28°	below: 14° above: 8°	800 mm: 141.5 kg 850 mm: 144 kg	SWL 300 kg 200 kg




Bariatric use

NB
The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Complies with IEC 60601-2-52

Standard frame colour and laminate bed end

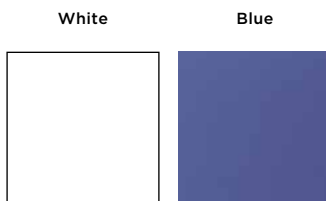

Please notice colors may vary slightly from those displayed above.

Tel: +44 61 487 70 70
hqueurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Neuhofweg 51
4147 Aesch BL
Switzerland



Yes, you can.®