

Kirilos and Atlante



Invacare **Atlante** - electrical version



Invacare **Atlante** - manual version

Compact and easy to maneuver -
the best choice for homecare



Compact in size and offering maximum maneuverability, these mobile lifters provide the ideal solution for the homecare environment. Carers will appreciate the easy handling and smooth operating functions of the **Atlante** and **Kirilos** models.

The **Atlante** lift comes in both an electric and hydraulic version; **Kirilos** is purely electric. Both models come equipped with a handle for easy leg opening.


Features and Options

Easy handling



Spreader bar

Spreader bar rotates 360 degree for easy positioning of the patient.

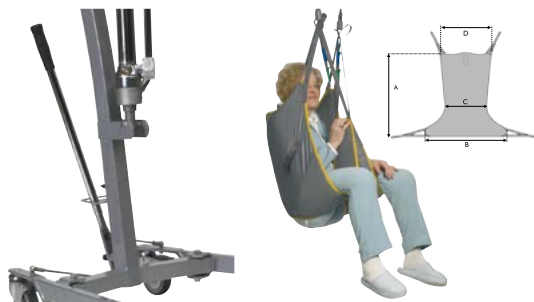

For easy leg opening

Equipped with a handle for comfortable leg opening.

Universal Standard sling, size M, solid

The lifter comes with a Universal Std solid sling (size M). Additional slings are available with or without head support and with or without commode opening.

The **Atlante** and **Kirilos** lifters are designed to make life easier for carers when handling patients. The lightweight, compact nature of the lifters makes them easy to maneuver even when space is limited. The spreader bar rotates 360 degrees, providing maximum handling opportunities for the carer. In addition the lifters can be dismantled for easy transportation and storage.

High safety

The control box of the **Kirilos** and the electric version of **Atlante** are equipped with emergency lowering and emergency stop. The low battery status is displayed by an LED light as well as an acoustic signal.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Overall length of base	Outside width: legs closed	Total product weight	Lifting range	Height of legs	Inside width: legs opened	Max. user weight	Number of cycles
Invacare Atlante hydraulic	1170 mm	700 mm	37 kg	610 - 1530 mm	150 mm	870 mm	140 kg	
Invacare Atlante	1170 mm	700 mm	33 kg	610 - 1530 mm	150 mm	870 mm	140 kg	40 cycles ¹
Invacare Kirilos	1200 mm	680 mm	33 kg	575 - 1490 mm	110 mm	900 mm	140 kg	40 cycles ¹

	Battery capacity	Turning diameter	Wheel diameter	Max. cleaning temperature
Invacare Atlante hydraulic		1230 mm	100 mm	95° C
Invacare Atlante	2,9 Ah	1230 mm	100 mm	95° C
Invacare Kirilos	2,9 Ah	1230 mm	75 mm	95° C

1. without battery charge with batteries at 50 % of full capacity

Frame colour

Grey



QR Code

Tel: +41 61 487 70 70
hqeuropa@invacare.com
www.invacare.eu.com

Invacare International GmbH
 Benkenstrasse 260
 4108 Witterswil
 Switzerland

