




Essential

Basic Cushion

RISK	Comfort
CUSHION WARRANTY	2" - 1 Year 3" - 2 Year
COMPONENT PARTS WARRANTY	1 Year
REPLACEMENTS	All parts
WEIGHT LIMIT	3" - 17 Stone (108 kg) 2" - 12 Stone (76 kg)
CUSHION HEIGHT	2" / 5.08 cm or 3" / 7.62 cm
CUSHION SIZE RANGE	8 size options available: 15x16 / 16x16 / 16x18 / 17x17 18x16 / 18x17 / 18x18 / 19x17" See cushion size chart

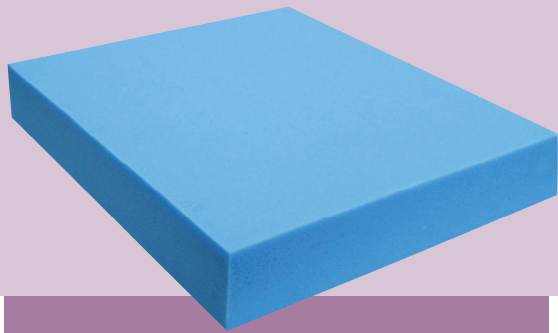

BASIC COMFORT CUSHION

FOR COMFORT

The Invacare Essential Basic Cushion provides basic comfort thanks to a highly resilient foam. The cushion is available in both 2" and 3" cushion depths, and comes with a two-way stretch, vapour permeable, ingress resistant black cover.

FEATURES

Cover features a zip

Available in two heights

Quality foam core, fire tested to Crib 5
Cover tested to BS ISO 7176-16:1997

Designed for purpose

Two way stretch, ingress resistant, vapour permeable cover.

Cover can be unzipped and removed from foam core. Cover is machine washable at temperature up to 80°C.

All covers have an anti-slip base.

BENEFITS

Enables easy inspection of foam core.

2"/ 5.08 cm and 3"/ 7.62 cm depth options

Offers comfort and assurance.

Delivers excellent value for money.

Helps to prevent a build up of shear and friction forces against the user. Allows fluid to 'roll off' fabric, therefore protecting clothing and foam core.

The cover material has been treated with a fluoro-carbon-based fluid repellent finish in order to prevent wetting of the cover and therefore reduce soiling and the risk of contamination.

Helps prevent slippage, ensuring safe use of the cushion.


* Essential nursing care is pivotal in pressure ulcer prevention. These cushions will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action planning based on vulnerability are fundamental factors in the prevention of pressure ulcers.



Essential

Visco Cushion

RISK	Comfort
CUSHION WARRANTY	2 Year
COMPONENT PARTS WARRANTY	1 Year
REPLACEMENTS	All Parts
WEIGHT LIMIT	17 Stone (108 kg)
CUSHION HEIGHT	3" / 7.62 cm
DIMENSIONS	8 size options available: 15x16 / 16x16 / 16x18 / 17x17 18x16 / 18x17 / 18x18 / 19x17" See cushion size chart

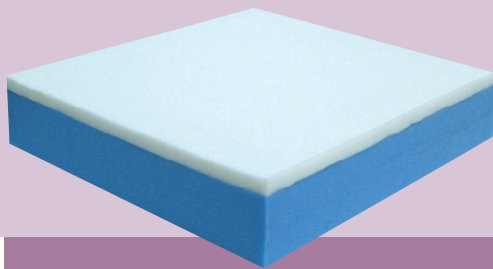

VISCO COMFORT CUSHION

FOR COMFORT

The Invacare Essential Visco Cushion is dual layered comfort, using a combination of high density, temperature sensitive, visco-elastic foam in conjunction with a highly resilient foam base to provide added support and longevity. The cushion has a two-way stretch, vapour permeable, ingress resistant black cover.

FEATURES

High density visco elastic foam top: weight and temperature sensitive

Single piece high density foam base.

Cover features a zip

Quality foam core, fire tested to Crib 5
Cover tested to BS ISO 7176-16:1997

Designed for purpose

Two way stretch, ingress resistant, vapour permeable cover.

Cover can be unzipped and removed from foam core. Cover is machine washable at temperature up to 80°C.

All covers have an anti-slip base.

BENEFITS

Conforms to the users shape, affording exceptional comfort.

Provides improved support and prolonged longevity.

Enables easy inspection of foam core.

Offers comfort and assurance.

Delivers excellent value for money.

Helps to prevent a build up of shear and friction forces against the user. Allows fluid to 'roll off' fabric, therefore protecting clothing and foam core.

The cover material has been treated with a fluoro-carbon-based fluid repellent finish in order to prevent wetting of the cover and therefore reduce soiling and the risk of contamination.

Helps prevent slippage, ensuring safe use of the cushion.

* Essential nursing care is pivotal in pressure ulcer prevention. These cushions will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action planning based on vulnerability are fundamental factors in the prevention of pressure ulcers.