

Dacapo Range

Quality Foam Mattresses




Flexibility & comfort



Invacare **Dacapo** is a broad range of mattresses designed to match customer's requirements within Invacare's homecare and hospital bed ranges. Selecting from high resilience, environmentally friendly, or high density visco-elastic foams within the 'Standard', 'Comfort' and 'Preventive' ranges, **Dacapo** has been developed to offer optimal cushioning and pressure relief without compromising on safety, hygiene and user friendliness.

Invacare **Dacapo Standard**

Dacapo Standard is a standard mattress comprehensively satisfies basic requirements for good cushioning. High quality, high resilient foam provides a resilient & firm supportive base.
10cm of high resilience foam.

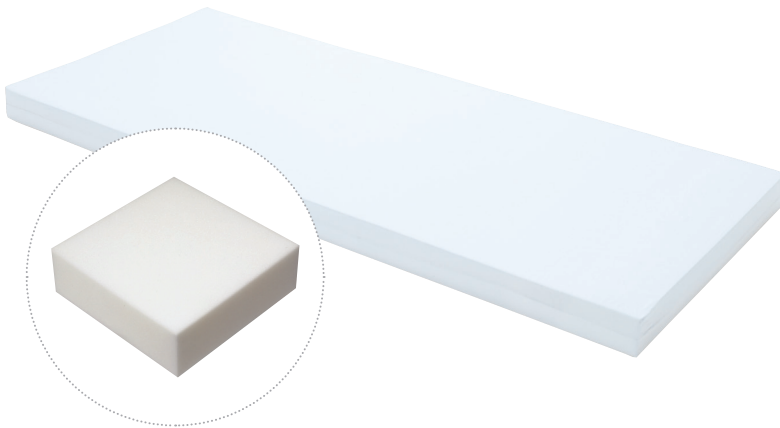

Dacapo Standard Cover

- ▶ 100% cotton
- ▶ Washable at 40 degrees

Invacare **Dacapo Basic**

Dacapo Basic is the first incontinence mattress in the series. It fully satisfies basic comfort & hygiene requirements for a sleeping surface. The high density, high resilient, foam provides a firm supportive base, strength & durability.
12cm of high resilience foam.




Dacapo Basic Cover

- ▶ 100% Fluid-proof
- ▶ Polyurethane coated polyester
- ▶ Breathable
- ▶ Two way stretch capacity
- ▶ Washable at 95 degrees

Invacare **Dacapo Heavy User Light**

Dacapo Heavy User Light is tailored to meet the needs of the larger user and offers a firm, comfortable cushioned support.

Uniquely constructed using two high resilience foams and a layer of visco elastic foam, this model can be reversed to adjust to the degree of firmness and support required by the user. The broader dimensions enjoyed by **Heavy User Light** make it the perfect complement for Invacare's Broad Home Care Bed models.

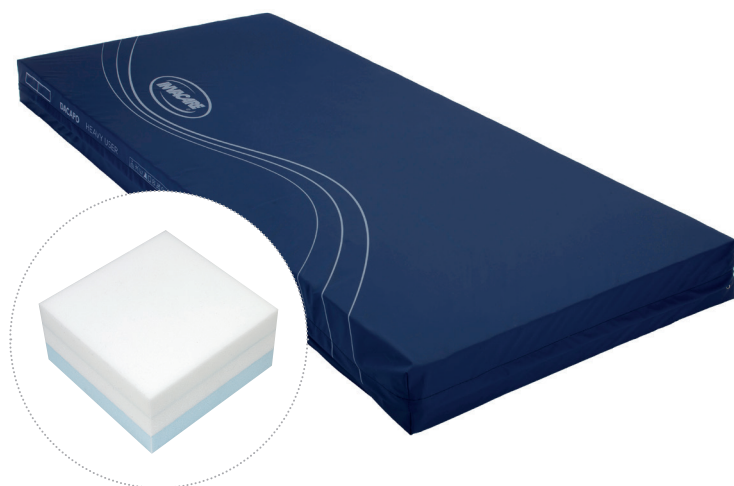

Dacapo Heavy User Light Cover


- ▶ 100% fluid-proof
- ▶ Polyurethane coated polyester
- ▶ Breathable
- ▶ Two way stretch capacity
- ▶ Washable at 95 degrees

Dacapo Heavy User Light is available in a variant that will offer support up to 350kg.

Invacare **Dacapo Comfort**

The **Dacapo Comfort** is a comfortable mattress comprised of a foam core whose density provides good recovery while the mattress' unique open cell structure, allows the skin to breathe. The mattress is constructed with great flexibility in order to provide a comfortable and supporting surface. The mattress is cut into various sized squares that are individually flexible. This profiling allows the contours of the body to immerse into the material so that optimum sleeping comfort is achieved.

12cm of environmentally friendly high resilience foam.

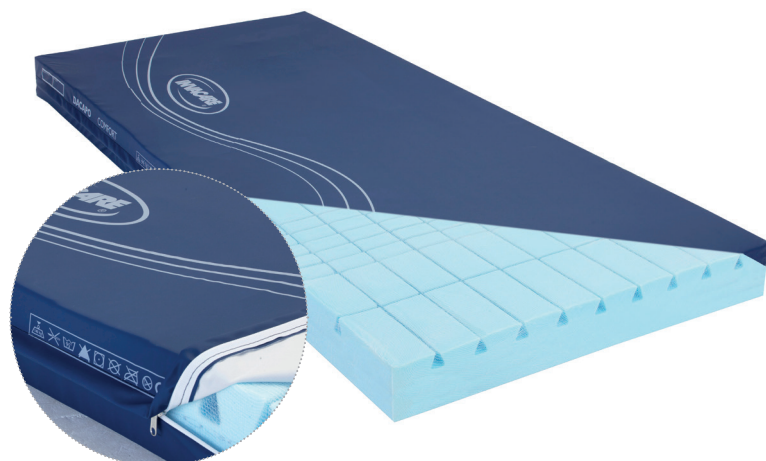


Dacapo Comfort Cover

- ▶ 100% fluid-proof
- ▶ Polyurethane coated polyester
- ▶ Breathable
- ▶ Two way stretch capacity
- ▶ Washable at 95 degrees

Invacare **Dacapo Comfort Plus**

The **Dacapo Comfort Plus** is a prevention mattress which satisfies the basic requirements for a comfortable sleep. A unique pressure relieving mattress consisting of various large keyhole cut squares that are individually flexible. The mattress has a good recovery and its unique open cell structure allows the body to breathe. The mattress offers a dynamic combination of comfort, support and stability while shear and friction forces are efficiently minimized due to the keyhole cut design. 14cm of environmentally friendly high resilience foam.

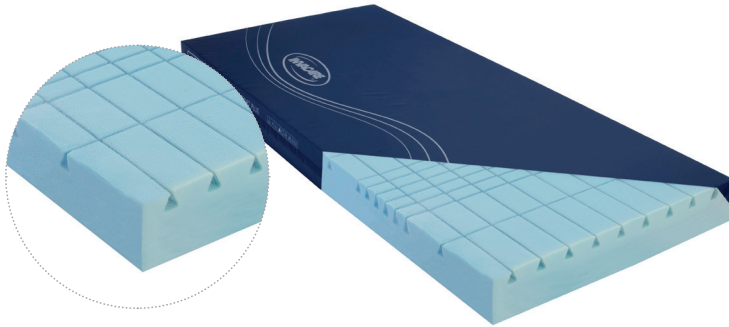


Dacapo Comfort Plus Cover

- ▶ 100% fluid-proof
- ▶ Polyurethane coated polyester
- ▶ Breathable
- ▶ Two way stretch capacity
- ▶ Washable at 95 degrees

Invacare **Dacapo Square**

Dacapo Square is a combination mattress with an environmentally- friendly foam top. The model's unique square-patterned body contour cut minimises surface pressure, affording the user outstanding pressure relief and body pressure equalisation. The foam's high volume weight gives it outstanding responsive & recovery properties, while its unique vapour permeable cell structure guarantees a well-aerated interior. In terms of user benefits, **Dacapo Square** excels in its ability to conform to every curve of the body and provide total body support.

Top (8 cm): Environmentally friendly high resilience foam. **Base (6 cm):** High resilience foam.

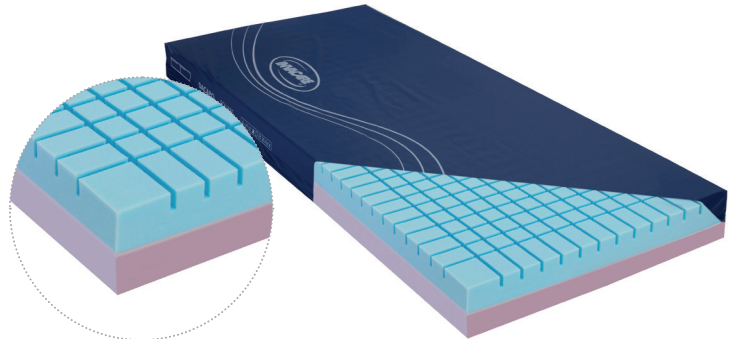

Dacapo Square Cover

- ▶ 100% Fluid proof
- ▶ Polyurethane coated polyester
- ▶ Breathable
- ▶ Two way stretch capacity
- ▶ Washable at 95 degrees

Invacare **Dacapo Combi**

Dacapo Combi - is a combination mattress with outstanding pressure relief properties. The top foam is composed of a self-regulating, heat sensitive material which works in conjunction with body temperature & movement. **Dacapo Combi** relaxes according to the user's body dynamics and softens as the user gently sinks into it, comfortably moulding & conforming to the body shape. The open cell structure negates heat & moisture build-up, while the base high resilient, foam reinforces the mattress strength & robustness.

Top (6 cm): Visco-elastic foam. **Base (8 cm):** High resilience foam.



Dacapo Combi Cover

- ▶ 100% fluid proof
- ▶ Polyurethane coated polyester
- ▶ Breathable
- ▶ Two way stretch capacity
- ▶ Washable at 95 degrees

Invacare **Dacapo 2Sided**

Dacapo 2Sided - is a combination mattress with outstanding pressure relief properties. One side is made from environmentally friendly organic foam. The high density of the foam allows for a good pressure distribution. The upper surface has been cut into squares, all being individually flexible. The other side is composed of visco-elastic foam which adapts to the user's body dynamics and temperature.




Dacapo 2Sided Cover

- ▶ 100% fluid proof
- ▶ Polyurethane coated polyester
- ▶ Breathable
- ▶ Two way stretch capacity
- ▶ Washable at 95 degrees

Invacare **Dacapo Premier**

The **Dacapo Premier** features a unique contoured insert pad which provides excellent patient comfort and optimises pressure reduction. The mattress cover is two-way stretch and helps to reduce shear and friction forces. The insert is supported by a foam u-core that facilitates patient transfer by providing effective side-wall support.

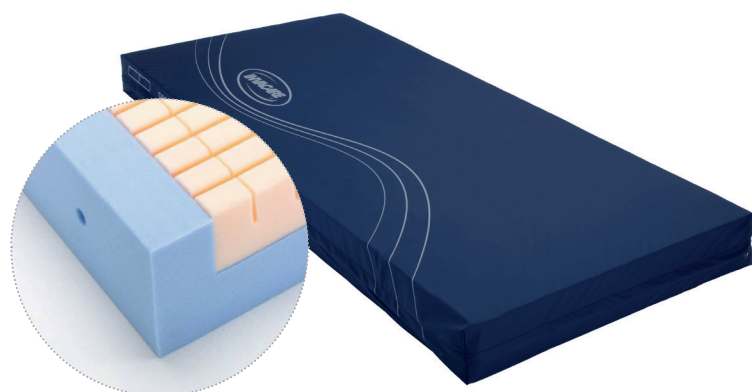


Dacapo Premier Cover

- ▶ 100% fluid proof
- ▶ Polyurethane coated polyester
- ▶ Breathable
- ▶ Two way stretch capacity
- ▶ Washable at 95 degrees

Invacare **Dacapo Top**

Dacapo Top - is a high quality, pressure relieving top layer ideally suited to basic mattresses. **Dacapo Top's** 'temperature smart' foam adapts to the user's body dynamics & temperature. On impact, the material relaxes and then moulds & conforms to the contours of the body leaving the user cradled in comfort. The foam's open cell structure affords a wellventilated interior which reduces the build up of heat & body moisture.

6 cm: Visco-elastic foam.

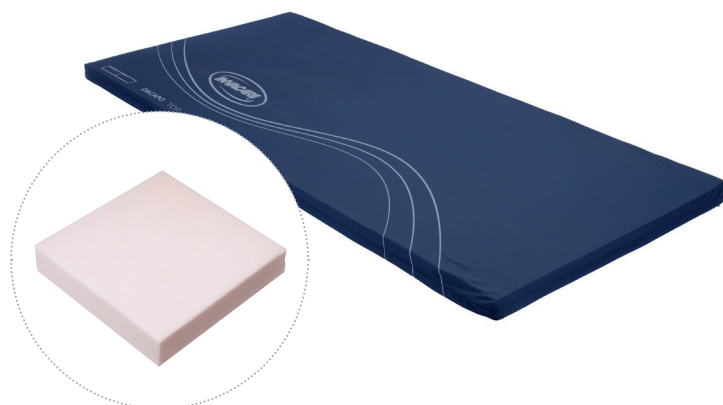


Dacapo Top Cover

- ▶ 100% fluid proof
- ▶ Polyurethane coated polyester
- ▶ Breathable
- ▶ Two way stretch
- ▶ Washable at 95 degrees

Invacare **Dacapo Top Pro**

The **Dacapo Top Pro** is a high quality pressure relieving mattress and has a foam core which has independent, geometrically-cut surface cells. The foam core has larger, more supportive cells at the centre and smaller cells at each end. This enables the mattress overlay to achieve even pressure redistribution by conforming to the shape and form of the body and cradles bony prominences.

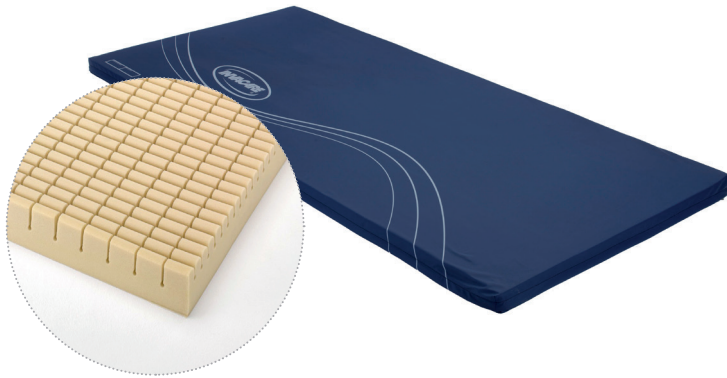


Dacapo Top Pro Cover

- ▶ 100% fluid proof
- ▶ Polyurethane coated polyester
- ▶ Breathable
- ▶ Two way stretch
- ▶ Washable at 95 degrees

Invacare **Dacapo Child**

Dacapo Child a single layer mattress constructed from an environmentally-friendly foam. The mattress is designed for the smaller, younger user and provides extra cushioning, ultra-durability and enhanced ventilation.

Top (8 cm): Environmentally friendly high resilience foam

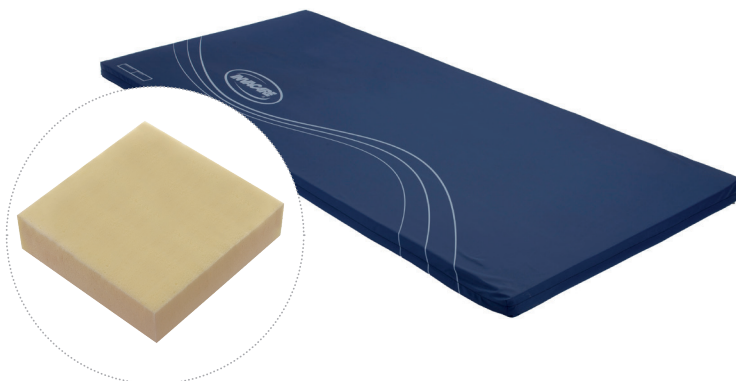


Dacapo Child Cover

- ▶ 100% fluid proof
- ▶ Coated polyester with polyurethane*
- ▶ Breathable
- ▶ Two way stretch capacity
- ▶ Washable at 95 degrees

*Also available with cotton cover option.

Invacare **Dacapo Kid**

The **Dacapo Kid** mattress is specifically designed for the **NordBed Kid**, available in 6 sizes to fit the bed platform variants. Based on the **Dacapo Child** design but with a single-layer high resilience foam, **Dacapo Kid** features square castellations which reduces shear and redistributes pressure as well as enhancing ventilation; ultimately providing a comfortable and durable solution for young children.

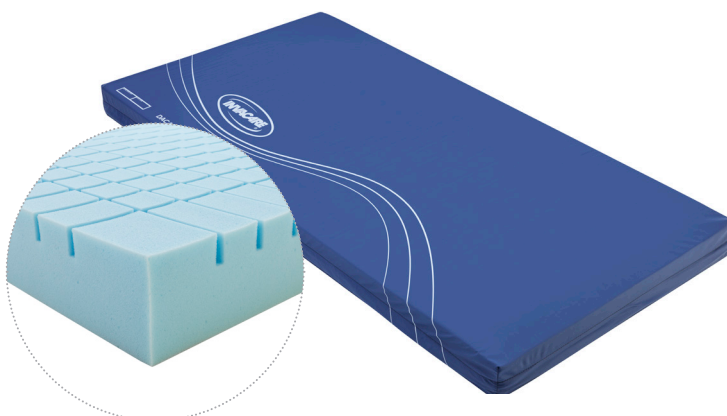

Dacapo Kid Cover

- ▶ Ingress resistant
- ▶ Breathable polyurethane cover
- ▶ Two way stretch capacity
- ▶ Washable at 95 degrees
- ▶ Fully concealed zip

Key features and options

◀ Cover options



Wipdown cover

The wipdown cover is composed of polyutherane coated polyester. The cover is 2 way stretch and gentle on the skin. The fabric is vapour permeable and waterproof. This cover is allergy tested and treated against dust mites.

Evacuation cover

The Evacuation cover is composed of 2 way stretch, vapour permeable, waterproof, polyurethane coated polyester. These covers are equipped with handles on all four sides. On the two long sides are straps, which in case of an emergency evacuation can be wrapped around the patient. This cover is allergy tested and treated against dust mites.

Cotton cover

The Cotton cover is composed of standard cotton and washable at 40 °C.



◀ High resilience foam

Polyurethane foams vary in cell structure and density. High Resilience foam has relatively uniform and dense cell structure that allows it to provide good support while loaded, and also maintain its shape while unloaded.



◀ Easy Handling

Optional handles further facilitate easy handling of the Dacapo range.



◀ Zipper

Discreetly zippered on three sides with cloth covering.

Technical data ▼


For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.



Product total weight ▼



Max. user weight ▼

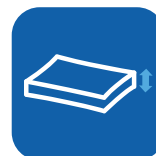

Height ▼

	Product total weight ▼	Max. user weight ▼	Height ▼
Standard	7 - 10.5 kg	100 kg	100 mm
Basic	8.5 - 12.5 kg	130 kg	120 mm
Heavy User Light	11 - 16.5 kg	250 kg	150 mm
Heavy User	12 - 17.5 kg	350 kg	150 mm
Comfort	11.5 - 17kg	130 kg	120 mm
Comfort Plus	13.5 - 19kg	140 kg	140 mm
Square	10 - 15 kg	140 kg	140 mm
Combi	8.5 - 13 kg	140 kg	140 mm
2Sided	11.5 - 17 kg	130 kg	150 mm
Premier	14 kg	180 kg	150 mm
Top	5 - 7.5 kg	130 kg	60 mm
Top Pro	4.5 - 5 kg	150 kg	70 mm
Child	3 kg	100 kg	80 mm
Kid	6 kg	70 kg	100 mm

Fire retardancy conforming to European Standard EN 597 – 1 and European Standard EN 597 – 2 is applicable for all Dacapo models with a polyurethane covered cover.

Tel: +41 61 487 70 70
hqeuropa@invacare.com
www.invacare.eu.com

Invacare International GmbH
 Benkenstrasse 260
 4108 Witterswil
 Switzerland

