

# Action<sup>®2</sup> Basic



Simple but without compromise



The **Action<sup>2</sup> Basic** is a simple but highly durable wheelchair featuring a reclining backrest. With **three standard sizes to choose from**, it is quick and easy to find the right size for your needs. In addition, the **chair weight has been kept low thanks to a balance of steel and aluminum components**. When removing the legrests, armrests and rear wheels the **Action<sup>2</sup> Basic** can be stored in even the smallest of spaces and will fit in most modern cars.

- **Improved comfort**
- **Proven reliability and quality**

## Improved comfort

With a simple handle the **backrest reclines in two quick steps**, providing **additional comfort** for those with hip flexion limitation and **accommodating individual body shape**.


## Proven reliability and quality

With years of experience and using **tried-and-tested technology**, the **Action<sup>2</sup> Basic** is a reliable partner. The one you can count on.

## Features and Options

### Elevating legrests

The optional Eco elevating legrests provide **extra comfort for users**.

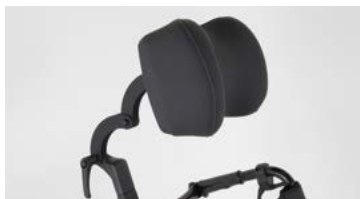

### Headrest

An optional headrest **supports the head for extra comfort** when the backrest is reclined.



### Anti-tip device

The anti-tip device provides **extra safety on slopes** or when backrest is reclined.

## Technical data

**For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.**

local Invacare website.								
	Seat width	Seat depth	Seat angle	Seat height	Backrest angle	Backrest height	Front frame angle	Length of footrests
Self Propulsion	405 / 455 / 505 mm	425 mm	3°	495 mm	0° / 7° / 0°-15°-30°	430 mm	70°	330 - 455 mm
	Armrest height	Driving unit width <sup>4</sup>	Total width, folded	Total length inc. Footrests <sup>4</sup>	Total length without footrests	Total height	Max. Safe slope <sup>1</sup>	Radius <sup>2</sup>
Self Propulsion	240 mm	600 - 708 mm (Seat width +210 mm)	280 mm	1116 - 1150 mm	860 mm	945 mm		
	Max. User weight	Total weight	Mass of the heaviest part <sup>3</sup>	1. Stability Values Max. slope parking brake: 7° Static stability downhill: 9° Static stability uphill: 21,5° Static stability sideways: 16,8°		2. Radius Mini turning radius: 900 mm Pivot width: 1380 mm		
Self Propulsion				3. Folded Dimensions Folded height: 635 mm Folded depth: 1065 mm Folded width: 225 mm		4. Overall dimensions In some configurations, overall width or length exceed the authorised limits, access to		
	125 kg	from 18 kg	from 10 kg					This wheelchair is not intended to be used

## Frame colour

Pearl Grey



## Upholstery

Black nylon



Please notice colors may vary slightly from those displayed above.

QR Code

Tel: +41 61 487 70 70  
hqeurope@invacare.com  
www.invacare.eu.com

Invacare International GmbH  
Benkenstrasse 260  
4108 Witterswil  
Switzerland



**Yes, you can.®**



© 2018 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. Action2 Basic - EU - 10/2018