



Azalea® Tall



No limits on comfort for
taller users



Rea Azalea Tall is specially designed to meet the needs of tall users who require a Tilt in Space wheelchair with a longer seat support. Adapted from the **Rea Azalea**, the **Rea Azalea Tall** boasts all the advantages of a reliable, tilting wheelchair and offers a unique weight-shifting mechanism. The **Rea Azalea Tall**: at the top of its class for comfort and stability. The improved Flex 3 backrest adds individual adaptability. 5 frame colours and three upholsteries are available.

Real compatibility

The **Rea Azalea Tall** has all the features of the **Rea Azalea** such as Tilt in Space weight-shifting for a really stable tilt. All accessories such as rear wheels for self-propelling or smaller wheels for transport and a different seat system also fit the **Rea Azalea Tall**.


Easy maneuverability

The adjustable push bar has been designed with the carer in mind, it offers a good ergonomic grip and can be adjusted to suit the needs of a carer.



Longer seat plate

An extra 50 mm of seat depth is designed to support longer legs.

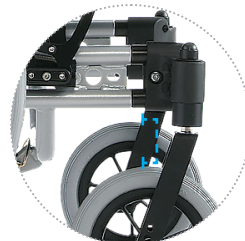

Chassis length

An extra 50 mm in length enables greater forward stability.






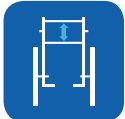






Seat height

An extra 50 mm of height allows for longer leg-rests.

For other Rea Azalea accessories, please consult the Rea Azalea brochure or visit www.invacarearea.com to see videos of them

Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

								
	Seat width ▼	Seat depth ▼	Seat height ▼	Backrest height ▼	Armrest height ▼	Length of footrests ▼	Seat angle ▼	Backrest angle ▼
Azalea Tall	390-450/ 440-500/ 490-550 mm (390-490/ 440-540/ 490-590 mm)	480 - 550 mm	500 mm	Backrest plate: 560 - 790 mm Tension adjustable backrest: 540 - 715 mm	240 - 360 mm	330 - 627 mm	-1° to 23° (+1° to 25°)	0° - 30°
								
	Driving unit width ▼	Total height ▼	Total length inc. Footrests ▼	Total weight ▼	Max. User weight ▼	Transport weight ▼	Crash tested ▼	Maximum slope for use of parking brake ▼
Azalea Tall	SW + 220 mm (SW + 230 mm)	1000 - 1300 mm	950 - 1085 mm	34 kg (SW 440 mm)	Max 135 kg	18,7 kg ¹	ISO 7176-19	7°

Frame colour ▼



Pearl Grey


Azurite Blue

Sand

Happy Red

Electric Green

Upholstery ▼

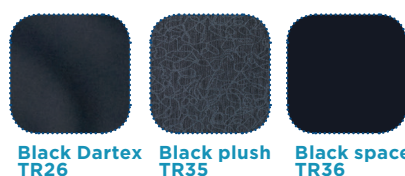

Black Dartex TR26

Black plush TR35

Black spacer TR36

1. Transport weight increases with 4,2 Kg when Power tilt & recline option is used



Tel: +41 61 487 70 70
hqeurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Benkenstrasse 260
4108 Witterswil
Switzerland

Rea, Rea design and DSS (Dual Stability System) design are registered trademarks of Invacare International. The **Rea Azalea** is designed with DSS (Dual Stability system). It provides the user with enhanced stability and safety when being tilted or the backrest is reclined.



Yes, you can.®