

# ***Flo-tech<sup>®</sup> Plus***



Deep contoured cushion with gel sac



Yes, you can.

The Invacare **Flo-tech Plus** cushion incorporates a fluid sac integrated in the cushion seat well. The positioning of the fluid sac targets the ischial tuberosities, helping to significantly reduce pressure in this vulnerable area. This delivers superior levels of protection for those patients identified as being at High Risk of developing pressure ulcers.

## Features and Options



### Ingress Resistant Cover

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle.



### Available in 80 sizes

The *Invacare Flo-tech Plus* is available in a range of 80 sizes, from 360 x 330 to 600 x 550 mm. (See size chart for details).



### Integral Sag Compensator option

Option of sag compensator to accommodate the sag in wheelchair seats, and ensure the cushion is correctly presented.



### Ultra-fresh™ added to all cushions for antimicrobial and odour protection

Antimicrobial treatment that has been added to all **Flo-tech** foam to help make it resistant to the growth of odour causing bacteria. Helps provide longer lasting freshness.

## Contoured seating surface

Offers improved stability and postural support. Also delivers a greater degree of pressure distribution compared to that of a flat surface. The castellations of the cushion allow slight movement which in turn reduces shear and friction forces. The deeper moulded foam offers greater protection to those who are considered clinically obese, or users who would benefit from a higher sitting position.

## Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

Various sizes	Height	Total product weight	Max. user weight <sup>1</sup>	Max. cleaning temperature
<b>Flo-tech Plus</b> A range of 80 sizes	103 mm	2.5 kg	152 kg 191 kg on sizes 560 x 560 mm, 560 x 610 mm	80°

Manufactured to European Standards:  
BS EN 1021-1 Cig, BS EN 1021-2 Match  
and BS 7176 (Crib 5)

1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.

## Colours

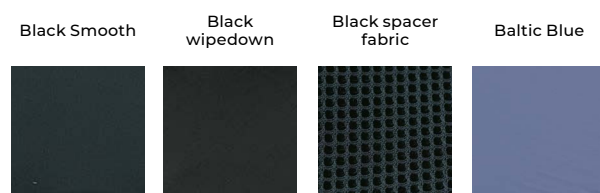

**Tel: +41 61 487 70 70**  
hqeurope@invacare.com  
www.invacare.eu.com

Invacare International GmbH  
Neuhofweg 51  
4147 Aesch BL  
Switzerland



**Yes, you can.®**