

Leo[®]



Safe, solid, secure



The *Invacare Leo* is a scooter designed for all those that value their independence and wish to get out and about unaided. Safety is a key feature of the **Leo**, but this does not detract from its stylish and sporty looks. **Leo** offers users the freedom and confidence to enjoy their essential daily outings and leisure excursions. With exactly high standards of design and safety engineering, the **Leo** undoubtedly sets a new market standard.



For a safe, powerful drive

The *Invacare Leo* offers a smooth, safe drive in a scooter that is easy to control and manoeuvre. The **Leo** intelligent lighting package includes indicators that automatically switch off and a rear brake light for safe driving. In addition to a visual warning for the battery level display, users are also alerted with an audible warning (beep) when the battery needs recharging.

Reliable

Users can trust the **Leo** to take them wherever they need to go. From a trip to the shops or a ride to the park, the **Leo** is rendered sturdy and reliable thanks to its exceptionally high quality electronic engineering and power. These elements are protected from corrosion and water, (splash protection) which safeguard the vital workings of the scooter against dirt and water. A two-step disengaging lever protects the scooter from free-wheeling if inadvertently knocked out of drive mode. The lighting system is designed so that if one of the lights goes, all other lights still continue to function.

Features and Accessories



Finite tiller movement

Users can easily adjust the tiller to suit their needs.



Brake light

Rear brake lighting warns other people that the scooter is slowing down, thus preventing collisions.



Comfort seating

Swivel functionality and sliding rails ensure optimal seating position.



Two-step disengaging lever

This unique concept prevents the scooter from accidentally freewheeling.



Compact

The **Leo** is a compact mini-scooter that can be easily dismantled for transportation during longer trips or weekend breaks. With a smooth drive and easy manoeuvrability, getting round tight corners in the local shopping centre, museum or leisure complex becomes refreshingly easy. Access to the battery and electronics for maintenance or servicing is much easier with the **Leo** as only one part needs to be removed.

Features and Accessories



Built in anti-splash guards

Protects electronics and the transaxle from dirt and water.



Disassembly into five parts

No tools are required to dismantle *Invacare Leo*



Crutch & cane holder

Keeps your crutch or cane stable during the drive.



Lockable front box

Secure box for personal valuables.



Seatbelt

Keeps the user well-positioned in the seat.



Large mirror

Can be mounted left or right



Heavy-duty storage cover

Keeps the scooter clean and dry.



Charging off-board

Simply remove batteries to charge with the off-board charger kit.



Seat suspension module

Adds extra comfort during driving.



Windscreen

Protects the user from the cold wind.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.


	Backrest width	Seat depth	Seat to floor height	Backrest height	Driving unit width	Total length	Empty weight
<i>Invacare Leo</i>							
4 wheels	470 mm	410 mm	460 - 560 mm	475 mm	590 mm	1220 mm	59 kg (excl. batteries) 83 kg (incl. batteries)
3 wheels	470 mm	410 mm	460 - 560 mm	475 mm	590 mm	1200 mm	54 kg (excl. batteries) 78 kg (incl. batteries)

	Max. User weight	Battery capacity	Turning diameter	Max. Climbable obstacle height	Ground clearance	Max. Safe slope	Motor capacity
<i>Invacare Leo</i>							
4 wheels	136 kg	2 x 36 Ah	2620 mm	60 mm	70 mm	10°	240 W
3 wheels	136 kg	2 x 36 Ah	2320 mm	60 mm	70 mm	10° 1	240 W

	Max. Driving distance	Max. Driving speed
<i>Invacare Leo</i>		
4 wheels	38 km	8 km/h
3 wheels	36 km	8 km/h


1. Seat in maximum fifth rearward locking position from the front.

Colour of shrouds

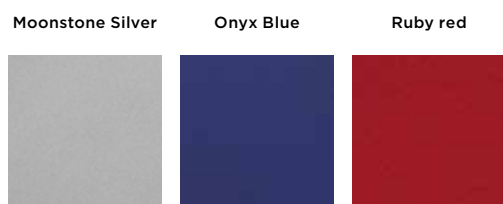

Upholstery

Black vinyl



Tel: +41 61 487 70 70
hqeuropa@invacare.com
www.invacare.eu.com

Invacare International GmbH
 Benkenstrasse 260
 4108 Witterswil
 Switzerland

