

Propad™ Premier Cushion




Profiled foam cushion designed
for comfort and support



Yes, you can.®

The *Invacare Propad Premier Cushion* has been designed to facilitate good pressure reduction when used with both wheelchairs and standard armchairs. This cushion is constructed to ensure improved support and prolonged longevity.

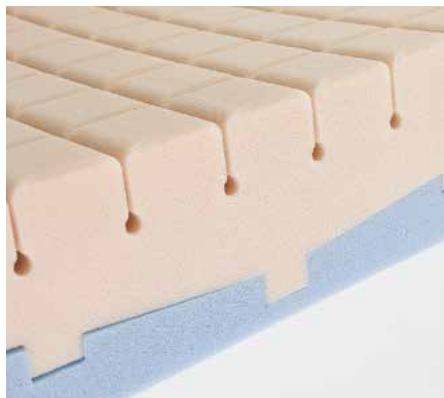

Features and Options



Cover helps to prevent a build up of shear and friction forces against the user

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base.



Available in a range of sizes

The *Invacare Propad Premier Cushion* is available in three standard sizes: 410 x 410, 430 x 430 and 460 x 410 mm. (See size chart for details).



Full length flap conceals the zip

The full length flap concealing the zip reduces the potential for fluid ingress into the cushion and protects the foam core. The cover can be removed and laundered.

Deep contoured dual-layer profiled foam cushion designed to deliver enhanced pressure care

The *Invacare Propad Premier Cushion* is a dual layer cushion. The lower layer is constructed from a firmer foam, shaped to deliver increased stability support and prevent bottoming out. The upper layer features a layer of castellated foam, each castellation delivering enhanced support and pressure area care to the user. The “keyhole” cut at the base of each castellations ensures good air circulation, reducing the potential for skin maceration.

Castellated seating surface

Offers improved stability and postural support; castellations reduce shear and friction forces on tissue.


Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Size range	Height	Total product weight	Max. user weight	Max. cleaning temperature
Propad Premier Cushion	406 x 406 mm 432 x 432 mm 406 x 457 mm	102 mm	1 kg	108 kg	80°

Manufactured to European Standards:
BS EN 597-1 Cig, BS EN 597-2 Match
and BS 7175 (Crib 5)

Colours

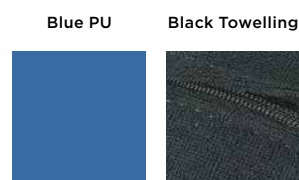

QR
Code

Tel: +44 61 487 70 70
hqeurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Neuhofweg 51
4147 Aesch BL
Switzerland



Yes, you can.®