

rea[®]

Azalea[®] Base



The perfect wheelchair base for
custom made seat units




Rea Azalea Base is designed for users seeking a custom-made seat combined with a first-class wheelchair base. The **Rea Azalea Base** is wider in the rear section of the seat, offering the possibility to use larger seating systems on smaller frames. As there are many balance options, it is now lighter to push. A key feature of the Base is the 'Tilt in Space' function - a weight shifting mechanism that ensures reliable stability for the user. All the advantages of our standard **Rea Azalea** models are included and is also cross-compatible with other **Rea Azalea** models (excluding trunk supports). Available in 5 colours.

Seat adaptation

Rea Azalea Base is specifically designed to accommodate custom-made seats. The backrest is angle adjustable including an increased recline which can be angled via a 'rigging screw' or an attendant operated gas piston.

The seat width can be adjusted from the arm and leg rests, and the depth changed via the telescopic seat frame. These features ensure that the **Rea Azalea** easily adapts to different-shaped seat moulds.

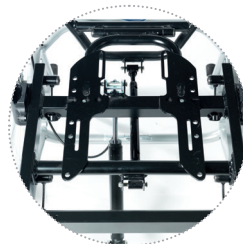

◀ Affixing attachments

The seat offers various attachment holes, plus greater 'open space' at the back.



◀ Harness attachment

An attachment for different kind of harnesses is also available as an option.



◀ Seat plate with turnable wedge

An additional +5° seat angle can be obtained, using the optional seat plate with a wedge.

For other Rea Azalea accessories, please consult the Rea Azalea brochure or visit www.invacarearea.com to see videos of them

Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

Seat width ▼	Seat depth ▼	Seat height ▼	Backrest height ▼	Armrest height ▼	Length of footrests ▼	Seat angle ▼	Backrest angle ▼	
Azalea Base	340/390/ 440/490 mm (+ 50/100 mm)	380 - 450 mm ¹ 430 - 500 mm	400/450 mm	600 ¹ / 700 mm	240 - 360 mm	330 - 627 mm	-1° to 23° (22 / 24") (+1° to 25°) (16")	0° - +30°
Driving unit width ▼	Total height ▼	Total length inc. Footrests ▼	Total weight ▼	Max. User weight ▼	Transport weight ▼	Crash tested ▼	Maximum slope for use of parking brake ▼	
Azalea Base	SW + 220 mm (SW + 230 mm)	1300 ¹ / 1400 mm	950 - 1035 mm	20 kg	Max 135 kg	15 kg ³	ISO 7176-19 ²	7°

Frame colour ▼



Pearl Grey

Azurite Blue

Sand

Happy Red

Electric Green



1. For a seat width 340 mm
2. Crash tested in specific configurations
3. Transport weight increases with 4,2 Kg if Power tilt & recline is chosen

Rea, Rea design and DSS (Dual Stability System) design are registered trademarks of Invacare International. The Rea Azalea is designed with DSS (Dual Stability system). It provides the user with enhanced stability and safety when being tilted or the backrest is reclined.

Tel: +41 61 487 70 70
hqeurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Benkenstrasse 260
4108 Witterswil
Switzerland



Yes, you can.®