



Invacare®

# Propad™ Mattress Overlay

High specification pressure  
redistributing overlay



The **Propad Overlay** is designed to sit on top of an existing conventional mattress, providing a high level of patient comfort with excellent pressure redistribution abilities. The overlay is lightweight, available in a range of sizes and suitable for those considered at 'High Risk'\* of developing pressure ulcers, making it ideal for the homecare environment.


## Effective Pressure Ulcer Prevention

The foam core has independent, geometrically-cut surface cells which significantly reduce shear and friction, whilst increasing the overall support surface area, thereby reducing peak pressure.

## Even Weight Distribution

The **Propad Overlay** has larger, more supportive cells at the centre and smaller cells at each enabling the mattress overlay to achieve even weight distribution by conforming to the shape and form of the body and cradle bony prominences.

## Features and Options

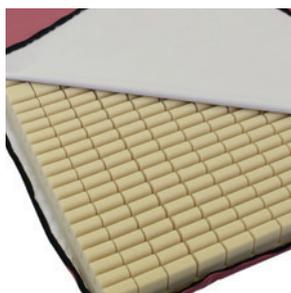


### Cover helps to prevent a build up of shear and friction forces

Cover is multi stretch, water-resistant and vapour-permeable.



### Available in a range of sizes

Single - 1875 x 876 x 85 mm,  
Double - 1875 x 1372 x 85 mm,  
King Size - 1980 x 1524 x 85 mm.

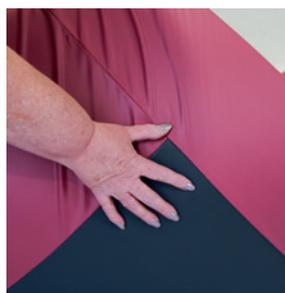

### Covers have an anti-slip base

Ensures safe use of the **Propad Overlay** and helps to prevent slippage.



## Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Size range	Height	Total product weight	Max. user weight	Max. cleaning temperature	
						Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match and BS 7175 (Crib 5)
<b>Propad Mattress Overlay</b>	Std sizes: 1875 x 876 mm <sup>1</sup> 1875 x 1372 mm <sup>2</sup> 1980 x 1524 mm <sup>3</sup>	85 mm	5.5 kg	108 kg	80°	1. Single 2. Double 3. King Size
	Additional sizes available on request					

## Colour

Violet Quartz



\* Essential nursing care is pivotal in pressure ulcer prevention. This product will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

QR  
Code

Tel: +41 61 487 70 70  
hqeurope@invacare.com  
www.invacare.eu.com

Invacare International GmbH  
Benkenstrasse 260  
4108 Witterswil  
Switzerland

