

Aero-tech 4



The Lighter Way to Pressure Care



Clinical Advantage

The *Invacare Aero-tech 4* is a four valve air cushion which is specifically designed for those considered to be at 'Very High Risk' of developing pressure ulcers but require improved stability. Individual interconnected air cells allow free movement of air within the cushion, which allows the level of air to be individually adjusted to body shape thus maximising immersion. This air fluidised movement minimises changes in peak pressures, therefore providing excellent protection against pressure damage.

The natural rubber material and vulcanised bonded base, offers excellent product longevity and may be suitable for those who experience allergic reactions to latex¹.


The *Invacare Aero-tech 4* addresses critical pressure points providing the client with a high level of skin protection. In addition the continual movement of air provides protection for the ischial tuberosities and coccyx area. This cushion provides an infinitely adjustable, low friction and low shear environment which can protect vulnerable clients against pressure damage. The four separate air chambers improve stability, and can be individually adjusted to accommodate mild postural requirements.

Features and Options



Air Valve Adjustment

Incremental air valve adjustment allows the level of air to be adjusted to suit individuals' needs for pressure reduction, stability and comfort.

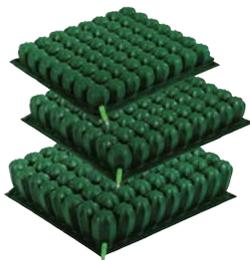

Custom Sizes Available

A range of cushion sizes are available, ensuring a perfect fit for all makes and models of wheelchair. (For size guide, please refer to order form).



Rounded Cell Shape

The rounded shape of the base of the cells ensures good immersion for users.



Optional Cover Plates

Rigid plastic cover plates are available in a range of sizes and can be used to suppress cells.



Bonded Base Layer

The vulcanized bonded base layer of the cushion ensures a robust seal, prolonging longevity.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Width	Length	Height	Min. / Max. user weight	Max. cleaning temperature
Aero-tech 4	380 - 470 mm	330 - 530 mm	100 mm	20 kg / 150 kg	60°

2 year warranty

1. The **Aero-tech 4** is free of latex allergen (ref: BGFA Bochum). Risk assessment must be carried out prior to use.

Colours



QR Code

Tel: +41 61 487 70 70
hqeuropa@invacare.com
www.invacare.eu.com

Invacare International GmbH
 Benkenstrasse 260
 4108 Witterswil
 Switzerland

