

**matrix**<sup>®</sup>  
SEATING SERIES

# Contour Visco<sup>NG</sup>



Conforms to the shape  
of the user's body



Yes, you can.


The *Invacare Matrx Contour Visco<sup>NG</sup>* cushion is an anatomically contoured 100% viscoelastic foam cushion. It provides excellent pressure redistribution properties to improve comfort and stabilise the pelvis. Its low-profile design promotes easy foot propulsion and enables access to lower tables.

### Rear Radius

▶ Contoured radius provides additional support behind the pelvis for optimal pelvic positioning. The cushion radius can fit underneath the backrest upholstery or a contoured rigid back support, filling the gap normally seen when rear squared cut cushions are used. Maximizing the surface area of contact will improve pressure redistribution and skin protection

### Waterfall front edge



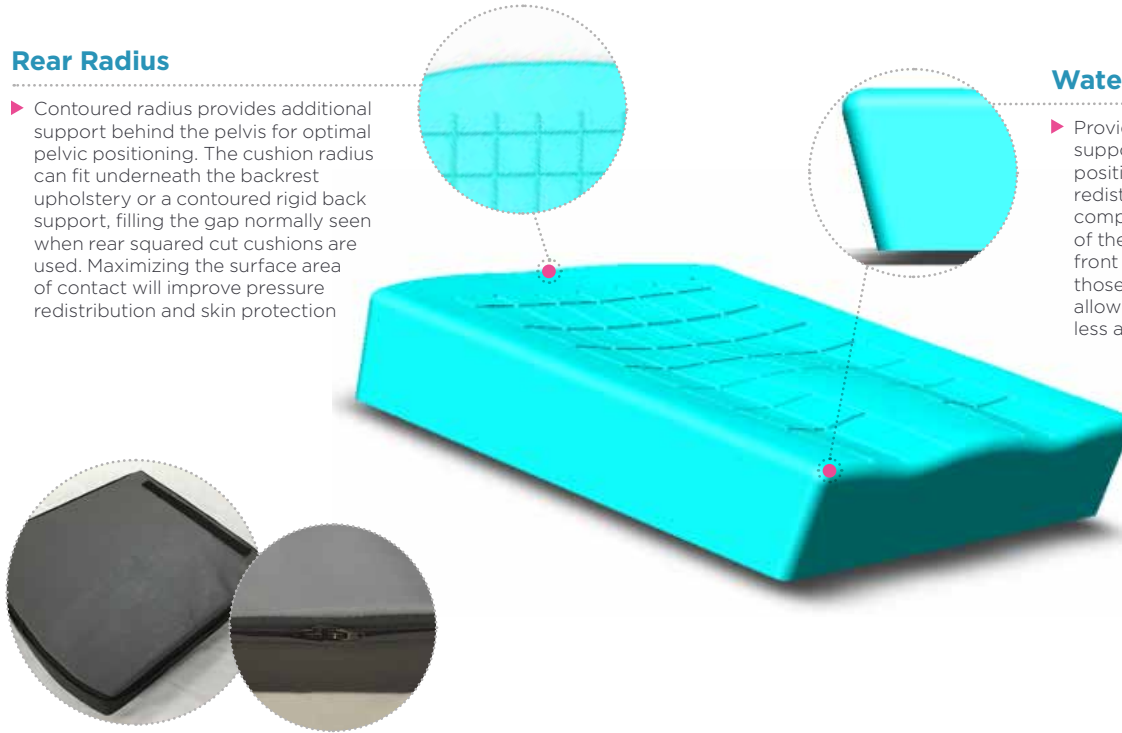





▶ Provides maximum thigh support for improved leg positioning and pressure redistribution, without compromising the popliteal area of the knee joint. The waterfall front edge also accommodates those with tight hamstrings to allow foot positioning at 90 or less and enables foot propulsion



▶ The *Matrx Contour Visco<sup>NG</sup>* provides superior immersion and envelopment. When sitting on the cushion the user will feel the viscoelastic foam slowly adapting to their body shape, increasing the surface of contact, and thus redistributing pressure. The *Matrx Contour Visco<sup>NG</sup>* is available in a wide range of sizes and weight capacities, providing the best fit to the user. It includes 2 covers with 2 zippers, easy to remove for cleaning and maintenance. The covers are incontinence proof and bi-elastic to allow efficient immersion of the user body in the foam cushion. The non-slip base ensures that the cushion remains in place and the carry-handle allows easy transfer of the cushion from one surface to another.

## Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

								
	Width ▼	Depth ▼	Height ▼	Height ▼	Total product weight ▼	Max. user weight ▼	Max. cleaning temperature ▼	
<b>Matrx Contour Visco<sup>NG</sup></b>	355 mm 380 mm 405 mm 430 mm 455 mm 480 mm 505 mm	405 mm 430 mm 455 mm 530 mm	A=100mm B=100mm	100mm	1.3kg	140Kg*	90°	

Cushion is manufactured to EN1021-1 Cig and EN1021-2 Match

\*Cushions with seat width 380 and 405 mm - minimum user weight is 40 kg  
Cushions with seat width 430 mm and above - minimum user weight is 60 kg

**Tel: +44 61 487 70 70**  
hqeurope@invacare.com  
www.invacare.eu.com

Invacare International GmbH  
Neuhofweg 51  
4147 Aesch BL  
Switzerland



**Yes, you can.®**