



Invacare® | **Flo-tech PT Back & Seating System**

A versatile system promoting good posture and exceptional comfort



## Invacare **Flo-tech PT Back & Seating System**

The *Invacare Flo-tech PT* System enables a wide range of posterior, anterior, lateral and angle adjustment, accommodating growth and individual user needs. The system can be applied to a wide range of wheelchair frames and can be used in conjunction with the **Flo-tech PT** Drop Seat Board and other accessories to create a fully adjustable seating system.

Made from **lightweight components** for easier propulsion



### Product Summary

**Flo-tech PT** features a modular seating system which can be configured and adjusted specific to individual needs to maximise user posture, comfort and support.

- Quick release back rest and drop seat board to allow the wheelchair frame to be folded for transportation and storage.
- Multi-adjustable system enabling adjustment of seat depth and angle, back rest angle, and the position of the back rest in relation to the axle.
- Wide range of sizes to fit wheelchair widths from 300 mm to 610 mm
- Crash tested to ISO16840-4
- Mounting brackets are constructed from high grade aluminium which is strong and lightweight
- Back rest available in Standard and Tall both of which can be upgraded to a Deep Contour Back
- New Spacer Fabric now comes as standard on back shell cushions for improved breath ability. This cover now also features a discreet back hardware concealing jacket
- The **Flo-tech PT Back & Seating System** standard package includes the back shell, castellated foam cushion with new spacer cover and your choice of Lock and Latch or Easy Clamp mounting system.



### Features and options



#### Foot Supports

Shoe Holders and Anklesures (High-Tops) are designed to optimise foot positioning and support.



### Drop Seat Boards

A fully adjustable seat board to provide a firm base of support for a wheelchair cushion.



### Headrests

A variety of multi-adjustable headrests and pads are available to provide additional comfort and support.



### Lateral Supports

Fixed and swing-away lateral brackets available with a choice of flat and curved pads.



### Knee Adductors

Pad and bracket combination which can be positioned at the side of the knee to prevent excessive abduction.



### Harnesses

A wide range of harnesses can be fitted to the Flo-tech PT back rest and wheelchair frame for additional anterior support and posture management.



### Pelvic Belts

Available in a variety of styles and fastening types, the pelvic belt range offers both postural support and pelvic stability.



### Abductors

Used to position the user's legs out of adduction for improved posture and stability in order to maintain a midline position.

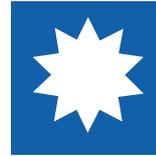


### Hip/Thigh Guides

Facilitate thigh alignment, hip abduction and helps to maintain central pelvic placement.

## Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Frame size	Back shell weight	Max. User weight	Max. cleaning temperature	Warranty
					
<b>Flo-tech PT Back &amp; Seating System</b>	300 - 610 mm	2.1 kg	Shell size 12 -20: 113 kg Shell size 22 -24: 190 kg	80°	3 years on all hardware and individual parts

Crash Testing Passed ISO 16840 - 4 and ISO 7176 - 19  
(in an **Action<sup>3</sup>** Wheelchair)

## Colour

Black spacer fabric



Please notice colors may vary slightly from those displayed above.

QR Code

Tel: +41 61 487 70 70  
[hqeuropa@invacare.com](mailto:hqeuropa@invacare.com)  
[www.invacare.eu.com](http://www.invacare.eu.com)

Invacare International GmbH  
Benkenstrasse 260  
4108 Witterswil  
Switzerland

Trusted Solutions,  
Passionate People

