

**AQUATEC®**

# **Ocean VIP Ergo**

## **Ocean Dual VIP Ergo**



**Modular personal care solutions  
for complex needs**

ERGONOMIC DESIGN | SAFETY | INDEPENDENCE



# Modular personal care solutions for those with more complex seating and positioning needs

The Aquatec Ocean VIP Ergo and Ocean Dual VIP Ergo are modular tilt in space shower chair commodes designed to suit a wide range of individuals and their postural needs. Both models benefit from a -5° to 40° seat tilt angle, whilst the Ocean Dual VIP has an additional backrest recline of 0° to 35° for those with more complex seating and positioning requirements.


## Highly customisable

As with all Ocean Ergo models, the VIP and Dual VIP can be tailored to each individual by making some simple adjustments or by adding some easy to fit, off the shelf accessories. Furthermore, their ergonomic design ensures ease of use for the carer whilst providing optimum safety and comfort to the user.

See full range of accessories online  
[www.invacare.eu.com](http://www.invacare.eu.com)


## NEW Contoured, ergonomic seat

- ▶ Based on the imprints people leave behind when seated to ensure a more upright and stable seating position.
- ▶ A 5° seat tilt to lift the knees into the natural squat position for toileting.
- ▶ Transfer supports at the front to allow the user to position themselves into a comfortable position for toileting and showering.

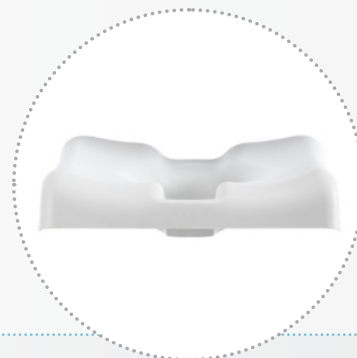

## NEW Ergonomic tilt

- ▶ The new pivot and slide mechanism means that during tilt, the chair now moves downwards and forwards which keeps the user's centre of gravity close to that of the chair. This makes the tilt operation lighter and more intuitive for the carer, creates a smoother and more comfortable experience for the user and has reduced the overall footprint of the product without affecting its stability.

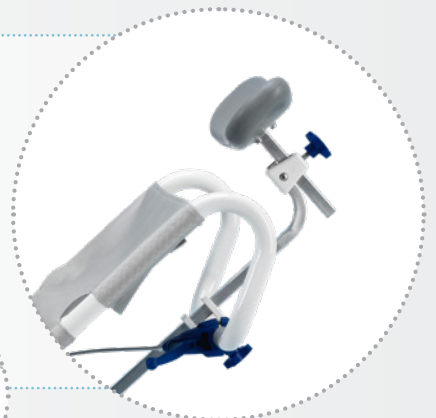

**NEW**

*FEATURES*

**-5° to 40°**

**1**

### **Extremely stable**

- ▶ Enhanced stability for additional safety and reassurance

**2**

### **Ergonomic tilt**

- ▶ Ergonomic tilt with weight shifting abilities

**3**


### **Innovative ergonomic seat**

- ▶ Contoured, ergonomic seat plate encourages a more upright seating position for enhanced independence and comfort

**4**

### **Smaller footprint**

- ▶ Optimised wheel base makes it easier to manoeuvre in small spaces

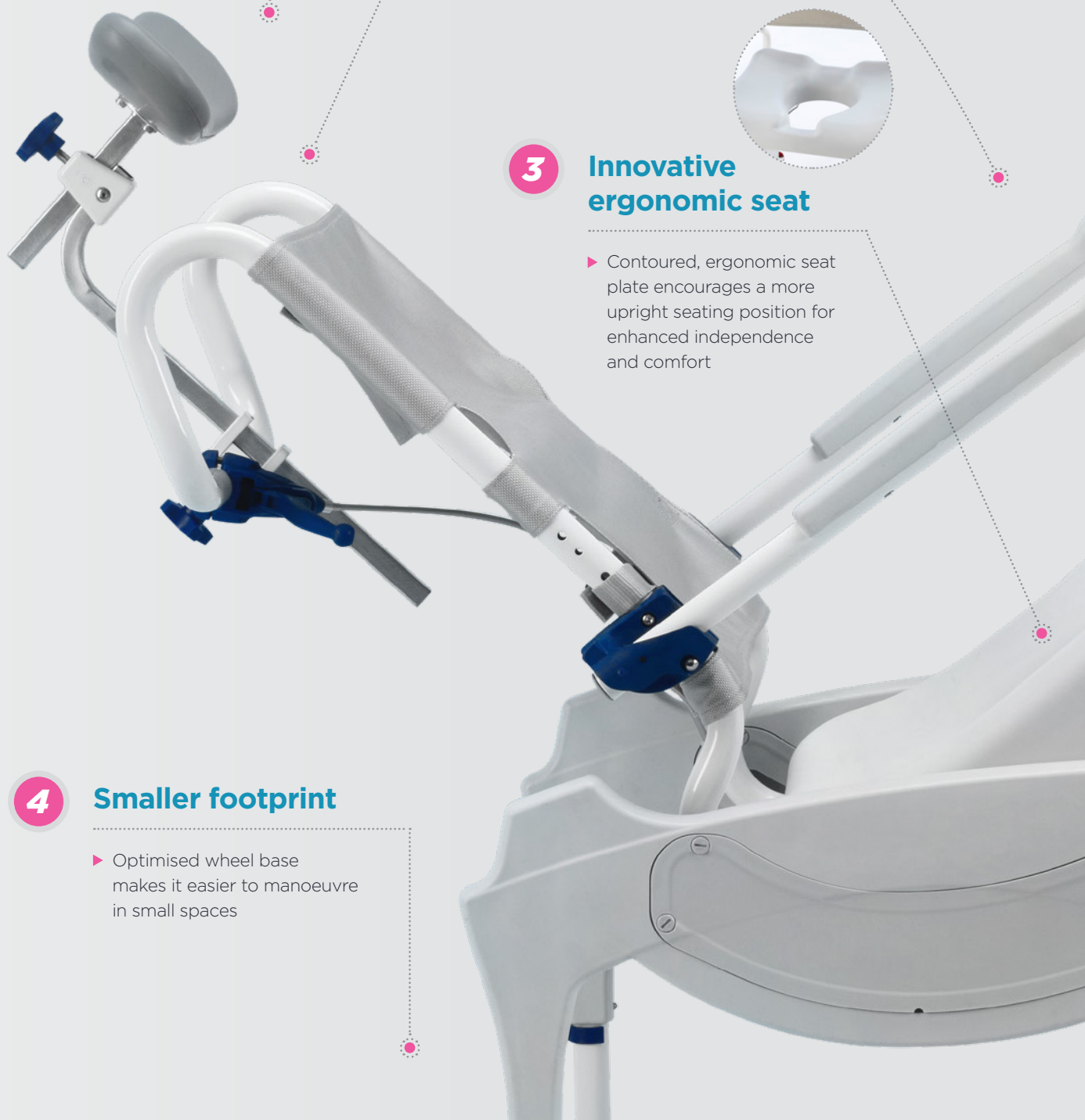

## Other features



▲ Stainless steel corrosion free frame.



▲ Extremely stable even at maximum tilt.




▲ Machine washable tension adjustable backrest.



▲ Armrest height, footrest height and seat height can be easily adjusted.



▲ 5° - 35° Backrest Recline  
**Dual VIP Ergo Only**

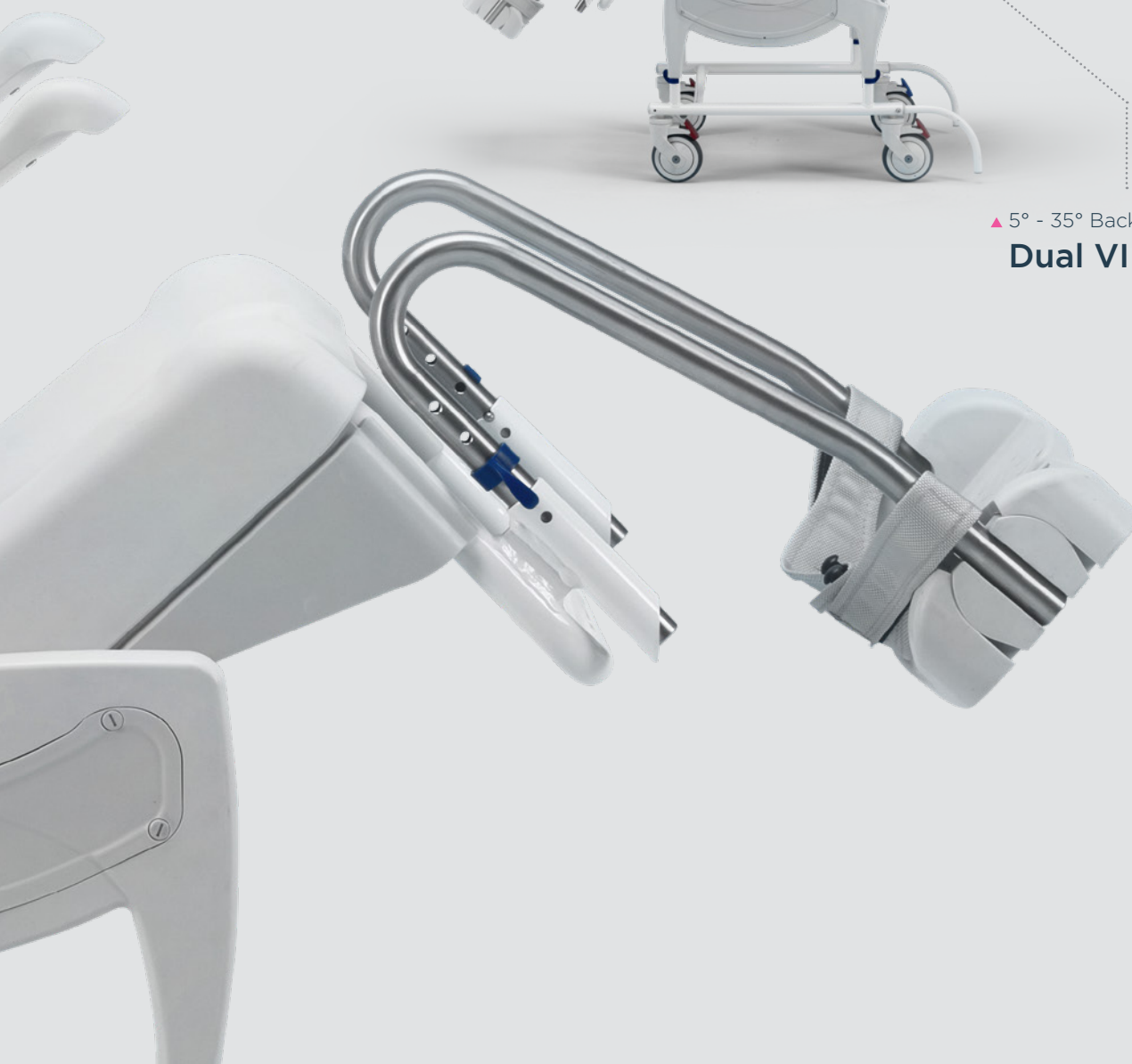



# Meet the Ocean VIP Ergo and Ocean Dual VIP Ergo

---

Modular tilt in space models to suit a wide range of individuals and their postural needs.



**Ocean VIP Ergo**

▲ -5°-40° seat tilt angle



**Ocean Dual VIP Ergo**

▲ -5°-40° seat tilt angle with an additional  
backrest recline of 5° - 35°

## Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.



Seat width ▼



Seat depth ▼



Total width ▼



Total depth ▼



Width between armrests ▼

<b>Ocean VIP Ergo</b>	480 mm	480 mm	565 mm	1000 mm	455 mm
<b>Ocean Dual VIP Ergo</b>	480 mm	480 mm	565 mm	1000 mm	455 mm



Tilt angle seat ▼



Tilt angle back ▼




Weight ▼



Load capacity ▼

<b>Ocean VIP Ergo</b>	-5° - 40°	-	23.5 kg	150 kg
<b>Ocean Dual VIP Ergo</b>	-5° - 40°	5° - 35°	28.3 kg	150 kg

## Frame colour ▼

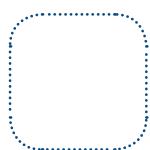

White

## Regulatory data ▼



**Tel: +41 61 487 70 70**  
[hqeuropa@invacare.com](mailto:hqeuropa@invacare.com)  
[www.invacare.eu.com](http://www.invacare.eu.com)

Invacare International GmbH  
 Benkenstrasse 260  
 4108 Witterswil  
 Switzerland



**Yes, you can.®**