

**Spin<sup>x</sup>**



**Are you sitting comfortably?**





The *Invacare Spin<sup>x</sup>* is designed to fit the user living an active social life and in need of body support and mobility. With an inbuilt postural frame and fully adjustable seat, every user gets individual fitting that supports them in their daily activities. Made with lightweight materials it has a dual folding mechanism making it suitable for storage and transportation.

### Great Seating

The *Invacare Spin<sup>x</sup>* allows erect posture and 90° angle legrests. The rear of the seat rail is specifically shaped to provide the appropriate additional space required by the pelvis, resulting in an ideal base for a stable and relaxed seating position.



### Attractive weight and easy to transport

The *Invacare Spin<sup>x</sup>*, thanks to its lightweight frame, is the key to its easy lifting. It weighs less than 13 kg in its minimal configuration and 7.5 kg when removable parts are detached.

Moreover, the dual folding mechanism, where the first stage is to fold the crossbars in the traditional way and then fold over the back posts at the base of the seat, makes for an incredibly compact and protected package, ideal for handling and stowage.

### Quality Proven Components

Many of the components used on the *Spin<sup>x</sup>*, such as armrests, legrests, and push handles, are currently being used on other wheelchairs in the Invacare product range. These items have a proven track record of quality and durability and are readily available and cost effective if refurbishment is required.

## Features and options



### Erect posture

The rear of the seat tube has a 5° angle providing sufficient room for the pelvis.



### Depth adjustable seat upholstery

For individual fitting and perfect support the front part of the seat upholstery is depth adjustable over 50 mm.



### Folding system

The single crossbar is integrated in the frame. The two largest seat widths have a double crossbar.



### External push handles

Height can be adjusted as desired over 100 mm by loosening the levers.

## Features and options



### One piece footplate

With the one piece footplate, increased stability as well as large surface for the position of the feet is achieved.



### Legrests assortment

① Elevating legrest, ② fixed 80° footrest and ③ fixed 90° footrest are available.



### Frame protection

With the frame protection scratches are avoided on the frame when transporting the wheelchair.



### Rear wheel support

Gives several balance and height options. It can also be turned to achieve a 3° or 0° camber.



### Backrest

The backrest is height adjustable and is angled to follow the body shape in the same way as the seat. An angle adjustable folding version is available as an option.



### Anti-tipper

Can easily be folded away by the foot. The device is also angled outwards, to give space for the feet for the attendant.



### Passive rearwheel fixation

For increased stability a passive rearwheel support can be offered as an option.






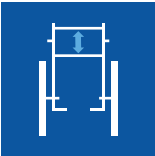















### Armrests

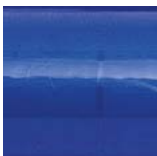

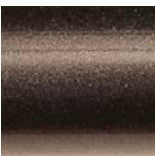



① Sideguard, ② mudguard, ③ flip up height adjustable long or short and ④ Rio height adjustable long or short are available.

## Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

Invacare <b>Spin<sup>x</sup></b>	<b>Seat width</b>  355* / 380 / 405 / 430 / 455 / 480 / 505 mm	<b>Seat depth</b>  400 / 450 / 500 mm (-50 mm)	<b>Seat angle</b>  0° / 3° / 6°	<b>Seat height</b>  410 / 430 / 455 / 480 / 505 / 530 mm	<b>Backrest angle</b>  ± 12°	<b>Backrest height</b>  350 - 410 / 360 - 400 / 410 - 490 mm	<b>Front frame angle</b>  80° / 90° / 0° - 80°	
Invacare <b>Spin<sup>x</sup></b>	<b>Length of footrests</b>  320 - 440 / 410 - 530 mm	<b>Armrest height</b>  190 - 290 mm	<b>Driving unit width</b>  Seat width + 210 mm	<b>Total width, folded</b>  290 mm	<b>Total length inc. Footrests</b>  800 - 1260 mm	<b>Total height</b>  760 - 1130 mm	<b>Max. User weight</b>  125 / 135** kg	
Invacare <b>Spin<sup>x</sup></b>	<b>Total weight</b>  13 kg	<b>Mass of the heaviest part</b>  7,5 kg	<div><div>COMPATIBLE</div><div></div></div>				<div><div>CRASH TESTED</div><div></div></div>	<div><div></div><div>1. Available as special 2. For seat widths 480 and 505 mm only</div></div>

## Frame colours

Blue	Red	Black	Grey	Lime green <sup>1</sup>	Orange <sup>1</sup>
					

## Upholstery

Black Nylon	Black Jemima
	

Please notice colors may vary slightly from those displayed above.

Tel: +41 61 487 70 70  
 hqeuropa@invacare.com  
 www.invacare.eu.com

Invacare International GmbH  
 Benkenstrasse 260  
 4108 Witterswil  
 Switzerland

Scan here to access the website.



© 2018 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. SpinX - EU - 10/2018