

ScanBeta[®] NG



JUNIOR

A safe, configurable bed that
changes with your child



Invacare **ScanBeta_{NG}** paediatric beds are available in three sizes and are perfect for children aged from 3 to 12 years old. **ScanBeta_{NG}** conforms to the very highest degrees of safety requirements, meeting the Nordic Requirements Specifications: there are no distances above 6 cm on the sides or at the mattress support. All beds in the range come complete with a handset that incorporates a lockable function to avoid unintentional use. Both the functionality of the mattress platform and the sides can be configured to meet individual requirements and can be up- or downgraded at a later stage, in line with a user's changing needs.



Stringent specifications

The Invacare **ScanBeta** in a standard configuration is rigorously tested against Nordic Requirement Specifications. Designed with safety as a priority, the bed meets all the current regulations and safety criteria.



Configuration

The mattress platforms for the three standard sizes have been constructed according to anthropometric analyses in order to maximize the child's comfort.

According to needs, the mattress platform can be configured to be either electrically or manually operated.

There is a variety of side options (fixed sides, opening panels etc.) in line with individual requirements. Upgrades are also available later on, as the child's needs change.

Ergonomic

The height is adjustable from 300 - 900 mm or 320 - 920 mm (depending on the size of the bed). When required it is very easy for a child to get in and out of bed. For caregivers, all the necessary features and functions for safety and comfort, such as central brake (option), height adjustability and lockable hand control are all in place.

Features and options



The optimized lengths of the sections are based on anthropometrical measurements of children. The 3 sizes of **ScanBeta_{NG}** Standard have been developed specifically for children aged between 3 and 12 years old.



A slim and ergonomic handset comes with an individual lockable function for increased safety.



Snap lock for quick and easy removal of the mattress support from the bed frame.

Features and options



ScanBeta_{NG} configured with 400mm lowering side rails

Bars, attractive beech finish, double-catch release mechanism for increased safety.



ScanBeta_{NG} configured with 400mm lowering side rails

Perspex, double-catch release mechanism for increased safety.



ScanBeta_{NG} configured with 400 mm panels

Bars, attractive beech finish, panels open to 180° and secure locking mechanism with catch release for increased safety.



ScanBeta_{NG} configured with 800 mm folding panels

Bars, attractive beech finish, panels open to 180°, minimizing space requirements.



Configuration of the mattress platform

Functionality of the mattress platform can be increased or decreased by adding or removing actuator(s) and rastofix. (See table on back page).



Support handle

Makes getting in and out of bed much easier.



Lifting Pole

The handle can be adjusted both in height and depth.



Upholstery

20 mm thick upholstery for the bed end and all versions of sides/side rails.



A safe locking mechanism which can only be opened via two different catches.



Stickers

Choose between different stickers for decoration of the bed ends – interior and exterior.



Possible configurations for the *Invacare ScanBeta_{NG}* mattress support

Fixed	
Motor on back rest and fixed leg section	
Motor on back rest and manual operated leg rest	
Motor on back rest and thigh section. Manual operated leg section	

* For each of the 3 sizes it is the same mattress platform being configured with actuator(s) and rastofix.

Sides and panels for *Invacare ScanBeta_{NG}*


	Fixed side	Lowering side rail	Panels	Foldable panels
400 mm height	✓	✓	✓	✓
800 mm height	✓		✓	✓

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local **Invacare website**.

	Width outside Width inside	Length outside Length inside	Height adjustment	Mattress support dimensions	Backrest angle	Thigh angle	Legrest lift
600 x 1400 mm	680 mm 620 mm	1480 mm 1420 mm	320 - 920 mm	620 - 70 - 210 - 460 mm	0 - 73°	0 - 36°	0 - 26°
700 x 1600 mm	780 mm 720 mm	1680 mm 1620 mm	300 - 900 mm	670 - 150 - 220 - 520 mm	0 - 73°	0 - 31°	0 - 24°
800 x 1800 mm	880 mm 820 mm	1880 mm 1820 mm	300 - 900 mm	750 - 190 - 270 - 550 mm	0 - 73°	0 - 27°	0 - 26°

Max. user weight

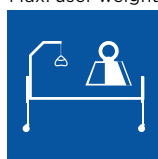

All sizes

SWL 100 kg
70 kg

Bespoke-build service



The *Invacare ScanBeta* can also be supplied in other dimensions. Please contact us if you would like more information about custom-made beds.

Note that *Invacare ScanBeta* beds made to custom requirements are not tested against Nordic Requirement Specifications.

Frame colour and wood work

Pearl white
(metal)

Beech (wood)

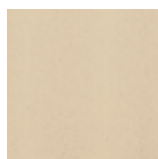

Please notice colors may vary slightly from those displayed above.

QR
Code

Tel: +41 61 487 70 70
hqeurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Benkenstrasse 260
4108 Witterswil
Switzerland



© 2018 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. ScanBetaNG - EU - 10/2018