



FOR THE ACTIVE WHEELCHAIR USER

**Ischial well - Two compartments**  
for superior pressure reduction in the IT-coccyx area & extra positioning possibilities

**Trochanter support**  
firmer density to provide lateral stability | redistributes pressure to help create an ischial well

**Ergonomically-shaped foam**  
for easy transfers

**SmartCell – Foam barrier**  
limits sliding forward

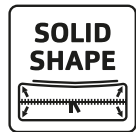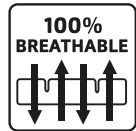

- ▲ Perfect for users who transfer frequently
- ▲ Excellent lateral and frontal stability
- ▲ Features a nifty storage pouch for personal items



The Vicair Active 02 is a pre-contoured cushion which consists of 4 individual compartments filled with Smartcells at the back of the cushion and an anatomically contoured foam inlay at the front. Striking a unique balance between preserving skin at high risk areas at the ischials and coccyx by redistribution of pressures, stabilizing the pelvis and facilitating easy transfers. For specific needs the filling of each of the 4 SmartCell filled compartments can be easily adjusted. Allowing you to change the density of each compartment benefitting users who: (1) Need correction or accommodation of an assymetric pelvis. (2) Want (partial) off-loading of the ischial and coccyx area. (3) Benefit from additional lateral support.

**PRODUCT PERFORMANCE**

**Skin protection / Anti decubitus**

|                            |       |
|----------------------------|-------|
| Pressure redistribution    | ▲▲▲▲△ |
| Shear force reduction      | ▲▲▲▲△ |
| Microclimate control       | ▲▲▲▲△ |
| Impact & vibration damping | ▲▲▲▲△ |

**Positioning**

|                           |       |
|---------------------------|-------|
| Lateral stability         | ▲▲▲▲▲ |
| Frontal stability         | ▲▲▲▲▲ |
| Sliding reduction         | ▲▲▲▲▲ |
| Allows manual adjustments | ▲▲▲▲△ |
| Increasing endurance      | ▲▲▲▲▲ |

**Maintenance & Safety**

|                            |        |
|----------------------------|--------|
| Durability                 | ▲▲▲▲▲  |
| Reliability                | ▲▲▲▲▲  |
| Functionality when misused | ▲▲▲△△△ |
| Easy to maintain           | ▲▲▲▲△  |
| Level of hygiene           | ▲▲▲▲△  |

**Functionality & Application**


|                          |       |
|--------------------------|-------|
| Comfort                  | ▲▲▲▲▲ |
| Ultra light weight       | ▲▲▲▲△ |
| Allows easy transfers    | ▲▲▲▲▲ |
| Available standard sizes | ▲▲▲▲▲ |
| Maximum user weight      | ▲▲▲▲△ |

**TECHNICAL SPECIFICATIONS**

|                  |                     |
|------------------|---------------------|
| Thickness        | 6 and 9 cm          |
| Average weight   | 1250 grams          |
| Max. user weight | 200 kg              |
| Machine washable | 60°C (without foam) |

**CUSHION SIZES**

- Standard (Sizes >53 cm thickness 9 only)
- On request

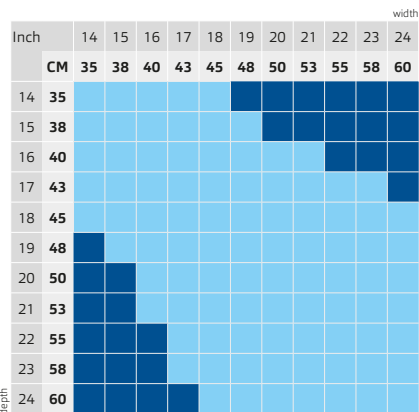

**AVAILABLE COVERS**

Not applicable. Vicair Active 02 only has one cover serving as inner and outer cover.