

ESPRIT *Action*[®]

DISCOVER. EXPLORE. EXPERIENCE.



NEW foldable & lightweight power chair

Featuring innovative gyro technology




Yes, you can.[®]

Advanced new design.

Our Esprit Action® has moved forward with technical precision.
Every part of its new design has made progressive steps.

An extensive customer-led design programme has made the difference. The new Esprit Action® is lightweight, foldable and uses integrated Alber powered wheels on a highly customisable

aluminium chassis from our popular Action® range. The results are unprecedented. It's finally time to enjoy the freedom and true spirit of the Esprit Action®.

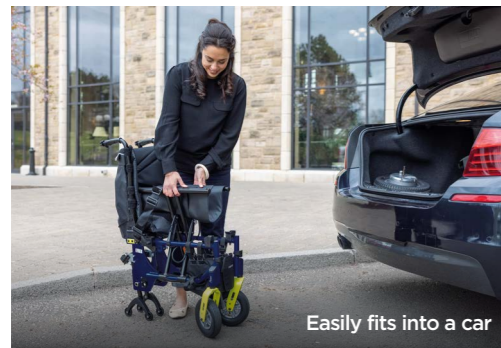

Easily fits into a car



Handles incline up to 6°



Foldable with removable parts

KEY FEATURES



Negotiates narrower spaces

Featuring new innovative gyro technology, the compact wheelbase gives enhanced manoeuvrability making the narrowest of environments easier to negotiate. The rear wheels are repositioned forward for better traction and a smaller turning radius.



Effortlessly handles inclines

Inclines of up to 6° are effortless, smooth and stable. The new gyro technology reduces acceleration according to the degree of the slope and quickly detects when the wheelchair is on a tilt.



Easier to transport

Designed using an Action® lightweight folding frame and removable powered rear wheels and battery, the Esprit Action® can be folded for transportation without ramps and adapted vehicles.



Cross compatibility

Cross-compatibility with the Action® range allows for quick adaptations. The Esprit Action® can be customised around specific needs quickly and cost effectively.

15kg*

*min. transportable weight

8km/h

speed option

14" Alber powered wheels



Better drive performance

14" rear wheels and pneumatic rear and front tyres give a smooth and comfortable drive. The new design also features a new 8km/h speed option for improved outdoor performances.

Options



Lightweight Alber swivel arm



Premium anti-tipper option provides jack-up function for quick release of wheels



Comfort & support

Optimal positioning is achieved through highly configurable options and accessories. The wheelchair can be discreetly fine-tuned to meet clinical needs.



Discover.

A new design

Explore.

New possibilities

Experience.

Even better performance

Introducing the next generation of Esprit Action®. An exciting new design using an innovative combination of technology, performance and carefully considered refinements.





ESPRIT Action®

Junior

The Esprit Action® Junior enables children to enjoy an active and independent lifestyle. The new Junior encourages kids to be kids, be happy, active and ooze confidence.

Whether it's connecting with friends online indoors or meeting up and playing outdoors, every child can actively participate and feel good about their Esprit Action®. Sharper visual aesthetics and a choice of three colour combinations, the Esprit Action® ensures children use their wheelchair without limitations. More possibilities, more fun.

Technical data

For more information about this product, including the user manual, please visit our website: www.invacare.eu

	Seat width	Seat depth	Seat height	Backrest height	Length of footrests	Seat angle	Backrest angle
Esprit Action®	380 - 605 mm	400 - 500 mm	460 or 485 mm	390 - 510 mm	340 - 470 mm	3°	0° - 30°
Esprit Action® Junior	305 - 380 mm	325 - 375 mm	460 or 485 mm	390 - 510 mm	200 - 330 mm	3°	0° - 30°
	Driving unit width	Total height	Total length inc. footrests	Total length without footrests	Total weight	Max. user weight	Battery capacity
Esprit Action®	Seat width + 180 mm	945 - 1000 mm	980 - 1115 mm	640 - 970 mm	40.4 - 45 kg ¹	125 kg	36 V (3 x 12V)/ 7.65 Ah
Esprit Action® Junior	Seat width + 180 mm	945 - 1000 mm	925 - 1075 mm	640 - 970 mm	40.4 - 45 kg ¹	80kg	36 V (3 x 12V)/ 7.65 Ah
	Turning diameter	Max. climbable obstacle height	Max. safe slope	Motor capacity	Driving range	Speed	Adjustable armrests
Esprit Action®	1430 - 1960 mm	30 mm	6°	6 km/h: 2 x 110 W 8 km/h: 2 x 150 W	15km	6 km/h or 8 km/h	+50 mm
Esprit Action® Junior	1430 - 1960 mm	30 mm	6°	6 km/h: 2 x 110 W 8 km/h: 2 x 150 W	15km	6 km/h or 8 km/h	+50 mm

¹Weight of heaviest part is the frame at 15kg

Frame colours

Adult



Slate Grey



Navy Blue

Adult and Junior



Anthracite Black



Dark Purple

Junior



Water Blue

Colour accents



Lime



Black

Tel: +41 61 487 70 70
hqeurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Neuhofweg 51
4147 Aesch BL
Switzerland



Yes, you can.®