




THE MOST STABLE CUSHION FOR WHEELCHAIR USERS

**Ischial well - Two compartments**  
for superior pressure reduction in the IT-coccyx area & extra positioning possibilities

**Trochanter support**  
firmer density to provide lateral stability | redistributes pressure to help create an ischial well

**Front compartments** provide leg alignment and stability | easily adjustable for individual positioning needs

**Mid height zip closure** for easy adjustment when client is seated | cushion keeps its shape, enhancing the user's skin protection

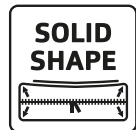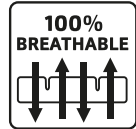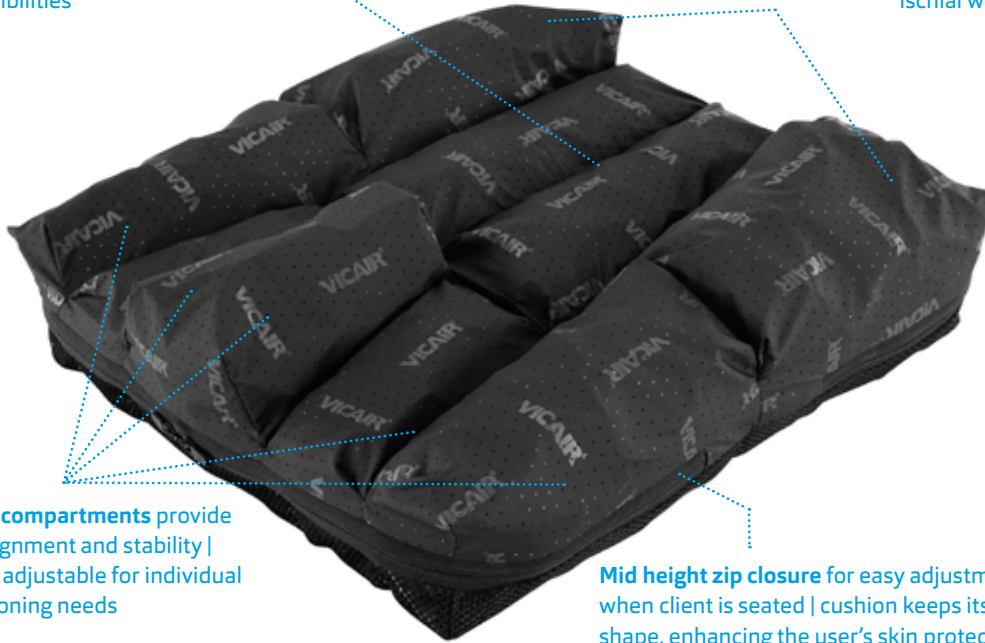

- ▲ **Nine adjustable compartments for optimal stability**
- ▲ **Expanding the user's reach**
- ▲ **Extra positioning of the thighs**

The Vicair Vector 02 is a pre-contoured cushion which consists of 9 individual compartments. Each compartment is designed to influence a specific area of the pelvis, striking a unique balance between preserving skin at the ischials and coccyx by re-distributing pressures, stabilizing the pelvis and positioning of the femurs. For specific needs the filling of each of the 9 compartments can be easily adjusted. Allowing you to change the density of each compartment benefitting users who: (1) Need correction or accommodation of an asymmetric pelvis. (2) Want (partial) off-loading of the ischial and coccyx area. (3) Need correction or accommodation for better thigh alignment. (e.g. abduction, one side foot propulsion, opening pelvic thigh angle).

**AVAILABLE COVERS**  
Vicair O2 Top Cover, Comfair Cover, Incotec Cover and Reversible Cover.

**PRODUCT PERFORMANCE**

**Skin protection / Anti decubitus**

Pressure redistribution	▲▲▲▲▲
Shear force reduction	▲▲▲▲▲
Microclimate control	▲▲▲▲▲
Impact & vibration damping	▲▲▲▲▲

**Positioning**

Lateral stability	▲▲▲▲▲
Frontal stability	▲▲▲▲▲
Sliding reduction	▲▲▲▲▲
Allows manual adjustments	▲▲▲▲▲
Increasing endurance	▲▲▲▲▲

**Maintenance & Safety**

Durability	▲▲▲▲▲
Reliability	▲▲▲▲▲
Functionality when misused	▲▲▲▲△
Easy to maintain	▲▲▲▲▲
Level of hygiene	▲▲▲▲▲

**Functionality & Application**

Comfort	▲▲▲▲▲
Ultra light weight	▲▲▲▲▲
Allows easy transfers	▲▲▲▲△
Available standard sizes	▲▲▲▲▲
Maximum user weight	▲▲▲▲△

**TECHNICAL SPECIFICATIONS**

Thickness	6 and 10 cm
Average weight	725 grams
Max. user weight	250 kg
Machine washable	60°C max.

**CUSHION SIZES**

- Standard (Sizes >53 cm thickness 10 only)
- On request

