

# Dacapo Ampla

Supporting plus-size patients



**NEW high specification foam mattress**

For effective pressure redistribution



**Yes, you can.®**

# An overview of challenges for plus-size patients


Worldwide obesity continues to increase in both adults and children. In most EU states, over half of adults according to Eurostat\* are declared overweight. Obesity causes a range of everyday challenges from self-care to limited mobility and leaves patients exposed to many adverse health problems.

The importance of access to adequately designed equipment that helps patients live a more active and independent lifestyle, has a direct impact on both mental health and clinical outcomes.

The need for comprehensive plus-size nursing care is increasingly required for global healthcare environments. Maintaining high levels of comfort and support to patients in a dignified way should inclusive to all.

**Our new Dacapo Ampla pressure redistributing mattress is the latest addition to our range of products specifically designed for plus-size patients.**

\* Source: [https://ec.europa.eu/eurostat/statistics-explained/index.php?title=Main\\_Page](https://ec.europa.eu/eurostat/statistics-explained/index.php?title=Main_Page)

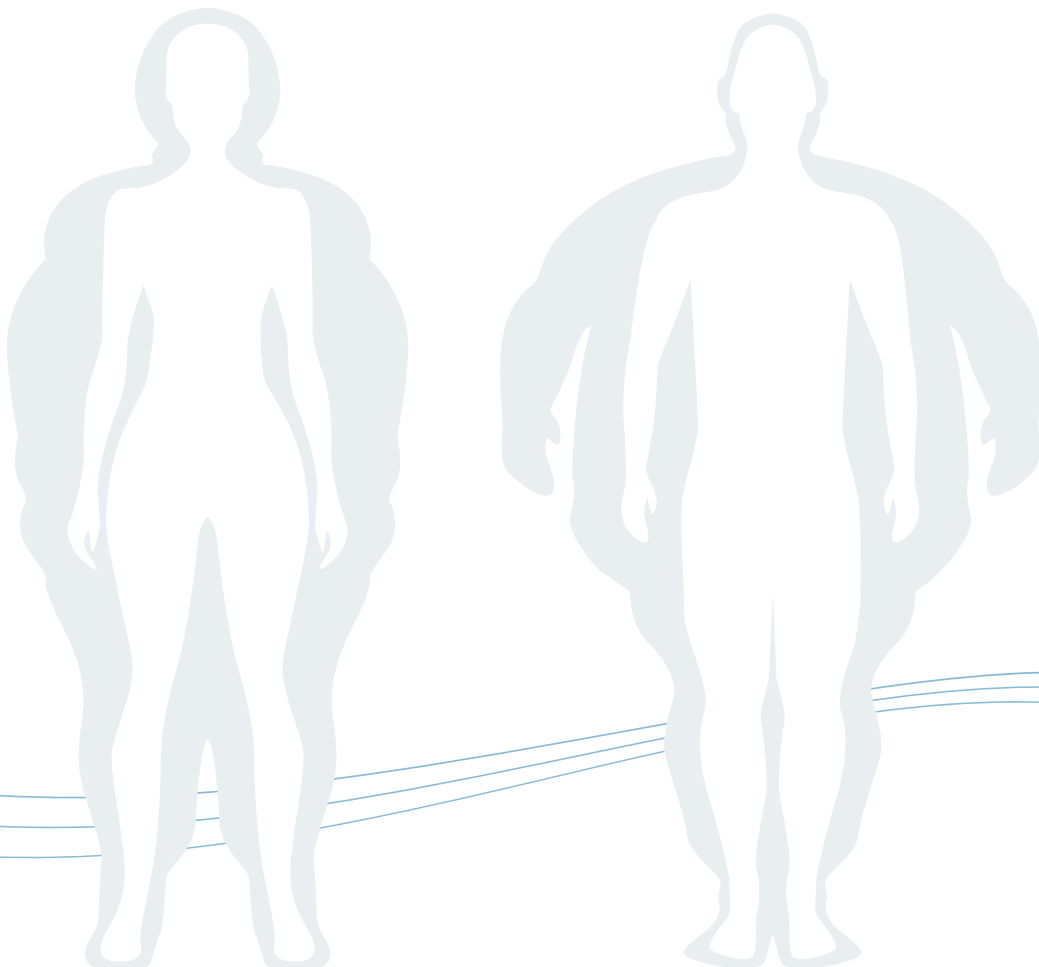

# Designed for plus-sized patients

Our new Dacapo Ampla high specification mattress provides effective pressure redistribution and microclimate management for plus-size patients.


The design focuses on achieving excellent clinical outcomes required by healthcare professionals including optimal immersion and envelopment, pressure redistribution and the reduction of shear and friction forces.

Maximum comfort in every position to support patient repositioning schedules and care giving routines such as hygiene and personal care.



**Delivered compressed and rolled** allowing 3x more mattresses to be stored and transported for a reduced carbon footprint

1

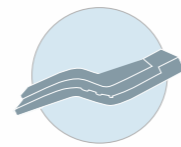

## Clinically effective

Presents outstanding pressure redistribution properties for effective pressure injury prevention.

2



## Demographically representative


Accommodates the growing need for larger, resilient and clinically effective products for plus-size patients.

3



## Quality assured

Delivering optimal quality. Our testing includes pressure mapping, independent durability testing, cytotoxicity testing, fire testing and accelerated aging.

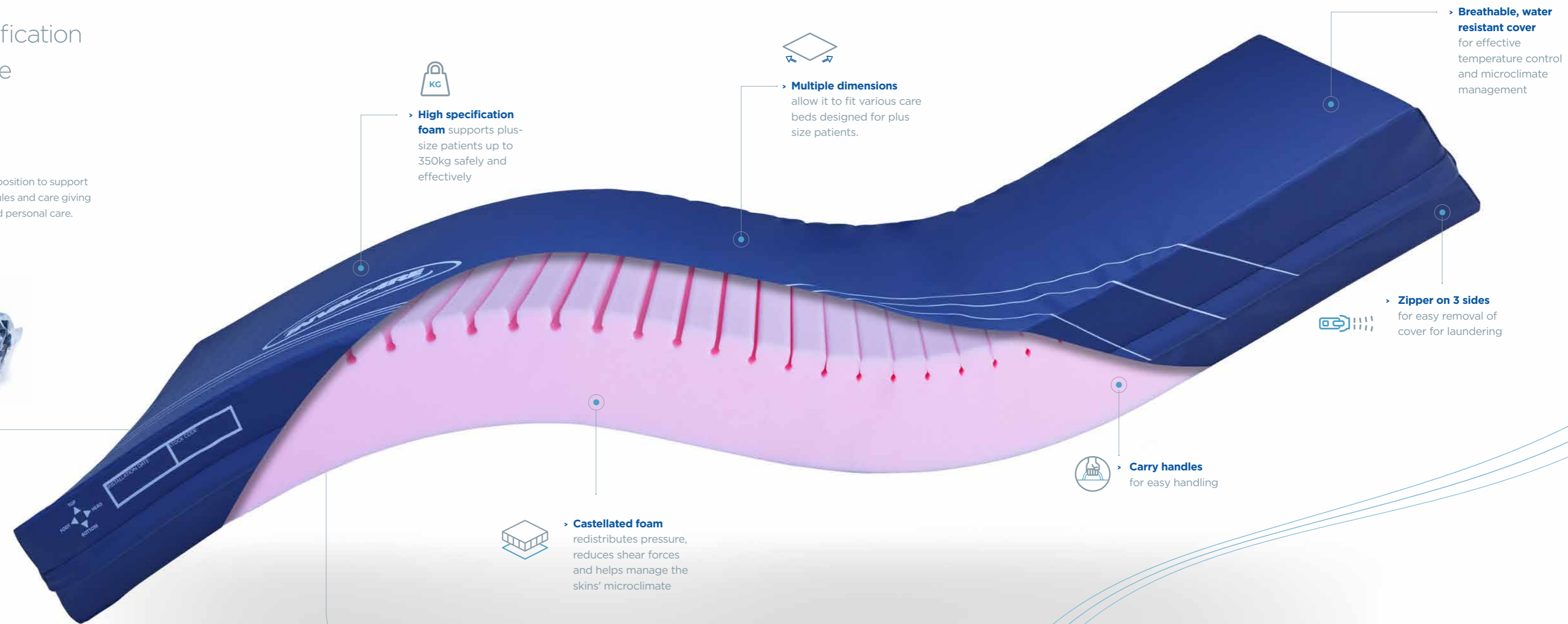

**High specification foam** supports plus-size patients up to 350kg safely and effectively



**Multiple dimensions** allow it to fit various care beds designed for plus size patients.



**Breathable, water resistant cover** for effective temperature control and microclimate management



**Zipper on 3 sides** for easy removal of cover for laundering



**Carry handles** for easy handling




**Castellated foam** redistributes pressure, reduces shear forces and helps manage the skins' microclimate

# Pressure Redistribution

Dacapo Ampla allows for optimal immersion and envelopment of the patient. This translates to a lower tissue interface pressure for vulnerable skin and facilitates a lower risk of developing pressure injuries.

The pressure mapping test results (illustrated right) demonstrates the excellent pressure redistribution of Dacapo Ampla.

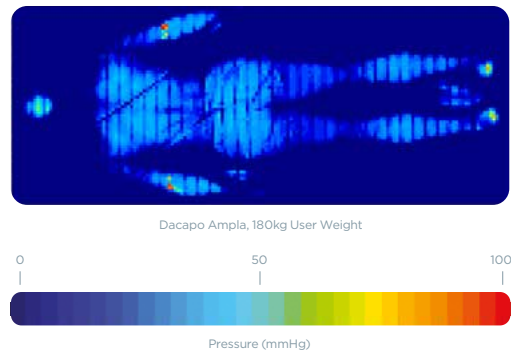


## Technical data ▼

For more information about this product, including the user manual, please visit our website: [www.invacare.eu](http://www.invacare.eu)



Width ▼

850, 900, 1030, 1150,  
1200 mm

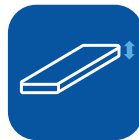

Height ▼

150 mm



Length ▼

2000, 2050 mm



Weight ▼

14 kg



User weight ▼

max. 350 kg



Clean temperature ▼

95 °C



Warranty ▼

2 years

**Tel: +44 61 487 70 70**  
[hqueurope@invacare.com](mailto:hqueurope@invacare.com)  
[www.invacare.eu.com](http://www.invacare.eu.com)

Invacare International GmbH  
 Neuhofweg 51  
 4147 Aesch BL  
 Switzerland



**Yes, you can.®**