

Medley[®] Ergo

More options. More choice.



A range of medical profiling beds
Trusted by healthcare professionals



Customised designs

The Medley Ergo range of medical beds fits perfectly at home and long-term care settings.

With more options, each design aims to achieve comfort, rest and sleep - while prioritising the needs of user and carer.

Choose from:

Medley Ergo

Medley Ergo Low



Safety

All Invacare beds meet IEC 60601-2-52 safety standards.



Quality

Our bed designs and specifications adhere to the highest standards.



Care-led

Each Invacare bed is designed with the knowledge of care settings in mind.



Comfort

Every Invacare bed promotes rest, dignity and independence.




Homely Design

Medley Ergo embodies comfort, warmth and personal touches.



* Medley Ergo shown with optional Select bed end panels and side rails. Invacare mattress sold separately.

More options. More choice.

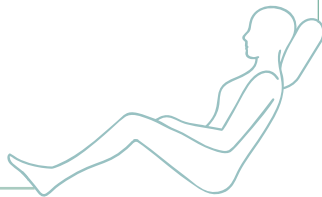

The Medley Ergo range is easy to use and handle as well as designed with options to maximise comfort and ensure practicality at the touch of a button.



Medley Ergo



Elevated leg section



Medley Ergo in tilt

01
Customised care
A range of options to meet various needs

11°
Supportive tilt
For positioning and pressure management.

Note the tilt function is only practised for users under medical supervision.

02
Quick & easy setup
Simple tool-free assembly.

03
Adjustable positioning
Backrest, thigh & legrest adjustment. Legrest elevates manually.

04
Height adjustable
Each model adjusts for comfortable working heights and optimal patient care.

Intuitive hand control
Clear visual display for ease of use. Optional lockable tilt available.

05
Homely design
Optional Select Beech wood bed end covers available to fit into any home.

08
Easy to clean
Suitable to wipe with mild detergents making it ideal for single or multiple use environments. See User Manual.

07
Interchangeable options
Healthcare environments can benefit from the versatility of using the range.

06
Range of side rails
Choose from Bella with integrated handles, Aria, Verso HC or Scala 2.

Medley Ergo shown with optional Select bed end covers and optional Bella side rails..

Medley® Ergo Low

Positioning the bed at a lower level increases security and reassurance for users and carers.
21 - 61 cm or 28 - 68 cm.



Medley® Ergo options

Customised care with a choice of options.



Transport trolley



Tech Support & After Sales

Ask our team. We are always here to help.

Tel: +41 61 487 70 70
support@invacare.com



Bella side rail in Beech wood



Aria side rail full length in steel



Scala 2 side rail collapsible in steel



Lifting pole with adjustable handle



Support handle 30 x 40 cm



Support handle 50 x 40 cm



Alternative handset with lockable tilt



Side rail cover for additional safety

Technical specifications



Width inside
90 cm



Width outside
102 cm



Length inside
208 cm



Length outside
224 cm



Backrest angle
0 - 70°



Tilt angle
11°



Height adjustment
40 - 80 cm
or 33 - 73cm



Heaviest part
17 kg



Mattress support
80 - 23 - 30
- 67 cm



Max user weight
145 kg



Total weight
64 kg



Thigh angle
0 - 33°

Medley Ergo Low

All technical details as above except:

Height adjustment:
21 - 61 cm
or 28 - 68 cm

Colours



Steel



Beech



Complies to the safety standard,
IEC 60601-2-52.



ISO 13485 certified

Tel: +41 61 487 70 70
globalcomms@invacare.com
www.invacare.eu.com

Invacare International GmbH
Neuhofweg 51
CH-4147 Aesch BL
SWITZERLAND

