



**matrx<sup>®</sup> PB**



**POSTURE BACK**

Classic Midline Support





# matrx® PB

## Clinician-Led Design

Introducing Matrx PB, our classic and popular approach to midline back support with our signature moulded foam.

### EXCEPTIONAL SUPPORT

Expertly contoured to provide increased midline postural support for people with moderate to high lateral trunk positioning needs.



#### PRECISION

Every specification developed by our clinician-led team aims for efficiency, quality and a precise fit.



#### COMFORT

Matrx signature moulded foam wrapped around the backshell for increased comfort and skin protection.



#### SUPPORT

Expertly contoured design to provide exceptional midline support.

## **m** Contoured for Full Support

Matrx® PB is a classic solution best suited to tilt and recline applications. Featuring full support across the shoulders and upper pelvis for people requiring additional torso support.




### DEPTH OPTIONS

13cm Contour




STANDARD

18cm Contour

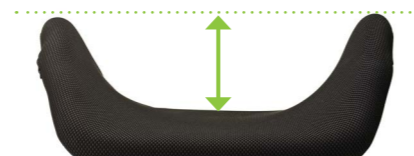

DEEP BACK

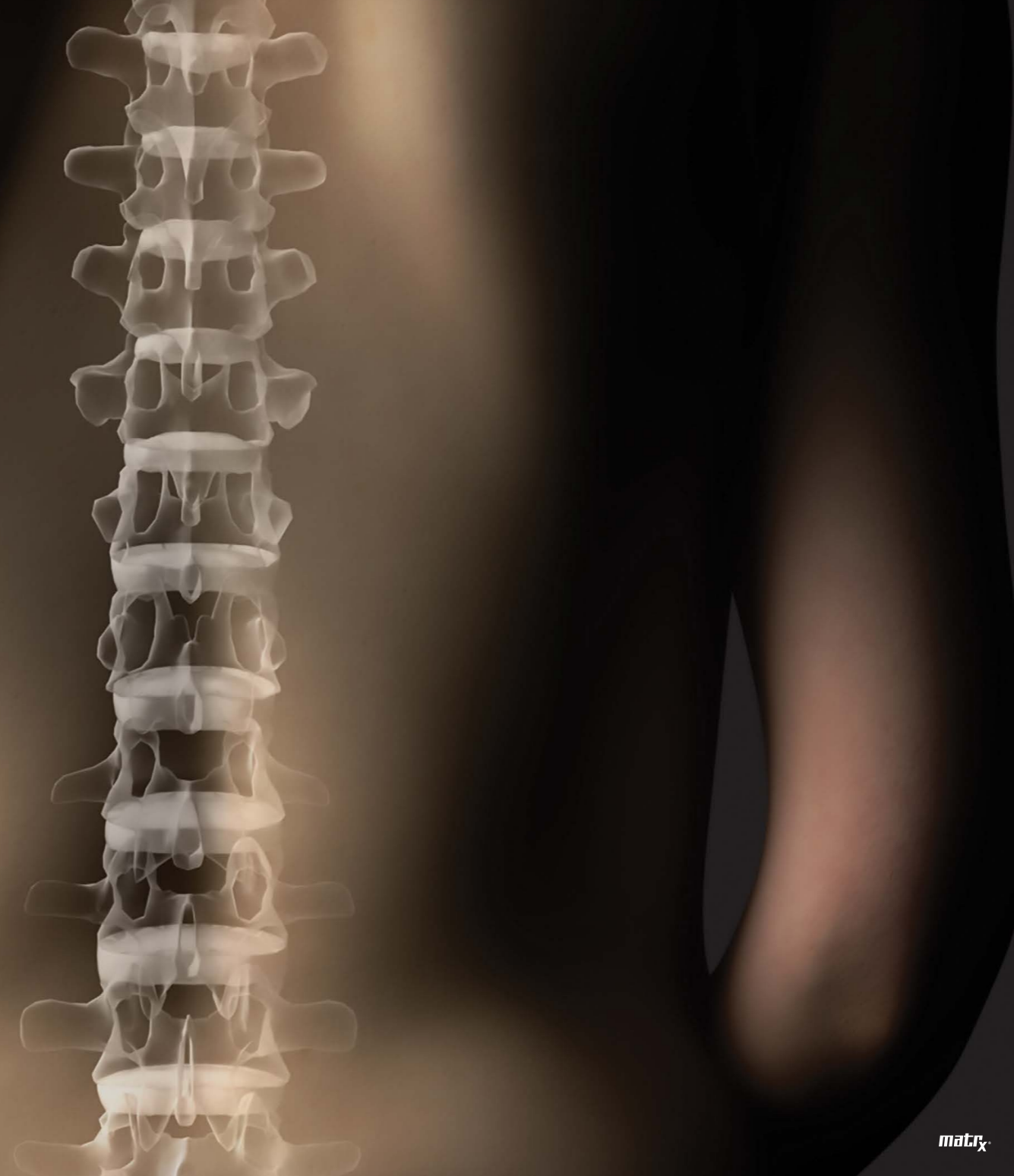



## **m** Take a Closer Look

Matrix® PB has a pronounced contoured back shell providing full lateral support for the upper body.

### **COVER**

Fully upholstered, moisture resistant and breathable

### **TWO DEPTHS**

Matrix PB 13cm or Matrix PB Deep 18cm options for required trunk support needs



### **BACK SHELL**

Lightweight ABS rigid back shell as standard or optional aluminium for heavy duty use

### **CONTOURED SHAPE**

Designed for comfortable midline support

### **MOULDED FOAM**


Wrapped around the shell for maximum comfort and skin protection. Made with Ultra-Fresh™ antimicrobial treatment for protection from odours and stains caused by bacteria, mould and mildew growth

## **M** Mounting & Hardware Options


The back supports are mounted onto the chair back canes and secured with the Easyset Mounting hardware. The flexible hardware is designed to match the needs of a range of body shapes and sizes, and can even be adjusted when the person is seated. It can also be positioned to avoid frame obstacles – such as arm supports and rigidizer bars etc.



### ADJUSTABLE WIDTH AND HEIGHT

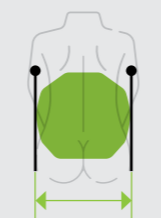

Fits chairs equal width to back model

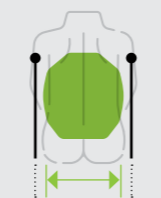

Fits chairs 2.5cm wider (or 5cm with optional extended pins) than back model

### STANDARD HARDWARE

#### EASYSET



- Simple and easy to install, and all adjustments can be made with the person sat in the wheelchair.
- Quick release for easy removal/portability.
- Highly adjustable in depth, angle (+/- 20°) and rotation.

### OPTIONAL HARDWARE

#### FIXED EASYSET



- Suitable for tilt/recline applications, manual chairs with limited hardware space or where back removal is not required.
- Permanent mount, requires tools for removal.

#### HEAVY-DUTY EASYSET



- Recommended for people with high tone or heavy intended use.

#### SPLIT



- Allows clearance of frame obstacles, quick release.







## Complete Range of Options



### PRIVACY FLAP

- Can be attached to the backrest with a pre-sewn loop strip located across the bottom rear of cushion cover.
- Available in 4 sizes.



### SHOULDER STRAP INTERFACE

For precise placement of Strap Guides.

## Back Sizes

		WIDTH				
		36cm	41cm	43cm	46cm	51cm
LENGTH	Fits W/C	36-41cm	41-46cm	43-48cm	46-51cm	51-56cm
	31cm	●	●	●	●	●
	41cm	●	●	●	●	●
	51cm	●	●	●	●	●

## Tech Data

- Angle adjustment ..... +/- 20°
- Depth adjustment ..... 10cm
- Cover cleaning ..... 60°C
- Max. User weight ..... 136kg std / 227kg HD





[invacare.eu.com](http://invacare.eu.com)

© 2025 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. EU – 01/25