

**matrx<sup>®</sup> MINI**



**MINI SERIES**

Supportive Paediatric Seating





## matrx<sup>®</sup> MINI

### Clinician-Led Design

The Matrx<sup>®</sup> Mini Series consists of a back support and cushion. It is designed specifically for children, providing exceptional support, comfort and maximum freedom of movement.

### LIGHTWEIGHT & RESILIENT

Helping children maximise their abilities. Lightweight, durable wheelchair seating system optimised for comfort and positioning.



#### PRECISION

The Mini Series has been carefully designed and constructed, with paediatrics in mind.



#### COMFORT

Ergonomically shaped for optimum support, immersion and comfort.



#### SUPPORT

Designed to promote good posture and prevent orthopedic changes.

## Optimum Positioning

Matrix Mini Series is comprised of our Matrix Mini Back and Matrix Mini PS cushion to give comfort, stability and improved positioning for children in their wheelchair.



## **M** Matrx® Mini PS

Matrx Mini PS is a high resilient wheelchair cushion made with our signature moulded foam. It is designed for mild positioning and targeted immersion around the ischial tuberosities.

### **MOULDED SHAPE**


Durable foam provides stability in the wheelchair. Includes Ultra-Fresh™ antimicrobial treatment for protection from odours and stains caused by bacteria, mould and mildew growth

### **ANATOMICALLY CONTOURED**

Maximised contact area for improved skin protection

### **WATERFALL FRONT EDGE**

Durable moulded high resilient foam provides stability in the wheelchair

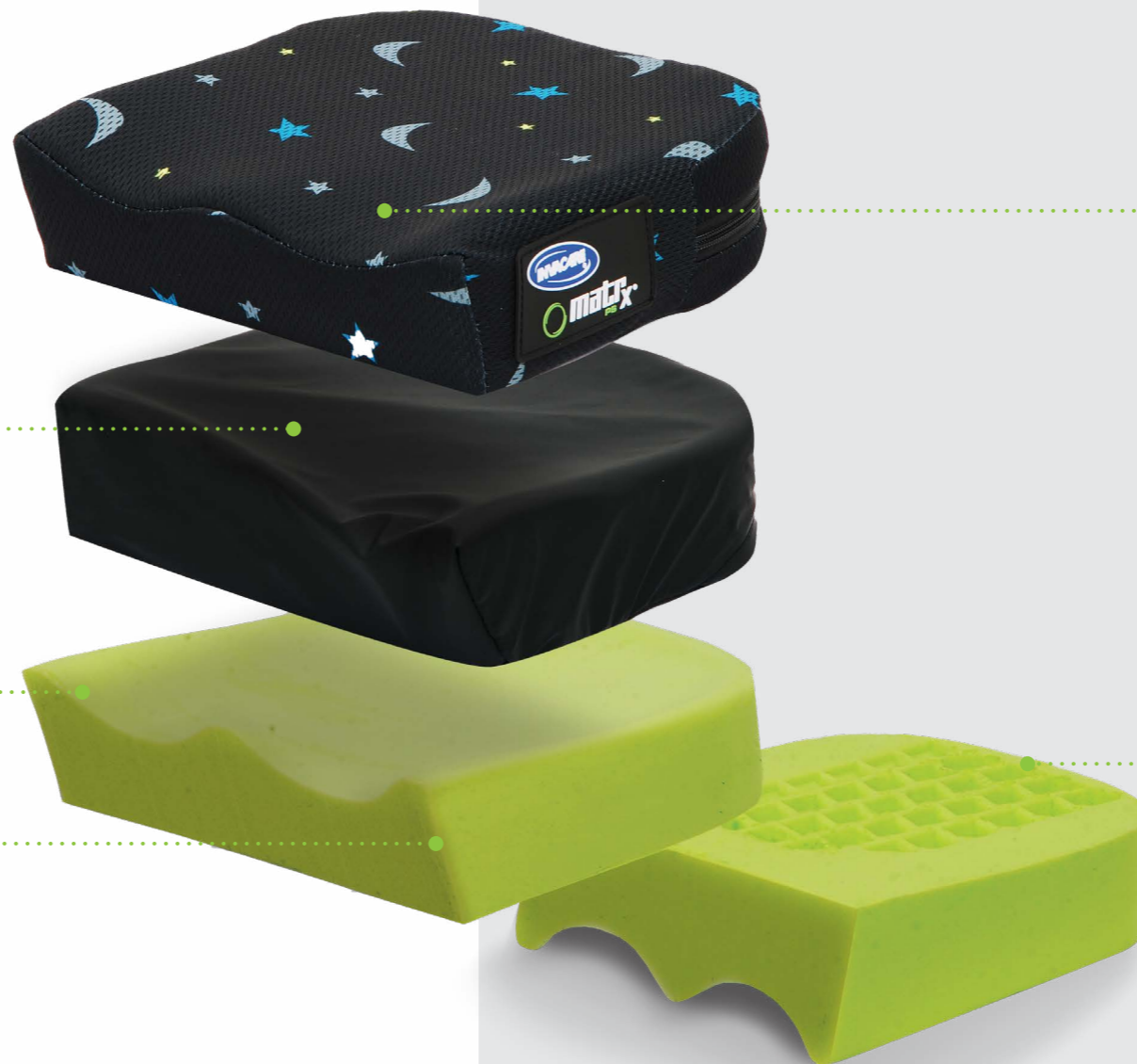

Action<sup>3</sup> Junior

IDEAL FOR ACTION JUNIOR WHEELCHAIRS

### **FUN DESIGN**

Moisture resistant and breathable cover with Kid\*a\*bra print

### **WAFFLED UNDERSIDE**

Provides ischial relief for improved immersion and comfort

## **M** Matrx® Mini Back

Matrx Mini Back together with the MiniSet hardware is a lightweight solution that can be installed and adjusted in minutes providing optimal positioning.

### SCAPULAR CUT OUTS

Allows freedom of upper extremity mobility for maximum function

### CONTOURED SHAPE

Designed to stabilise the hips and trunk to optimise positioning and function

### MINISET HARDWARE

Quick release hardware allows adjustments to accommodate growth



Action<sup>3</sup> Junior

IDEAL FOR ACTION JUNIOR WHEELCHAIRS

### FUN DESIGN

Moisture resistant and breathable cover with Kid\*a\*bra print

### LIGHTWEIGHT BACKSHELL


Durable aluminium, anatomically contoured

## **M** Mounting & Hardware Options


Adjustable pins secure Matrix Mini Back to the MiniSet hardware mounted on the wheelchair. Hardware can be positioned to avoid frame obstacles such as arm supports and can be adjusted while the child is seated.



### ADJUSTABLE WIDTH AND HEIGHT

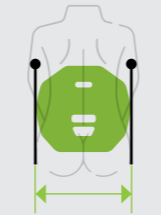

Fits chairs equal width to back model

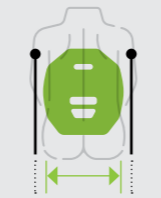

Fits chairs 2.5cm wider. Optional Extended mounting pins available to fit chairs up to 5cm wider

### MINI ELAN HEADREST

Infinitely adjustable through 6 adjustment points. A full range of Elan positioning pads are available.



### HEADREST

Mini backs are pre-drilled to allow attachment of headrests.



### MINI FIXED LATERALS

Height adjustable.



### CHEST HARNESS INTERFACE

Allows attachment of chest support straps to mini backs.



### MINISET HARDWARE

100mm depth adjustment for growth and easy to reposition with +/- 20° recline and rotation. Fits 25mm, 22mm & 19mm diameter back canes.



### MINI SWING-AWAY LATERALS

Height, width and depth adjustable. Latch release mechanism allows pads to be rotated out of the way for side transfers.



### CHEST HARNESS FOR MINI BACK

Allows attachment of chest support straps to mini backs.








## Back sizes

		WIDTH				
		26cm	28cm	31cm	33cm	36cm
LENGTH	Fits W/C	26-31cm	28-33cm	31-36cm	33-38cm	36-41cm
	18cm	●	●	●	●	●
	26cm	●	●	●	●	●
	31cm	●	●	●	●	●
	36cm	●	●	●	●	●







## Cushion sizes

		WIDTH				
		26cm	28cm	31cm	33cm	36cm
DEPTH	26cm	●	●	●	●	●
	31cm	●	●	●	●	●
	36cm	●	●	●	●	●

## Back tech data

-  Contour depth ..... 8cm
-  Depth adjustment ..... 10cm
-  Angle adjustment ..... +/- 20°
-  Total product weight ..... 400g and up
-  Max. user weight ..... 68kg

## Cushion tech data

-  Thickness at Abductor/adductor ... 8.9cm
-  Thickness at leg trough ..... 5cm
-  Total product weight ..... 150g and up
-  Max. user weight ..... 68kg
-  Max. cleaning temp(Outer cover) ... 60°C
-  Surface cleaning  
(Inner liner) ..... with damp cloth





[invacare.eu.com](http://invacare.eu.com)

© 2025 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. EU – 01/25